



ConvaTec

Ostomy Care

2022

OSTOMY

Digital Product Guide

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)


HELPING TO CREATE HEALTHY BONDS



Together, we can help your patients create a healthy bond with their ostomy.

ConvaTec's skin-friendly and secure pouch connections empower your patients with the security to go out, live life, and connect with others. In addition to providing a comprehensive portfolio of ostomy product solutions, we are committed to elevating the experience of living with an ostomy through our me+™ program, which provides your patients with the education, support and resources they need.

Because, when your patient has a healthy bond with their ostomy, they have a healthy bond to their life.

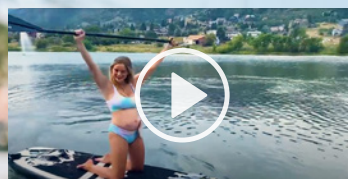

WATCH NOW

to see how ConvaTec is creating Healthy Bonds

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

TABLE OF CONTENTS

me+™ SUPPORT	1
• me+™ Program Overview	3
• me+™ Virtual Support	5
• me+™ Style	7
EDUCATION	9
<i>Ostomy Overview</i>	<i>11</i>
• Introduction to Ostomy Care	11
• Choosing the Right Ostomy System	13
<i>Guide to Ostomy Systems</i>	<i>15</i>
• A Guide to Adhesives	15
• A Guide to Skin Barriers	19
• Moldable Skin Barriers	21
• Accordion Skin Barriers	23
• The Purpose of Convexity	25
• Ostomy Pouch Options	29
• Common Pouch Features	31
TWO-PIECE SYSTEMS	33
• Natura™ Durahesive® Skin Barriers	35
• Natura™ Stomahesive® Skin Barriers	39
• Natura™ + Pouches	42
• Esteem Synergy® Skin Barriers	49
• Esteem Synergy® Pouches	51
ONE-PIECE SYSTEMS	53
• Esteem™ + Pouches	55
• ActiveLife® Pouches	65

ACUTE KITS	69
• Two-Piece Kits	71
• One-Piece Kits	79
• Unit-Dose Kits	80
PEDIATRIC PRODUCTS	81
• Extra Small Systems	83
• Two-Piece Systems	85
• One-Piece Systems	88
ACCESSORIES	89
• Skin Preparation	91
• Pouching Support	93
• Added Security	97
• Pouch Removal	99
• Irrigation	101
WOUND & FISTULA POUCHES	103
RESOURCES	107
• Discontinued Products	109
• Patient Ostomy Resources & Support	117
• Quick Links Guide	118
• Medicare Utilization Chart	119
• Product Code Index	120
• References	126

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

THERE IS A NEED FOR ONGOING PATIENT SUPPORT

Most common ostomy self-care challenges are related to *leakage, skin issues, frequency of changes and difficulty caring for an ostomy*¹

Following an interaction with a **me+™ WOC nurse:**

98%

of survey respondents said they would *recommend this service*²

97%

of those surveyed had a better understanding of how to *care for their ostomy* after a me+™ virtual visit (telehealth)³

90%

of those surveyed said they *avoided an additional visit* to a healthcare facility²



THE me+™ PROGRAM SUPPORTS PATIENTS

The me+™ program is a comprehensive support program designed to complement clinical care by providing patients with support, inspiration, insights and ideas throughout their entire ostomy journey.


me+™ program members can:

- ✓ Request free ostomy samples
- ✓ Speak to a WOC nurse or product specialist
- ✓ Receive complimentary virtual telehealth services
- ✓ Access online resources, education and guided recovery tools
- ✓ Connect to our Language Line Services, which offers translation support in over 240 foreign languages and dialects – at bedside or home

Health Care Professionals

TO ENROLL YOUR PATIENTS in the me+ program

call us at 877-585-0470 or email concierge@convatec.com

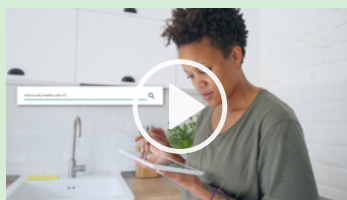

LEARN HOW

our me+™ program can provide the right ostomy support for your patient

TABLE OF CONTENTS

me+™ SUPPORT

me+™ SUPPORT

- me+™ Program Overview
- me+™ Virtual Support
- me+™ Style

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)



Our goal is to *partner with you and your patients to deliver best-in-class care* throughout the patient journey:

✓ Pre-Surgery / Surgery

- Access online resources at mepluscare.com
- Request me+™ patient Starter Kit and education materials
- Enroll in me+ program
- Learn about the me+ recovery series, developed to aid rehabilitation after surgery

✓ Weeks Immediately After Surgery

- Consult with me+ WOC nurses and product specialists
- Personalized Starter Kit
- me+ welcome email and support
- me+ recovery phase 1 exercise booklet*

✓ First Year & Beyond

- Explore online resources and tools
- Get inspired by the me+ community
- Sample new products
- Continue me+ recovery series exercises*

Throughout their journey, your patients can call 1-800-422-8811 for live support from me+™ WOC nurses and product specialists

A me+™ Concierge team member can assist Healthcare professionals with:

- ✓ Personalized support
- ✓ Enrolling your patients into the me+™ program
- ✓ Providing a me+ Starter Kit and free product samples for your patients

TO LEARN MORE contact your me+ concierge directly at 877-585-0470

OSTOMY CARE *has met its* APP



KEY FEATURES

- ✓ Nurse Login
- ✓ Product Selector
- ✓ Ship Samples
- ✓ Saved Kits
- ✓ me+™ Virtual Support
- ✓ Enroll Patients in me+™

Scan the code or use the buttons to download the ConvaTec Ostomy Nurse Solutions App



SCAN NOW



OstomyNurseSolutions.com

TABLE OF CONTENTS

me+™ SUPPORT

me+™ SUPPORT

- me+™ Program Overview
- me+™ Virtual Support
- me+™ Style

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

*The information in the series is not intended to provide medical advice. Patients and users should speak with their health care professional before engaging in physical activity following ostomy surgery. In the event they experience pain or discomfort, they should stop the activity immediately and contact their health care professional.



me+™ OFFERS VIRTUAL TELEHEALTH NATIONWIDE

The me+™ program is the first of its kind to offer Telehealth in all 50 states. Our complimentary virtual telehealth services provides visual support and live demonstrations through our telehealth partner, Corstrata™.

Interested in learning more about the me+™ virtual support option?



WATCH THE VIDEO

or call us today at
1-800-422-8811 to see
if a telehealth consultation
is right for your patient.

*"It's truly allowed me to have a more positive outlook,
which was certainly suffering before the visit."*

- Michael, me+™ member and telehealth user

me+™ virtual support is not intended to be a substitute for professional medical advice, diagnosis, or treatment. Always seek the advice of your physician or other qualified health provider with any questions you may have regarding a medical condition.



Telehealth visits will be managed through ConvaTec partner, Corstrata™. Only visits scheduled through a ConvaTec referral are complimentary.



You Can Feel Confident Recommending me+™ Virtual Support Because It Is:

- ✓ Provided by Board-Certified Ostomy Nurses
- ✓ Available in all 50 states
- ✓ HIPAA Compliant, Private & Secure
- ✓ Available in English and Spanish

TABLE OF CONTENTS

me+™ SUPPORT

me+™ SUPPORT

- me+™ Program Overview
- me+™ Virtual Support
- me+™ Style

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)



LOOK GREAT AND FEEL CONFIDENT

Ostomysecrets® offers a full line of stylish and functional undergarments, swimwear and skin care designed to support your ostomy system.

Who says living with an ostomy can't be stylish, sleek and even sexy? With me+™ style, you don't have to change what you wear or hide your ostomy pouch under bulky, unflattering clothes. Ostomysecrets.com is a one-of-a-kind clothing site created by people with ostomies, for people with ostomies. From swimwear to women's thong underwear to men's athletic boxers, it's a full line of apparel you won't find anywhere else. Ostomysecrets.com gives you everything you need to conceal and support your ostomy pouch, without sacrificing your personal style.

HAVE QUESTIONS?

Call an Ostomysecrets® specialist at 1-877-613-6246, chat online at ostomysecrets.com or send us an email at info@ostomysecrets.com.

SHOP NOW

and **SAVE 15% OFF YOUR FIRST ORDER** with code **STYLE15**.

TABLE OF CONTENTS

me+™ SUPPORT

me+™ SUPPORT

- me+™ Program Overview
- me+™ Virtual Support
- me+™ Style

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

EDUCATION

Feel more confident supporting your ostomy patients.

Take advantage of our educational resources to learn about ostomy basics, importance of skin care, and more to help you elevate the patient experience.

- ✓ An Introduction to Ostomy Care
- ✓ Choosing the Right Ostomy System
- ✓ A Guide to Adhesives
- ✓ A Guide to Skin Barriers
- ✓ Moldable Skin Barriers
- ✓ Accordion Skin Barriers
- ✓ The Purpose of Convexity
- ✓ The Characteristics of Convexity
- ✓ Ostomy Pouch Options
- ✓ Common Pouch Features

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care
- Choosing the Right Ostomy System
- A Guide to Adhesives
- A Guide to Skin Barriers
- Moldable Skin Barriers
- Accordion Skin Barriers
- The Purpose of Convexity
- The Characteristics of Convexity
- Ostomy Pouch Options
- Common Pouch Features

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

An Introduction to OSTOMY CARE

There are 3 primary types of ostomy surgeries. Understanding the fundamental difference is important to select the right product and provide the right support.



TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care

TWO-PIECE SYSTEMS

- Choosing the Right Ostomy System

- A Guide to Adhesives

ONE-PIECE SYSTEMS

- A Guide to Skin Barriers

- Moldable Skin Barriers

- Accordion Skin Barriers

ACUTE KITS

- The Purpose of Convexity

- The Characteristics of Convexity

PEDIATRIC PRODUCTS

- Ostomy Pouch Options

- Common Pouch Features

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

COLOSTOMY (large intestine)

A colostomy is a surgical procedure that creates an opening on the abdomen through which the large intestine is redirected to bypass diseased or damaged parts of the large intestine. The large intestine is brought to the surface of the skin, and stitched to the abdomen to form a new exit for the stool from the bowel.

Things to know about a Colostomy

- ✓ Primarily formed due to colon or rectal cancer or trauma
- ✓ Can be temporary or permanent
- ✓ Stool tends to be semi-formed to formed and predictable 1-2 times per day

ILEOSTOMY (small intestine)

An ileostomy is a surgical procedure that creates an opening on the abdomen through which the end of your small intestine is redirected to bypass diseased or damaged parts of the digestive system. The small intestine is brought to the surface of the skin, and stitched to your abdomen to form a new exit for waste from your bowels.

Things to know about an Ileostomy

- ✓ Primarily formed due to inflammatory bowel disease (Crohn's, Ulcerative Colitis)
- ✓ Can be temporary or permanent
- ✓ Stool is liquid or pasty and constant throughout the day
- ✓ The small intestine contains digestive enzymes and acids, and the output may cause skin irritation if it comes in contact with the skin

UROSTOMY (urinary tract)

A urostomy is a surgical procedure that creates an opening on the abdomen through which urine is redirected away from a bladder that is damaged or diseased. After the bladder is removed, a piece of small intestine (ileum) is separated from the rest of the bowel. This piece is used as a conduit for urine to pass out of the body. The ureters are attached to this new passageway and one end is brought up to an opening made on your abdomen.

Things to know about a Urostomy

- ✓ Primarily formed due to cancer or inflammation of the bladder
- ✓ Is always permanent
- ✓ Urine flows through the stoma as it is produced which means output is constant
- ✓ Urine may cause skin irritation if it comes in contact with the skin

Choosing the Right OSTOMY SYSTEM

The ostomy pouch collects stoma output. The pouch comes as part of a two-piece or one-piece system.



2 PIECE SYSTEM



The two-piece system includes a skin barrier and a separate pouch that securely join together. This option allows the pouch to be removed to empty or change, without removing the skin barrier as frequently.

Designed for patients

- ✓ who want to change their pouch more frequently than their skin barrier
- ✓ with an ileostomy who need to empty/change the pouch frequently
- ✓ who exercise, travel and play sports

1 PIECE SYSTEM



The skin barrier and pouch are provided as an integrated solution that cannot be separated. It offers one simple and convenient step to attach or remove, providing a fresh start every time. This option is thin and flexible, and designed for more frequent system changes.

Designed for patients

- ✓ post-operatively – removes any added pressure on the abdomen caused when connecting a skin barrier and pouch
- ✓ who prefer to change their ostomy system frequently
- ✓ who prefer a flexible and low profile ostomy system

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care

TWO-PIECE SYSTEMS

• Choosing the Right Ostomy System

- A Guide to Adhesives
- A Guide to Skin Barriers
- Moldable Skin Barriers
- Accordion Skin Barriers

ONE-PIECE SYSTEMS

- The Purpose of Convexity

ACUTE KITS

- The Characteristics of Convexity

PEDIATRIC PRODUCTS

- Ostomy Pouch Options
- Common Pouch Features

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

A Guide to ADHESIVES

The adhesive is an important component of the skin barrier. Our product portfolio features a range of adhesive technologies to enable you to select an appropriate solution that meets the needs of your patients.



STOMAHESIVE® (standard wear)

Designed for people with mostly solid output, our unique Stomahesive® technology is a hydrocolloid adhesive that gently adheres to the skin, protecting it from the damaging effects of stoma effluent.

We also offer a Modified Stomahesive® option, which features a touch of mineral oil that makes the barrier more flexible during wear, and removal easier on the skin.

Designed for patients

- ✓ with semi-formed/formed output
- ✓ who have sensitive skin
- ✓ with frequent changes
- ✓ who prefer a standard wear time

DURAHESIVE® (extended wear)

Designed for people with mostly liquid output, our unique Durahesive® technology is a hydrocolloid adhesive that swells up upon hydration to protect the stoma. This special effect called turtlenecking, helps create a more secure seal without harming the skin.

We also offer Durahesive® Plus, which features an enhanced tack, strength and cohesion designed for extended wear time or aggressive output.

Designed for patients

- ✓ with semi-liquid/liquid output, including urine
- ✓ who want or need longer wear time

ESTEEM™ + Soft & Flex CONVEX

Adhesive designed with a high initial tack, to aid in the initial attachment to the skin. The hydrocolloid barrier features a thin tapered edge to help prevent the edge catching on clothing, and a thicker center to help absorption around the stoma.⁴

Designed for patients

- ✓ with both solid and liquid output
- ✓ who prefer a standard wear time

TAPE COLLARS
are available in two forms:

Acrylic: Tape collar designed with good tack
Hydrocolloid: Designed for those with sensitivity to other tape collars

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care
- Choosing the Right Ostomy System
- **A Guide to Adhesives**
- A Guide to Skin Barriers
- Moldable Skin Barriers
- Accordion Skin Barriers
- The Purpose of Convexity
- The Characteristics of Convexity
- Ostomy Pouch Options
- Common Pouch Features

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

A Guide to ADHESIVES




ABOUT OUR ADHESIVES

ConvaTec's proprietary adhesives are integral in all of our ostomy skin barriers. Our hydrocolloid adhesives are designed for use in high-moisture environments and have decades of proven product performance.

The two charts below will help determine which adhesive may be best for your patient.

Adhesives Feature Guide

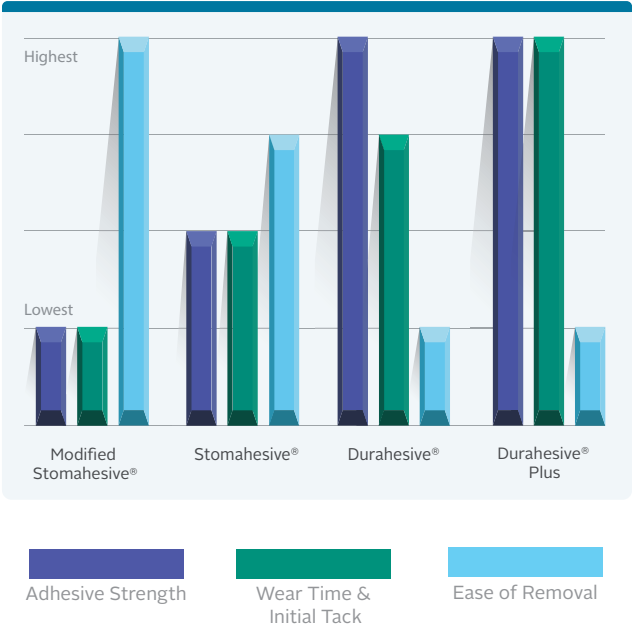

Our formulations are based on:

- Over four decades of proven and trusted ostomy solutions
- Leading adhesive technology, formulated after 15 years of research and with over 400 formulations
- Unmatched support designed to dramatically improve life with an ostomy

Adhesives Selection Guide

Adhesive Type	Formed Stool	Semi-Formed Stool	Semi-Liquid Stool	Liquid Stool	Urine	System Type	
Durahesive® Plus			✓	✓	✓	1-Piece Only	Extended Wear
Durahesive®		✓	✓	✓	✓	1-Piece and 2 Piece	
Esteem™+ Flex/Soft Convex Line	✓	✓	✓	✓	✓	1-Piece Only	Standard Wear
Stomahesive®	✓	✓	✓	✓		1-Piece and 2-Piece	
Modified Stomahesive®	✓	✓				1-Piece Only	

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care

TWO-PIECE SYSTEMS

- Choosing the Right Ostomy System

A Guide to Adhesives

ONE-PIECE SYSTEMS

- A Guide to Skin Barriers

- Moldable Skin Barriers

- Accordion Skin Barriers

ACUTE KITS

- The Purpose of Convexity

- The Characteristics of Convexity

PEDIATRIC PRODUCTS

- Ostomy Pouch Options

- Common Pouch Features

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[LEARN MORE](#) about ConvaTec Adhesives


[Request Product Sample](#)

A Guide to SKIN BARRIERS

The skin barrier is the piece of the pouching system that sticks to the body. It should fit snugly around the stoma to help prevent leakage and have an opening customized to stoma size and shape.



MOLDABLE

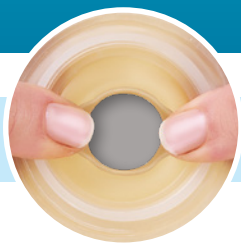

ConvaTec Moldable Technology™ skin barrier with Rebounding Memory Technology™ creates an elastic-like seal that rebounds to fit most shapes and sizes, and then “turtlenecks” upward, to enhance the seal and protect the skin.⁵

Designed for patients

- ✓ who have large, round, oval or irregular shaped stomas
- ✓ with dexterity issues who find cutting challenging

CUT-TO-FIT




Cut-to-Fit barriers can be cut to match the shape and size of individual stomas. Using scissors, one can cut a circle, oval or any shape needed to fit the stoma.

Designed for patients

- ✓ with a stoma size in-between pre-cut sizes
- ✓ with an irregular stoma shape

PRE-CUT

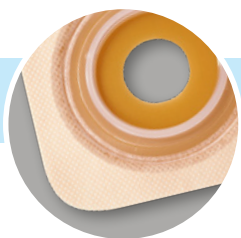

Pre-cut barriers have a uniform, circular stoma hole and are ready to apply right out of the box. They are available for the most common stoma sizes.

Designed for patients

- ✓ with an evenly round stoma that fits the available pre-cut size
- ✓ with dexterity issues who find cutting challenging

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care

TWO-PIECE SYSTEMS

- Choosing the Right Ostomy System

- A Guide to Adhesives

• A Guide to Skin Barriers

- Moldable Skin Barriers
- Accordion Skin Barriers

ONE-PIECE SYSTEMS

- The Purpose of Convexity

ACUTE KITS

- The Characteristics of Convexity

PEDIATRIC PRODUCTS

- Ostomy Pouch Options

- Common Pouch Features

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

MOLDABLE

Skin Barriers

ConvaTec Moldable Technology™ is the **LEADING MOLDABLE SKIN BARRIER** that new patients wear leaving the hospital⁶. The skin barrier is designed to hug the stoma for a more custom and secure fit.

97% strongly agree that it is *easier to apply* Moldable than Cut-to-Fit⁷

100% strongly agree that it is *faster to apply* Moldable than Cut-to-Fit⁷

CLINICALLY PROVEN


to prevent and improve peristomal skin complications⁵

Up to 55% of ostomates have skin issues⁸

- **95% of people** who started on ConvaTec Moldable Technology™ kept their skin healthy⁵
- **86% of people** who switched to ConvaTec Moldable Technology improved their skin condition⁵


See our Moldable Technology in action and

request a FREE sample now

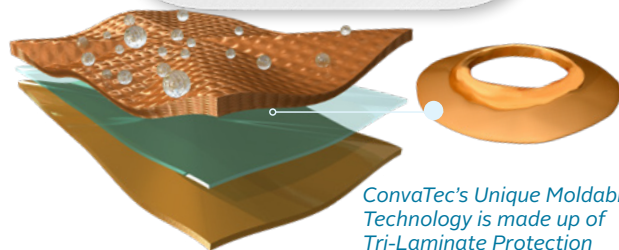


LEAK-LOCK™

First-of-its kind system provides a snug “turtleneck” effect for maximum leakage control and added security

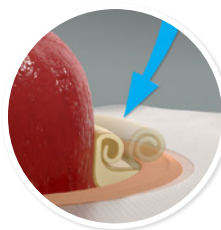

PATENTED REBOUNDED MEMORY TECHNOLOGY™

Customizable so that it molds to the stoma for a perfect fit every time

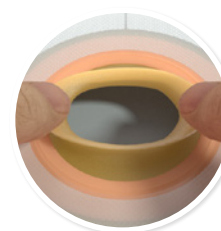

STAYS ON

Roll, mold, and hold to the stoma's shape and size. Barriers are reinforced with our acclaimed adhesives, to keep it secure and give your patients peace of mind

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care

TWO-PIECE SYSTEMS

- Choosing the Right Ostomy System
- A Guide to Adhesives

ONE-PIECE SYSTEMS

- A Guide to Skin Barriers
- **Moldable Skin Barriers**

ACUTE KITS

- Accordion Skin Barriers
- The Purpose of Convexity
- The Characteristics of Convexity

PEDIATRIC PRODUCTS

- Ostomy Pouch Options
- Common Pouch Features

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

ACCORDION

Skin Barriers

The Accordion Flange is an option on a range of Natura® skin barriers. The accordion flange lifts to provide generous finger room so coupling is easy and comfortable.



PRODUCT FEATURES

- 2X more finger space than any other lifting flange⁹
- Simple to measure and cut or mold
- Range of flat, convex, Cut-to-Fit and Moldable options
- Compatible with SUR-FIT® Natura™ and Natura™+ pouch option

Designed for patients

- ✓ with dexterity issues as it provides finger room so coupling is easy and comfortable
- ✓ designed to prevent pressure or trauma if abdominal tenderness is present due to post-operative stoma surgery/incisions

HOW IT WORKS

- Designed to help minimize abdominal pressure and pain when attaching the pouch.
- Connects with an audible 'click' and physical 'snap' to confirm the pouch is firmly attached.

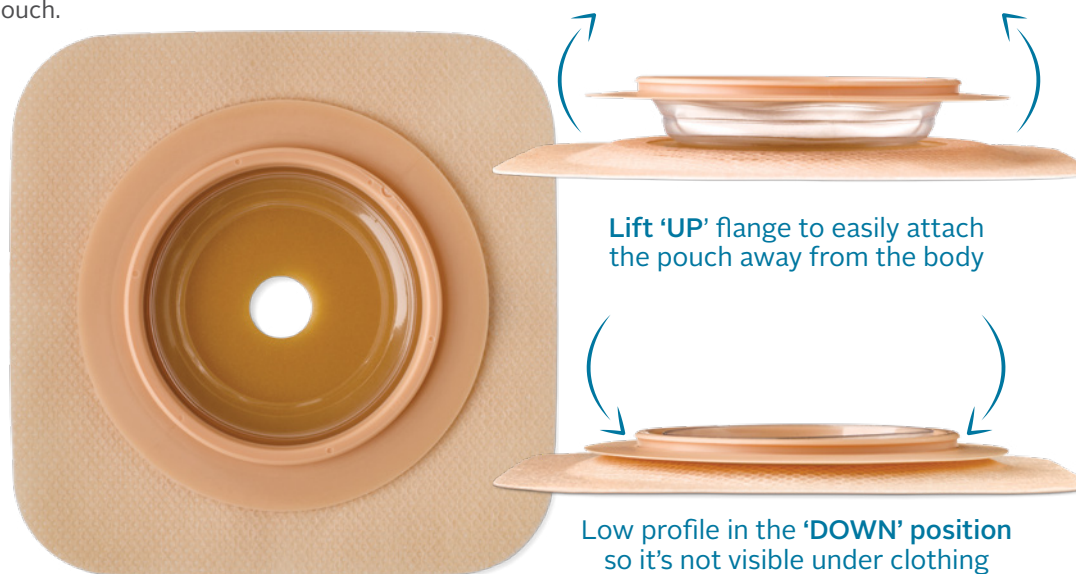


TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care

TWO-PIECE SYSTEMS

- Choosing the Right Ostomy System

- A Guide to Adhesives

ONE-PIECE SYSTEMS

- A Guide to Skin Barriers

- Moldable Skin Barriers

- **Accordion Skin Barriers**

ACUTE KITS

- The Purpose of Convexity

- The Characteristics of Convexity

PEDIATRIC PRODUCTS

- Ostomy Pouch Options

- Common Pouch Features

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

The Purpose of CONVEXITY



WHAT IS CONVEXITY?

All bodies are different. Convex wafers are designed to help manage irregular stomas and abdomens to help reduce leakage and simplify stoma management.

- Convexity is the outward curving of the baseplate or skin barrier.^{10,11}
- Convexity applies a uniform pressure to the immediate peristomal skin, improving the projection of the stoma for drainage.¹²
- Convexity helps maintain a secure seal between the pouch and the skin.¹³

GROWING CONVEXITY NEED

Data shows the need for convex skin barriers growing year on year.¹⁴ This is due to a number of factors including:


- Changing surgical techniques/outcomes of stoma shape & profile
- Increased obesity rates

WHEN TO USE CONVEXITY



Flush Stomas

Located at the same level as the abdomen skin surface, and does not protrude¹⁵

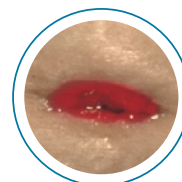

Retracted Stomas

Stoma pulled into the abdomen below skin level or having little or no spout¹⁶



Skin Folds

A stoma which is deep in a skin fold/crease¹³



Soft or Flaccid Abdomens

An abdomen requiring firm level of support to stabilize the pouching system¹⁵

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care

TWO-PIECE SYSTEMS

- Choosing the Right Ostomy System

- A Guide to Adhesives

ONE-PIECE SYSTEMS

- A Guide to Skin Barriers

- Moldable Skin Barriers

- Accordion Skin Barriers

ACUTE KITS

- **The Purpose of Convexity**

- The Characteristics of Convexity

PEDIATRIC PRODUCTS

- Ostomy Pouch Options

- Common Pouch Features

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

The Characteristics of **CONVEXITY**



RIGIDITY

Firm/Rigid Convexity

Typically used on soft or flaccid abdomens.


Soft Convexity

Typically used on firm or toned abdomens.

DEPTH


Height of the convex plateau.

✓ Deep ✓ Shallow

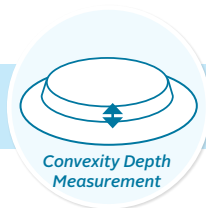

PLATEAU WIDTH

Distance across the raised flat surface from edge-to-edge

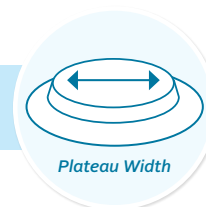


Why Plateau Width Matters

Successful convexity must press down on the peristomal skin enough to make a good seal around the base of the stoma and help with protrusion.¹⁴

- ✓ **Convexity Pressure:** A narrow plateau width may be used to create pressure directly around the stoma with lesser pressure as the plateau width gets wider.
- ✓ **Plateau Width:** A wider plateau width may be used to flatten out skin creases away from the stoma.

Plateau Width Guidance for Esteem™ + *Flex* and *Soft Convex*

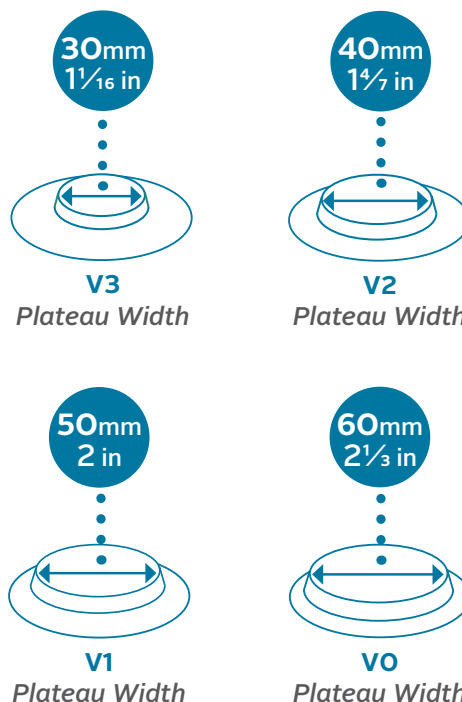

Choosing the right convexity solution for different body profiles can be challenging. Please feel free to consult with one of our me+™ Specialists or WOC Nurses for further guidance.

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care

TWO-PIECE SYSTEMS

- Choosing the Right Ostomy System

ONE-PIECE SYSTEMS

- A Guide to Adhesives
- A Guide to Skin Barriers
- Moldable Skin Barriers
- Accordion Skin Barriers

ACUTE KITS

- The Purpose of Convexity

PEDIATRIC PRODUCTS

- **The Characteristics of Convexity**
- Ostomy Pouch Options
- Common Pouch Features

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

Ostomy Pouch OPTIONS

ConvaTec offers a wide range of pouching options across both 2 Piece and 1 Piece Systems.



CLOSED

A closed pouch doesn't require draining. Once they are full, you simply remove and replace with a new pouch.

Additionally, ideal during physical activity, intimate moments, swimming, and travel.



Designed for patients

- ✓ who have thick, formed stool that is predictable (1-2 times per day)

DRAINABLE

Drainable pouches allow you to drain the output periodically throughout the day so that you do not have to change pouches as frequently.



Designed for patients

- ✓ who have loose, liquid stool or semi-formed stool that can be unpredictable or constant throughout the day

UROSTOMY

A Urostomy pouch has a faucet tap, soft tap, or fold over tap at the bottom that allows for draining urine as frequently as needed.



Designed for patients

- ✓ with a Urostomy

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care
- Choosing the Right Ostomy System
- A Guide to Adhesives
- A Guide to Skin Barriers
- Moldable Skin Barriers
- Accordion Skin Barriers
- The Purpose of Convexity
- The Characteristics of Convexity
- **Ostomy Pouch Options**
- Common Pouch Features

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

Common POUCH FEATURES



POUCH COLOR



Transparent (*Clear front of pouch*)

Clear pouch front for easy output assessment and determining when a pouch needs to be changed or emptied. Soft and quiet material on the body side of the pouch.

Opaque (*Tan-colored pouch*)

Soft and quiet material cover both the front and the body side of the pouch for more discretion and comfort.



Easy-View Window

Some options come with an easy-view window for accurate positioning and observation.

FILTER or NO FILTER



Filters help release and deodorize gas and keep the closed or drainable pouch flat. Note that the filter must be covered with a filter cover while swimming, bathing or showering.

DRAINABLE POUCH CLOSURE



InvisiClose™ Clipless Tail Closure

Self closing drainable pouch closure which provides an easy closure option that does not require a separate clip for added security.



Lock-it Pocket™

Folds inside out to hide away the secured drainage tail—providing the flexibility of a drainable pouch, with the feel and comfort of a closed-end pouch.

Traditional drainable pouches are closed with a separate clip.

UROSTOMY POUCH DRAINAGE



The Accuseal™ Tap with valve is a faucet type tap designed for easy manipulation and emptying. Other drainage options include a Soft Tap or Fold Over Tap.

All Urostomy pouches include an anti-reflux valve which is designed to help prevent backflow.

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

EDUCATION

- An Introduction to Ostomy Care
- Choosing the Right Ostomy System
- A Guide to Adhesives
- A Guide to Skin Barriers
- Moldable Skin Barriers
- Accordion Skin Barriers
- The Purpose of Convexity
- The Characteristics of Convexity
- Ostomy Pouch Options
- **Common Pouch Features**

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

TWO-PIECE


Systems



NATURA™

The Natura™ brand (inclusive of Natura™, Sur-Fit Natura™, and Natura™+) is an extensive range of skin barriers and pouches that connect together with a plastic ring called a flange in a simple snap-off, snap-on design.



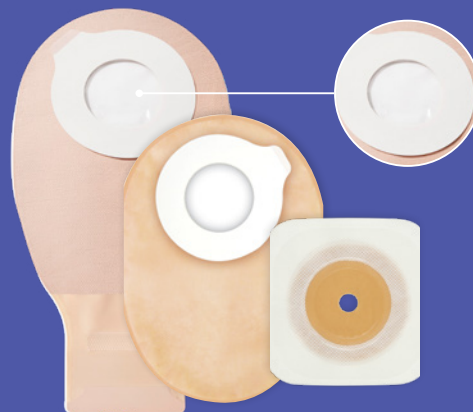



Two-piece systems include a **skin barrier** and a **separate pouch** that securely join together. There are two brands that make up ConvaTec's Two-Piece Ostomy Systems: Natura™ + and Esteem Synergy® +.






ESTEEM SYNERGY®

The Esteem Synergy®+ brand are skin barriers and pouches that connect together with Adhesive Coupling Technology™. This is a flat, flexible adhesive foam ring on the pouch that sticks to a flat 'landing zone' on the skin barrier.

Color Key & Sizing Guide

COLOR MATCH						
FLANGE SIZE	1 1/4" (32mm)	1 1/2" (38mm)	1 3/4" (45mm)	2 1/4" (57mm)	2 3/4" (70mm)	4" (100mm)

Color Key & Sizing Guide

COLOR MATCH			
COUPLING SIZE	Small Stoma Opening	Medium Stoma Opening	Large Stoma Opening



When ordering a Two-Piece System, **make sure you select the same flange/coupling size for the skin barrier and pouch.** Use the colored dots in the table to help.

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

TWO-PIECE SYSTEMS

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- Natura™ + Pouches
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

NATURA™ DURAHEESIVE® SKIN BARRIERS

Compatible with Natura™ + and SUR-FIT® Natura™ pouches



SUR-FIT® NATURA™ SKIN BARRIER

Cut-to-Fit

- Durahesive® (Extended Wear)
- Quantity = 10 (*Box of 5)

HCPCS A4409

(413157 and 401905 are HCPCS A4410)

Flange Size	Cutting Range	Acrylic Tape Collar WHITE	Acrylic Tape Collar TAN	Without Tape Collar
● 1 1/4" (32mm)	1/2 - 3/4"	413159	413164	413153
● 1 1/2" (38mm)	1/2 - 1"	413160	413165	413154
● 1 3/4" (45mm)	1/2 - 1 1/4"	413161	413166	413155
● 2 1/4" (57mm)	1/2 - 1 3/4"	413162	413167	413156
● 2 3/4" (70mm)	1/2 - 2 1/4"	413163	413168	413157
● 4" (100mm)	1/2 - 3 1/2"	N/A	N/A	401905*



NATURA™ SKIN BARRIER

Cut-to-Fit | Accordion Flange

- Durahesive® (Extended Wear)
- Quantity = 10

HCPCS A4409

Flange Size	Cutting Range	Acrylic Tape Collar WHITE
● 1 3/4" (45mm)	1/2 - 13/16" (13mm-21mm)	421454
● 2 1/4" (57mm)	1/2 - 1 5/16" (13mm-33mm)	421458
● 2 3/4" (70mm)	1/2 - 1 3/4" (13mm-45mm)	421462



NATURA™ SKIN BARRIER

Moldable | Accordion Flange

- Durahesive® (Extended Wear)
- Quantity = 10

HCPCS A4409

Flange Size	Stoma Size	Hydrocolloid Tape Collar TAN
● 2 1/4" (57mm)	1/2 - 7/8" (13-22mm)	421039
● 2 1/4" (57mm)	7/8 - 1 1/4" (22-33mm)	421040
● 2 3/4" (70mm)	1 1/4 - 1 3/4" (33-45mm)	421041

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

TWO-PIECE SYSTEMS

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- Natura™ + Pouches
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

NATURA™ SKIN BARRIER

**Cut-to-Fit | Convex |
Accordion Flange**

- Durahesive® (Extended Wear)
- Convex Wafer Depth = 6.65mm
- Quantity = 10

HCPCS A4407

Flange Size	Cutting Range	Acrylic Tape Collar TAN	Hydrocolloid Tape Collar TAN
● 1 3/4" (45mm)	1/2 - 7/8" (13-22mm)	421638	421639
● 2 1/4" (57mm)	1/2 - 1 3/8" (13-35mm)	421640	421641
● 2 3/4" (70mm)	1/2 - 1 7/8" (13-48mm)	421642	421643

NATURA™ SKIN BARRIER

Moldable

- Durahesive® (Extended Wear)
- Quantity = 10

HCPCS A4409

Flange Size	Stoma Size	Acrylic Tape Collar WHITE	Hydrocolloid Tape Collar WHITE
● 1 3/4" (45mm)	1/2 - 7/8" (13-22mm)	413417	411800
● 1 3/4" (45mm)	7/8 - 1 1/4" (22-33mm)	413418	411802
● 2 1/4" (57mm)	1 1/4 - 1 3/4" (33-45mm)	413419	411804
● 2 3/4" (70mm)	1 1/4 - 1 3/4" by 1 3/4 - 2 1/8" (33x45mm by 45x56mm) Oval shape	413420	411806



NATURA™ SKIN BARRIER

Moldable | Convex

- Durahesive® (Extended Wear)
- Quantity = 10

HCPCS A4407

Flange Size	Stoma Size	Acrylic Tape Collar WHITE	Convex Wafer Depth
● 1 3/4" (45mm)	1/2 - 7/8" (13-22mm)	404592	6.35mm
● 1 3/4" (45mm)	7/8 - 1 1/4" (22-33mm)	404593	7mm
● 2 1/4" (57mm)	1 1/4 - 1 3/4" (33-45mm)	404594	7mm

SUR-FIT® NATURA™ SKIN BARRIER

Pre-cut | Convex

- Durahesive® (Extended Wear)
- Convex Wafer Depth = 6.65mm
- Quantity = 10

HCPCS A4407

Flange Size	Stoma Size	Acrylic Tape Collar WHITE
● 1 3/4" (45mm)	1/2" (13 mm)	413177
● 1 3/4" (45mm)	5/8" (16 mm)	413178
● 1 3/4" (45mm)	3/4" (19 mm)	413179
● 1 3/4" (45mm)	7/8" (22 mm)	413180
● 1 3/4" (45mm)	1" (25 mm)	413181
● 1 3/4" (45mm)	1 1/8" (28 mm)	413182
● 1 3/4" (45mm)	1 1/4" (32 mm)	413183
● 1 3/4" (45mm)	1 3/8" (35 mm)	413184
● 2 1/4" (57mm)	1 1/2" (38 mm)	413185
● 2 1/4" (57mm)	1 5/8" (41 mm)	413186
● 2 1/4" (57mm)	1 3/4" (45 mm)	413187
● 2 1/4" (57mm)	2" (50 mm)	413188



TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE
SYSTEMS

TWO-PIECE SYSTEMS

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- Natura™ + Pouches
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

ONE-PIECE
SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

NATURA™ STOMAHESIVE® SKIN BARRIERS

Compatible with Natura™ + and SUR-FIT® Natura™ pouches



SUR-FIT® NATURA™ SKIN BARRIER


Cut-to-Fit

- Stomahesive® (Standard Wear)
- Quantity = 10 (*Box of 5)

HPCPS A4414

(401577 and 401906 are HPCPS A4415)

	Flange Size	Cutting Range	Acrylic Tape Collar WHITE	Acrylic Tape Collar TAN	Without Tape Collar
●	1 1/4" (32mm)	1/2 - 3/4"	125257	125262	401573
●	1 1/2" (38mm)	1/2 - 1"	125258	125263	401574
●	1 3/4" (45mm)	1/2 - 1 1/4"	125259	125264	401575
●	2 1/4" (57mm)	1/2 - 1 3/4"	125260	125265	401576
●	2 3/4" (70mm)	1/2 - 2 1/4"	125261	125266	401577
●	4" (100mm)	1/2 - 3 1/2"	N/A	N/A	401906*

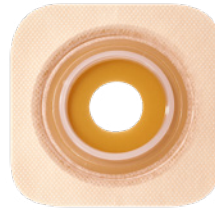

SUR-FIT® NATURA™ SKIN BARRIER

Pre-cut

- Stomahesive® (Standard Wear)
- Quantity = 10

HPCPS A4414

	Flange Size	Stoma Size	Acrylic Tape Collar TAN
●	1 3/4" (45mm)	1/2" (13 mm)	125267
●	1 3/4" (45mm)	5/8" (16 mm)	125268
●	1 3/4" (45mm)	3/4" (19 mm)	125269
●	1 3/4" (45mm)	7/8" (22 mm)	125270
●	1 3/4" (45mm)	1" (25 mm)	125271
●	1 3/4" (45mm)	1 1/8" (28 mm)	125272
●	1 3/4" (45mm)	1 1/4" (32 mm)	125273
●	1 3/4" (45mm)	1 3/8" (35 mm)	125274
●	2 1/4" (57mm)	1 1/2" (38 mm)	125275
●	2 1/4" (57mm)	1 5/8" (41 mm)	125276
●	2 1/4" (57mm)	1 3/4" (45 mm)	125277



NATURA™ SKIN BARRIER

Cut-to-Fit | Accordion Flange

- Stomahesive® (Standard Wear)
- Quantity = 10

HPCPS A4414

	Flange Size	Cutting Range	Acrylic Tape Collar WHITE
●	1 3/4" (45mm)	1/2 - 13/16" (13-21mm)	421456
●	2 1/4" (57mm)	1/2 - 1 5/16" (13-33mm)	421460
●	2 3/4" (70mm)	1/2 - 1 3/4" (13-45mm)	421464

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

TWO-PIECE SYSTEMS

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- Natura™ + Pouches
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)



NATURA™ SKIN BARRIER

Moldable | Accordion Flange

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Flange Size	Stoma Size	Hydrocolloid Tape Collar TAN
● 2 1/4" (57mm)	1/2 - 7/8" (13-22mm)	421033
● 2 1/4" (57mm)	7/8 - 1 1/4" (22-33mm)	421034
● 2 3/4" (70mm)	1 1/4 - 1 3/4" (33-45mm)	421035



NATURA™ SKIN BARRIER

Moldable

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Flange Size	Stoma Size	Acrylic Tape Collar WHITE	Hydrocolloid Tape Collar WHITE
● 1 3/4" (45mm)	1/2 - 7/8" (13-22mm)	413421	411801
● 1 3/4" (45mm)	7/8 - 1 1/4" (22-33mm)	413422	411803
● 2 1/4" (57mm)	1 1/4 - 1 3/4" (33-45mm)	413423	411805
● 2 3/4" (70mm)	1 1/4 - 1 3/4" by 1 3/4 - 2 1/8" (33x45mm by 45x56mm) Oval shape	413424	411807

**WATCH NOW**

to learn how to apply
a Natura™ Moldable
Skin Barrier

NATURA™ + POUCHES

Compatible with SUR-FIT® Natura™ and Natura™ skin barriers



NATURA™+ POUCH

**Drainable | InvisiClose™
Tail Closure**

- Pouch Length = 12" (*pouch length = 14")
- Quantity = 10

HCPCS A5063 | No Filter

	Flange Size	Transparent	Opaque
●	1 3/4" (45mm)	416418	421738
●	2 1/4" (57mm)	416421	421739
●	2 3/4" (70mm)	416424	421740
●	4" (100mm)	416472*	N/A

HCPCS A4425 | Filter (Filter covers included)

	Flange Size	Transparent	Opaque w/ Easy-View Window	Opaque
●	1 1/2" (38mm)	416414	421746	416415
●	1 3/4" (45mm)	416416	421747	416417
●	2 1/4" (57mm)	416419	421748	416420
●	2 3/4" (70mm)	416422	421749	416423

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

**TWO-PIECE
SYSTEMS****TWO-PIECE SYSTEMS**

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- **Natura™ + Pouches**
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

**ONE-PIECE
SYSTEMS**

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

SUR-FIT® NATURA™ POUCH

Drainable | Clipped Closure

- Pouch Length = 10"
- 2 Tail Closures per box
- Quantity = See table below



HPCPS A5063 | No Filter

	Flange Size	Opaque Qty = 10	Opaque Qty = 20
●	1 1/4" (32mm)	401505	404030
●	1 1/2" (38mm)	401506	404031
●	1 3/4" (45mm)	401507	404032
●	2 1/4" (57mm)	401508	404033
●	2 3/4" (70mm)	401509	404034

SUR-FIT® NATURA™ POUCH

Drainable | Clipped Closure

- Pouch Length = 12"
(*pouch length = 14" ; **pouch length = 10")
- 2 Tail Closures per box
- Quantity = See table below



HPCPS A5063 | No Filter

	Flange Size	Transparent Qty = 10	Opaque Qty = 10	Opaque Qty = 20	Opaque w/ 1 comfort panel Qty = 10	Opaque w/ 1 comfort panel Qty = 20
●	1 1/4" (32mm)	401510	401500	404025	401932	404014
●	1 1/2" (38mm)	401511	401501	404026	401933	404015
●	1 3/4" (45mm)	401512	401502	404027	401934	404016
●	2 1/4" (57mm)	401513	401503	404028	401935	404017
●	2 3/4" (70mm)	401514	401504	404029	401936	404018
●	4" (100mm)	401515*	N/A	N/A	N/A	N/A

Continued ➡

HPCPS A4425 | Filter (Filter covers included)

	Flange Size	Transparent Qty = 10	Opaque Qty = 20	Opaque Qty = 20**
●	1 1/4" (32mm)	N/A	411489	411288**
●	1 1/2" (38mm)	411264	411490	411289**
●	1 3/4" (45mm)	411265	411491	411290**
●	2 1/4" (57mm)	411266	411492	411291**
●	2 3/4" (70mm)	411267	411493	411292**

SUR-FIT® NATURA™ POUCH

Drainable | High-Output

- Anti-Reflux Valve
- Outlet with spout and cap for high output drainage
- Pouch Length = 14"
- 1 Tail Closure, 1 Night Drainage Adapter, and 5 Replacement Filters per box
- Quantity = 5



HPCPS A4413 | Filter

	Flange Size	Transparent	Opaque
●	1 3/4" (45mm)	420695	401557
●	2 1/4" (57mm)	420696	401558
●	2 3/4" (70mm)	420697	401559

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

TWO-PIECE SYSTEMS

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- Natura™ + Pouches
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample



NATURA™+ POUCH

Closed-End

- Pouch Length = 8"
- Quantity = 30

HCPCS A5054 | No Filter

	Flange Size	Opaque w/ Easy-View Window
●	1 1/4" (32mm)	421677
●	1 1/2" (38mm)	421679
●	1 3/4" (45mm)	421681
●	2 1/4" (57mm)	421683
●	2 3/4" (70mm)	421893

HCPCS A4419 | Filter (Filter covers included)

	Flange Size	Opaque w/ Easy-View Window	Opaque
●	1 1/4" (32mm)	421676	N/A
●	1 1/2" (38mm)	421678	416403
●	1 3/4" (45mm)	421680	416406
●	2 1/4" (57mm)	421682	416409
●	2 3/4" (70mm)	421798	416412

NATURA™+ POUCH

Closed-End

- Pouch Length = 6"
- Quantity = 30

HCPCS A4419 | Filter (Filter covers included)

	Flange Size	Opaque w/ Easy-View Window
●	1 3/4" (45mm)	421796
●	2 1/4" (57mm)	421797

SUR-FIT® NATURA™ POUCH

Closed-End | Flange Cap

- Quantity = 25

HCPCS A5055 | Filter (Filter covers included)

	Flange Size	Opaque w/ 1 comfort panel
●	1 3/4" (45mm)	401909



TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE
SYSTEMS

TWO-PIECE SYSTEMS

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- **Natura™ + Pouches**
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

ONE-PIECE
SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample



SUR-FIT® NATURA™ POUCH

Urostomy | Accuseal™ Tap with Valve

- Anti-Reflux Valve
- Pouch Length = 10" (*pouch length = 9")
- 2 Accuseal Night Drainage Adapters per box
- Quantity = 10

HCPCS A4432

	Flange Size	Transparent	Opaque	Opaque*
●	1 1/4" (32mm)	401542	401551	401547*
●	1 1/2" (38mm)	401543	401552	401548*
●	1 3/4" (45mm)	401544	401553	401549*
●	2 1/4" (57mm)	401545	401554	401550*
●	2 3/4" (70mm)	401546	401555	N/A

NATURA™+ POUCH

Urostomy | Soft Tap

- Anti-Reflux Valve
- Pouch Length = 10" (*pouch length = 9")
- 2 Soft Tap Drainage Adapters per box
- Quantity = 10

HCPCS A5073

	Flange Size	Transparent	Transparent*
●	1 1/4" (32mm)	N/A	421074*
●	1 3/4" (45mm)	413437	N/A
●	2 1/4" (57mm)	413438	N/A
●	2 3/4" (70mm)	413439	N/A

SUR-FIT® NATURA™ POUCH

Urostomy | Fold Over Tab

- Anti-Reflux Valve
- Pouch Length = 10" (*pouch length = 9")
- 2 Universal Night Drainage Adapters per box
- Quantity = 10

HCPCS A5073

	Flange Size	Transparent	Transparent*
●	1 1/4" (32mm)	401533	401538*
●	1 1/2" (38mm)	401534	401539*
●	1 3/4" (45mm)	401535	401540*
●	2 1/4" (57mm)	401536	401541*
●	2 3/4" (70mm)	401537	N/A



TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

TWO-PIECE SYSTEMS

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- **Natura™ + Pouches**
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES


RESOURCES

[Request Product Sample](#)

ESTEEM SYNERGY® SKIN BARRIERS

Compatible with Esteem Synergy®+ Pouches

ESTEEM SYNERGY® SKIN BARRIER

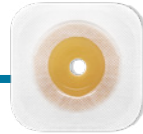

Cut-to-Fit

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Coupling Size	Cutting Range	Acrylic Tape Collar WHITE
Small Stoma Opening	up to 1 3/8" (35mm)	405456
Medium Stoma Opening	up to 1 7/8" (48mm)	405457
Large Stoma Opening	up to 2 3/8" (61mm)	405458

HCPCS A4415 (*HCPCS is A4414)

Coupling Size	Cutting Range	Without Tape Collar
Small Stoma Opening	up to 1 3/8" (35mm)	405467*
Medium Stoma Opening	up to 1 7/8" (48mm)	405468
Large Stoma Opening	up to 2 3/8" (61mm)	403948

ESTEEM SYNERGY® SKIN BARRIER

Moldable

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Coupling Size	Stoma Size	Hydrocolloid Tape Collar WHITE
Medium Stoma Opening	1/2 - 7/8" (13-22mm)	411657
Medium Stoma Opening	7/8 - 1 1/4" (22-33mm)	411658
Large Stoma Opening	1 1/4 - 1 3/4" (33-45mm)	411659

ESTEEM SYNERGY® SKIN BARRIER

Moldable | Convex

- Durahesive® (Extended Wear)
- Convex Wafer Depth = 6.65mm
- Quantity = 10



HCPCS A4407

Coupling Size	Stoma Size	Acrylic Tape Collar WHITE
Small Stoma Opening	1/2 - 7/8" (13-22mm)	409268
Medium Stoma Opening	7/8 - 1 1/4" (22-33mm)	409269
Large Stoma Opening	1 1/4 - 1 3/4" (33-45mm)	409270

ESTEEM SYNERGY® SKIN BARRIER

Pre-Cut

- Stomahesive® (Standard Wear)
- Quantity = 10

HCPCS A4414

Coupling Size	Stoma Size	Acrylic Tape Collar WHITE
Small Stoma Opening	5/8" (16 mm)	405473
Small Stoma Opening	3/4" (19 mm)	405474
Small Stoma Opening	7/8" (22 mm)	405475
Small Stoma Opening	1" (25 mm)	405476
Small Stoma Opening	1 1/8" (28 mm)	405477
Small Stoma Opening	1 1/4" (32 mm)	405478
Small Stoma Opening	1 3/8" (35 mm)	405479
Medium Stoma Opening	1 1/2" (38 mm)	405480
Medium Stoma Opening	1 5/8" (41 mm)	405481
Medium Stoma Opening	1 3/4" (45 mm)	405482

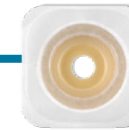


TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

TWO-PIECE SYSTEMS

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- Natura™ + Pouches
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

ESTEEM SYNERGY® POUCHES

Compatible with Esteem Synergy® Skin Barriers

ESTEEM SYNERGY® + POUCH

Drainable | InvisiClose™ Tail Closure

- Pouch Length = 12"
- Quantity = 10



HCPCS A5063 | No Filter

Coupling Size	Transparent
● Large Stoma Opening	416903

HCPCS A4425 | Filter (Filter covers included)

Coupling Size	Transparent	Opaque
● Small Stoma Opening	416795	416796
● Medium Stoma Opening	416798	416799
● Large Stoma Opening	416901	416902

ESTEEM SYNERGY® + POUCH

Closed-End

- Pouch Length = 8"
- Quantity = 30



HCPCS A4419 | Filter (Filter covers included)

Coupling Size	Transparent	Opaque
● Small Stoma Opening	416774	416775
● Medium Stoma Opening	416778	416779
● Large Stoma Opening	416782	416783

ESTEEM SYNERGY® + POUCH

Closed-End

- Pouch Length = 6"
- Quantity = 30



HCPCS A4419 | Filter (Filter covers included)

Coupling Size	Transparent	Opaque
● Small Stoma Opening	416776	416777
● Medium Stoma Opening	N/A	416781

ESTEEM SYNERGY® POUCH

Urostomy | Accuseal™ Tap with Valve

- Anti-Reflux Valve
- Pouch Length = 10" (*pouch length = 9")
- 2 Accuseal Night Drainage Adapters per box
- Quantity = 10



HCPCS A4432

Coupling Size	Transparent	Transparent*	Opaque
● Small Stoma Opening	405451	405446*	405448
● Medium Stoma Opening	405452	405447*	405449
● Large Stoma Opening	405453	N/A	405450

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

TWO-PIECE SYSTEMS

- Natura™ Durahesive® Skin Barriers
- Natura™ Stomahesive® Skin Barriers
- Natura™ + Pouches
- Esteem Synergy® Skin Barriers
- Esteem Synergy® Pouches

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

ONE-PIECE

Systems

TWO BRANDS
make up ConvaTec's
ONE-Piece Ostomy Systems



Esteem™ +
Pouches

Active Life®
Pouches

ESTEEM™

The Esteem™ + brand is an all-in-one solution combining the skin barrier and pouch in a single unit. Designed to be easy and convenient, it's one simple step to attach or remove. Available in multiple adhesive options, growing convexity options, and with ConvaTec's latest pouch technologies.

ActiveLife®

Active Life® is a legacy one-piece system brand that also combines the skin barrier and pouch in a single unit.

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

**ONE-PIECE
SYSTEMS**

ONE-PIECE SYSTEMS

- Esteem™ + Pouches
- ActiveLife® Pouches

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

ESTEEM™ + POUCHES



ESTEEM™+ POUCH

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Durahesive® Plus (*Extended Wear*)
- Acrylic Tape Collar
- Pouch Length = 12"
- Quantity = 10

HCPCS A5056 | Filter (Filter covers included)

Cutting Range

Transparent

Opaque

3/4 - 2 1/2" (19-64mm)

416976

416975



ESTEEM™+ POUCH

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 12"
- Quantity = 10

HCPCS A5061 | No Filter

Cutting Range

Transparent

Opaque w/ Easy-View Window

13/16 - 2 3/4" (20-70mm)

416721

421853

HCPCS A4424 | Filter (Filter covers included)

Cutting Range

Transparent

Opaque w/
Easy-View Window

Opaque

13/16 - 2 3/4" (20-70mm)

416718

421858

416719



ESTEEM™+ POUCH

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (*Standard Wear*)*
- Without Tape Collar
- Pouch Length = 14"
- Quantity = 10

HCPCS A5061 | No Filter

Cutting Range

Transparent

Up to 4" (100mm)

416908

*Double-layer adhesive (Stomahesive body-side, Durahesive top-side) which combines the skin friendly aspects of Stomahesive with the turtlenecking protection of Durahesive.



ESTEEM™+ POUCH

Cut-to-Fit | Drainable | Clipped Closure

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 12"
- 2 Tail Closures per box
- Quantity = 10

HCPCS A4424 | Filter (Filter covers included)*

Cutting Range

Transparent

13/16 - 2 3/4" (20-70mm)

421848

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE
SYSTEMS

ONE-PIECE SYSTEMS

- Esteem™+ Pouches
- ActiveLife® Pouches

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample



ESTEEM™+ POUCH

Pre-Cut | Drainable | InvisiClose™ Tail Closure

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 12"
- Quantity = 10

HPCPS A5061 | No Filter

Stoma Size	Opaque w/ Easy-View Window
1" (25mm)	421854
1 3/16" (30mm)	421855
1 3/8" (35mm)	421856
1 9/16" (40mm)	421857

HPCPS A4424 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque w/ Easy-View Window	Opaque
1" (25mm)	416724	421859	416725
1 3/16" (30mm)	416728	421860	416729
1 3/8" (35mm)	416732	421861	416733
1 9/16" (40mm)	416736	421862	416737



ESTEEM™+ POUCH

Pre-Cut | Drainable | Convex | InvisiClose™ Tail Closure

- Durahesive® (*Extended Wear*)
- Acrylic Tape Collar
- Convex Wafer Depth = 5mm
- Pouch Length = 12"
- Quantity = 10

HPCPS A5057 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque
3/4" (19mm)	416738	416739
7/8" (22mm)	416740	416741
1" (25mm)	416742	416743
1 1/8" (28mm)	416744	416745
1 1/4" (32mm)	416746	416747
1 3/8" (35mm)	416748	416749
1 1/2" (38mm)	416750	416751



ESTEEM™+ POUCH

Moldable | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 12"
- Quantity = 10

HPCPS A5061 | No Filter

Stoma Size	Transparent
13/16 - 1 3/16" (20-30mm)	413518
1 3/16 - 1 9/16" (30-40mm)	413519
1 9/16 - 1 15/16" (40-50mm)	413520

HPCPS A4424 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque w/ Easy-View Window
13/16 - 1 3/16" (20-30mm)	413515	413521
1 3/16 - 1 9/16" (30-40mm)	413516	413522
1 9/16 - 1 15/16" (40-50mm)	413517	413523



WATCH NOW

to learn how to apply
an Esteem™+ One-Piece
Moldable pouch

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ONE-PIECE SYSTEMS

- Esteem™+ Pouches
- ActiveLife® Pouches

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)



ESTEEM™+ FLEX CONVEX POUCH

Cut-to-Fit | Drainable | Convex

- Esteem™ + Flex Convex Adhesive (Standard Wear)
- Without Tape Collar
- Convex Wafer Depth = 7mm
- Safe Seal Clipless Closure (1 Tail Closure per box)
- Pouch Length = 8 2/3"
- Includes Belt Tabs
- Quantity = 10

HCPCS A4389 | Filter (Filter covers included)

Cutting Range	Plateau Size	Opaque w/ Easy-View Window
3/4 - 1" (20-25mm)	V3	421622
3/4 - 1 3/8" (20-35mm)	V2	421619
3/4 - 1 2/3" (20-43mm)	V1	421615
3/4 - 2 6/25" (20-57mm)	V0	422159


Pre-Cut | Drainable | Convex

HCPCS A4389 | Filter (Filter covers included)

Stoma Size	Plateau Size	Opaque w/ Easy-View Window
1" (25mm)	V2	421620
1 1/6" (30mm)	V2	421621
1" (25mm)	V1	421616
1 1/6" (30mm)	V1	421617
1 3/8" (35mm)	V1	421618



WATCH how Esteem™+ Soft Convex is changing lives.

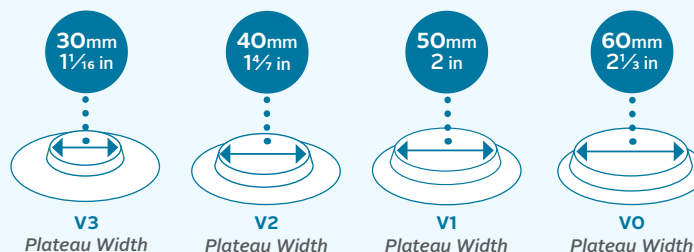Plateau Width Guidance for
Esteem™+ Flex and Soft Convex

ESTEEM™+ SOFT CONVEX POUCH

Cut-to-Fit | Drainable | Convex

- Esteem™ + Soft Convex Adhesive (Standard Wear)
- Without Tape Collar
- Convex Wafer Depth = 4mm
- Safe Seal Clipless Closure (1 Tail Closure per box)
- Pouch Length = 8 2/3"
- Includes Belt Tabs
- Quantity = 10

HCPCS A4389 | Filter (Filter covers included)

Cutting Range	Plateau Size	Opaque w/ Easy-View Window
3/8 - 1 1/8" (10-28mm)	V3	422365
3/5 - 1 3/5" (15-40mm)	V2	422364
3/4 - 1 7/8" (20-47mm)	V1	422363

Pre-Cut | Drainable | Convex

HCPCS A4389 | Filter (Filter covers included)

Stoma Size	Plateau Size	Opaque w/ Easy-View Window
1" (25mm)	V2	422545
1 1/6" (30mm)	V2	422546
1" (25mm)	V1	422539
1 1/6" (30mm)	V1	422540
1 3/8" (35mm)	V1	422541



ESTEEM™+ POUCH

Cut-to-Fit | Closed-End

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 8"
- Quantity = 30

HCPCS A4416 | Filter (Filter covers included)

Cutting Range	Transparent	Opaque w/ Easy-View Window	Opaque
13/16 - 2 3/4" (20-70mm)	416700	421821	416701


ESTEEM™+ POUCH

Cut-to-Fit | Closed-End

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 6"
- Quantity = 30

HCPCS A4416 | Filter (Filter covers included)

Cutting Range	Opaque w/ Easy-View Window
13/16 - 2 3/4" (20-70mm)	421816

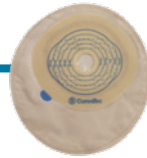

ESTEEM™+ POUCH

Pre-Cut | Closed-End

- Modified Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 8"
- Quantity = 30

HCPCS A4416 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque w/ Easy-View Window	Opaque
1" (25mm)	416703	421687	416704
1 3/16" (30mm)	416706	421688	416707
1 3/8" (35mm)	416709	421689	416710
1 9/16" (40mm)	416712	421822	416713
2" (50mm)	N/A	421823	416715



ESTEEM™+ POUCH

Moldable | Closed-End

- Stomahesive® (*Standard Wear*)
- Without Tape Collar
- Pouch Length = 8"
- Quantity = 30

HCPCS A4416 | Filter (Filter covers included)

Stoma Size	Opaque w/ Easy-View Window
13/16 - 1 3/16" (20-30mm)	413509
1 3/16 - 1 9/16" (30-40mm)	413510



TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ONE-PIECE SYSTEMS

- Esteem™+ Pouches
- ActiveLife® Pouches

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)



ESTEEM™+ POUCH

Urostomy | Cut-to-Fit |
Accuseal™ Tap with Valve

- Durahesive® (*Extended Wear*)
- Acrylic Tape Collar
- Anti-Reflux Valve
- Pouch Length = 10"
- Includes Belt Tabs
- 2 Accuseal Night Drainage Adapters per box
- Quantity = 10

HCPCS A4428

Cutting Range	Transparent
1/2 - 1 3/4" (13-45mm)	421631



ESTEEM™+ FLEX CONVEX POUCH

Urostomy | Cut-to-Fit | Convex

- Esteem™ + *Flex* Convex Adhesive (*Standard Wear*)
- Without Tape Collar
- Convex Wafer Depth = 7mm
- Soft Tap Closure
- Pouch Length = 7 2/3"
- Includes Belt Tabs
- 2 Esteem™ + *Flex* Convex Night Drainage Adapters per box
- Quantity = 10

HCPCS A4392

Cutting Range	Plateau Size	Opaque w/ Easy-View Window
3/8 - 1" (10-25mm)	V3	421628
3/8 - 1 3/8" (10-35mm)	V2	421626
3/8 - 1 2/3" (10-43mm)	V1	421623

Urostomy | Pre-Cut | Convex

HCPCS A4392

Stoma Size	Plateau Size	Opaque w/ Easy-View Window
9/16" (15mm)	V3	421629
9/16" (15mm)	V2	421627
9/16" (15mm)	V1	421624
4/5" (20mm)	V1	421625



ESTEEM™+ SOFT CONVEX POUCH

Urostomy | Cut-to-Fit | Convex

- Esteem™ + *Soft* Convex Adhesive (*Standard Wear*)
- Without Tape Collar
- Convex Wafer Depth = 4mm
- Soft Tap Closure
- Pouch Length = 7 2/3"
- Includes Belt Tabs
- 2 Esteem™ + *Soft* Convex Night Drainage Adapters per box
- Quantity = 10

HCPCS A4392

Cutting Range	Plateau Size	Opaque w/ Easy-View Window
3/8 - 1 1/8" (10-28mm)	V3	422562
3/5 - 1 3/5" (15-40mm)	V2	422557
3/4 - 1 7/8" (20-47mm)	V1	422551

Urostomy | Pre-Cut | Convex

HCPCS A4392

Stoma Size	Plateau Size	Opaque w/ Easy-View Window
4/5" (20mm)	V3	422564
1" (25mm)	V3	422565
4/5" (20mm)	V2	422558
1" (25mm)	V2	422559
1 1/6" (30mm)	V2	422560
1" (25mm)	V1	422552
1 1/6" (30mm)	V1	422553

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE
SYSTEMS

ONE-PIECE SYSTEMS

- Esteem™ + Pouches
- ActiveLife® Pouches

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

ACTIVELIFE® POUCHES

ACTIVELIFE® POUCH

Cut-to-Fit | Drainable | Clipped Closure

- Stomahesive® (Standard Wear)
- Acrylic Tape Collar
- Pouch Length = 12"
- 1 Tail Closure per box
- Quantity = See table below

HCPCS A5061 | No Filter

Cutting Range	Transparent Qty = 1	Transparent Qty = 10	Transparent Qty = 20
3/4 - 2 1/2" (19-64mm)	025433	022771	125330

ACTIVELIFE® POUCH

Cut-to-Fit | Drainable | Clipped Closure

- Durahesive® Plus (Extended Wear)
- Acrylic Tape Collar
- Pouch Length = 12"
- 1 Tail Closure per box
- Quantity = 10

HCPCS A4388 | No Filter

Cutting Range	Transparent	Opaque
3/4 - 2 1/2" (19-64mm)	400598	400599

ACTIVELIFE® POUCH

Pre-Cut | Drainable | Clipped Closure

- Stomahesive® (Standard Wear)
- Acrylic Tape Collar
- Pouch Length = 12" (*pouch length = 10")
- 1 Tail Closure per box
- Quantity = See table below

HCPCS A5061 | No Filter

Stoma Size	Transparent Qty = 10	Opaque Qty = 10	Opaque Qty = 20	Opaque* Qty = 10	Opaque* Qty = 20
3/4" (19mm)	022764	022757	125338	022750*	125331*
1" (25mm)	022765	022758	125339	022751*	N/A
1 1/4" (32mm)	022766	022759	125340	022752*	125333*
1 1/2" (38mm)	022767	022760	125341	022753*	125334*
1 3/4" (45mm)	022768	022761	125342	022754*	125335*
2" (50mm)	022769	022762	125343	022755*	125336*
2 1/2" (64mm)	022770	022763	125344	022756*	125337*

ACTIVELIFE® POUCH

Pre-Cut | Drainable | Convex | Clipped Closure

- Durahesive® (Extended Wear)
- Acrylic Tape Collar
- Pouch Length = 12"
- 1 Tail Closure per box
- Quantity = See table below

HCPCS A4390 | No Filter

Stoma Size	Transparent Qty = 5
3/4" (19mm)	175777
7/8" (22mm)	175778
1" (25mm)	175779
1 1/8" (28mm)	175780
1 1/4" (32mm)	175781
1 3/8" (35mm)	175782
1 1/2" (38mm)	175783
1 3/4" (45mm)	175784
2" (50mm)	175785

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ONE-PIECE SYSTEMS

- Esteem™ + Pouches
- ActiveLife® Pouches

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

ACTIVELIFE® POUCH

Cut-to-Fit | Stoma Cap

- Stomahesive® (Standard Wear)
- Without Tape Collar
- Quantity = 30

HCPCS A5055 | Filter (Filter covers included)

Cutting Range	Opaque
3/4 - 2" (19 - 50mm)	175611



ACTIVELIFE® POUCH

Urostomy | Pre-Cut | Accuseal™ Tap with Valve

- Durahesive® (Extended Wear)
- Acrylic Tape Collar
- Anti-Reflux Valve
- Pouch Length = 8"
- Includes Belt Tabs
- 2 Accuseal Night Drainage Adapters per box
- Quantity = 10

HCPCS A4428

Stoma Size	Transparent
3/4" (19mm)	650828
7/8" (22mm)	650829
1" (25mm)	650830
1 1/8" (28mm)	650831
1 1/4" (32mm)	650832
1 1/2" (38mm)	650833



ACTIVELIFE® POUCH

Urostomy | Cut-to-Fit | Fold Over Tap

- Stomahesive® (Standard Wear)
- Without Tape Collar
- Pouch Length = 9"
- Includes Belt Tabs
- 1 Universal Night Drainage Adapter per box
- Quantity = 10

HCPCS A5071

Cutting Range	Transparent
3/4 - 1 3/4" (19-45mm)	064927



ACTIVELIFE® POUCH

Urostomy | Pre-Cut | Convex | Fold Over Tap

- Durahesive® (Extended Wear)
- Acrylic Tape Collar
- Pouch Length = 8"
- Includes Belt Tabs
- 1 Universal Night Drainage Adapter per box
- Quantity = See table below

HCPCS A4393

Stoma Size	Transparent Qty = 5	Transparent Qty = 10
1/2" (13mm)	175790	125361
5/8" (16mm)	175791	125362
3/4" (19mm)	175792	125363
7/8" (22mm)	175793	125364
1" (25mm)	175794	125365
1 1/8" (28mm)	175795	125366
1 1/4" (32mm)	175796	125367
1 3/8" (35mm)	175797	125368
1 1/2" (38mm)	175798	125369



TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ONE-PIECE SYSTEMS

- Esteem™+ Pouches
- ActiveLife® Pouches

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

ACUTE KITS

ConvaTec's range of Acute Kits includes Post-Operative/Surgical kits in Sterile and Non-Sterile combinations, and Unit Dose kits.

Information on the components included in the Acute Kits can be found below the product information tables.

Color Key & Sizing Guide







COLOR MATCH						
FLANGE SIZE	1 1/4" (32mm)	1 1/2" (38mm)	1 3/4" (45mm)	2 1/4" (57mm)	2 3/4" (70mm)	4" (100mm)

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

- ACUTE KITS**
- Two-Piece Kits
 - One-Piece Kits
 - Unit-Dose Kits

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

POST-OPERATIVE/ SURGICAL KIT TWO-PIECE SYSTEM

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Double Layer Skin Barrier† (Standard Wear)
- Without Tape Collar
- Pouch Color: Transparent
- Pouch Length = 12" (*pouch length is 14")
- 1 Low-Pressure Adaptor
- Quantity = 5 Kits

HCPCS N/A

	Flange Size	Non-Sterile Kit	Sterile Kit	Ostomy Loop Rod
●	1 3/4" (45mm)	416919	416920	None
●	2 1/4" (57mm)	416923	416924	None
●	2 3/4" (70mm)	416927	416928	65mm
●	4" (100mm)	416929	416930	90mm

*Double-layer adhesive (Stomahesive body-side, Durahesive top-side) which combines the skin friendly aspects of Stomahesive with the turtlenecking protection of Durahesive.

Above kits include the following components:

Kit	Comparable to Skin Barrier in Kit* Stomahesive®	Comparable to Skin Barrier in Kit* Durahesive®	Pouch In Kit	Low Pressure Adaptor	Ostomy Loop Rod
416919/ 416920	401575	413155	416418	401993	None
416923/ 416924	401576	413156	416421	401994	None
416927/ 416928	401577	413157	416424	401995	022355
416929/ 416930	401906	401905	416472	401996	022356

*The double-layer skin barrier in the kit is not available in the retail market. Comparable skin barriers are in the chart above.

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Moldable | Drainable | InvisiClose™ Tail Closure

- Durahesive® (Extended Wear)
- Hydrocolloid Tape Collar (White)
- Pouch Color: Transparent
- Pouch Length = 12"
- 1 Low-Pressure Adaptor
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416917
●	2 1/4" (57mm)	416921
●	2 3/4" (70mm)	416925

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit	Low Pressure Adaptor
416917	411802	416418	401993
416921	411804	416421	401994
416925	411806	416424	401995

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

ACUTE KITS

- Two-Piece Kits
- One-Piece Kits
- Unit-Dose Kits

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Moldable | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (Standard Wear)
- Hydrocolloid Tape Collar (White)
- Pouch Color: Transparent
- Pouch Length = 12"
- 1 Low-Pressure Adaptor
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416918
●	2 1/4" (57mm)	416922
●	2 3/4" (70mm)	416926

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit	Low Pressure Adaptor
416918	411803	416418	401993
416922	411805	416421	401994
416926	411807	416424	401995

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Moldable | Accordion Flange | Drainable | InvisiClose™ Tail Closure

- Durahesive® (Extended Wear)
- Hydrocolloid Tape Collar (Tan)
- Pouch Color: Transparent
- Pouch Length = 12"
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	2 1/4" (57mm)	421045
●	2 1/4" (57mm)	421046
●	2 3/4" (70mm)	421047

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit
421045	421039	416421
421046	421040	416421
421047	421041	416424

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

- Two-Piece Kits
- One-Piece Kits
- Unit-Dose Kits

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

NATURA™ SKIN BARRIER & NATURA™ + POUCH KIT

Moldable | Accordion Flange | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (Standard Wear)
- Hydrocolloid Tape Collar (Tan)
- Pouch Color: Transparent
- Pouch Length = 12"
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	2 1/4" (57mm)	421042
●	2 1/4" (57mm)	421043
●	2 3/4" (70mm)	421044

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit
421042	421033	416421
421043	421034	416421
421044	421035	416424

NATURA™ SKIN BARRIER & SUR-FIT® NATURA™ + POUCH KIT

Moldable | Urostomy | Accuseal™ Tap with Valve

- Durahesive® (Extended Wear)
- Hydrocolloid Tape Collar (White)
- Anti-Reflux Valve
- Pouch Color: Transparent
- Pouch Length = 10"
- 1 Low Pressure Adaptor
- 1 Accuseal Night Drainage Adapter per box
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416932
●	2 1/4" (57mm)	416940
●	2 3/4" (70mm)	416948

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit	Low Pressure Adaptor
416932	411802	401544	401993
416940	411804	401545	401994
416948	411806	401546	401995

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

- Two-Piece Kits
- One-Piece Kits
- Unit-Dose Kits

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

NATURA™ SKIN BARRIER & SUR-FIT® NATURA™ + POUCH KIT

Moldable | Urostomy | Accuseal™ Tap with Valve

- Stomahesive® (Standard Wear)
- Hydrocolloid Tape Collar (White)
- Anti-Reflux Valve
- Pouch Color: Transparent
- Pouch Length = 10"
- 1 Low Pressure Adaptor
- 1 Accuseal Night Drainage Adapter per box
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416934
●	2 1/4" (57mm)	416942
●	2 3/4" (70mm)	416950

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit	Low Pressure Adaptor
416934	411803	401544	401993
416942	411805	401545	401994
416950	411807	401546	401995

NATURA™ SKIN BARRIER & SUR-FIT® NATURA™ + POUCH KIT

Cut-to-Fit | Urostomy | Accuseal™ Tap with Valve

- Double Layer Skin Barrier† (Standard Wear)
- Acrylic Tape Collar (White)
- Anti-Reflux Valve
- Pouch Color: Transparent
- Pouch Length = 10"
- 1 Low Pressure Adaptor
- 1 Accuseal Night Drainage Adapter per box
- Quantity = 5 Kits

HCPCS N/A | Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	416938
●	2 1/4" (57mm)	416946
●	2 3/4" (70mm)	416954

**Double-layer adhesive (Stomahesive body-side, Duraheasive top-side) which combines the skin friendly aspects of Stomahesive with the turtlenecking protection of Duraheasive.*

Above kits include the following components:

Kit	Comparable to Skin Barrier in Kit* Stomahesive®	Comparable to Skin Barrier in Kit* Duraheasive®	Pouch In Kit	Low Pressure Adaptor
416938	401575	413155	401544	401993
416946	401576	413156	401545	401994
416954	401577	413157	401546	401995

**The double-layer skin barrier in the kit is not available in the retail market. Comparable skin barriers are in the chart above.*

CONVATEC LOOP OSTOMY ROD

Sterile

- Quantity = 10

HCPCS A4421

	Flange Size	Size	Product Code
●	2 3/4" (70mm)	65mm	022355
●	4" (100mm)	90mm	022356



TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS**ACUTE KITS**

- Two-Piece Kits
- One-Piece Kits
- Unit-Dose Kits

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

POST-OPERATIVE/ SURGICAL KIT ONE-PIECE SYSTEM

ESTEEM™ + POUCH KIT

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Double Layer Skin Barrier* (Standard Wear)
- Without Tape Collar
- Color: Transparent
- Pouch Length = 14"
- Quantity = 5 Kits

HCPCS N/A

Cutting Range	Non-Sterile Kit	Sterile Kit
Can be enlarged to accommodate stoma sizes up to 4" (100mm)	416915	416916

*Double-layer adhesive (Stomahesive body-side, Durahesive top-side) which combines the skin friendly aspects of Stomahesive with the turtlenecking protection of Durahesive.

UNIT-DOSE KITS

SUR-FIT® NATURA™ SKIN BARRIER & SUR-FIT® NATURA™ + POUCH KIT

Cut-to-Fit | Urostomy | Accuseal™ Tap with Valve

- Durahesive® (Extended Wear)
- Acrylic Tape Collar (White)
- Pouch Color: Transparent
- Pouch Length = 10"
- 1 Accuseal Night Drainage Adapter per box
- Quantity = 5 Kits

HCPCS N/A | Non-Sterile

	Flange Size	Kit Product Code
●	1 3/4" (45mm)	401923
●	2 1/4" (57mm)	401924

Above kits include the following components:

Kit	Skin Barrier in Kit	Pouch In Kit
401923	413161	401544
401924	413162	401545

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

ACUTE KITS

- Two-Piece Kits
- One-Piece Kits
- Unit-Dose Kits

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)

PEDIATRIC





Products



Little Ones® is a pediatric Two-piece and One-piece system designed for babies smaller than 10 lbs, and children over 10 lbs.

The pouches have soft-welded edges and smaller pouch sizes designed to provide extra comfort for little ones.

Color Key & Sizing Guide

COLOR MATCH				
FLANGE/ COUPLING SIZE	0 - 9/10" (0-23mm)	1/5 - 1 1/4" (4.8-31mm)	1 1/4" (32mm)*	1 3/4" (45mm)*

*Indicates Flange size measurement

Color Key & Sizing Guide


 When ordering a Pediatric Two-Piece System, ***make sure you select the same flange/coupling size for the skin barrier and pouch.*** Use the colored dots in the table to help.

TABLE OF CONTENTS
me+™ SUPPORT
EDUCATION
TWO-PIECE SYSTEMS
ONE-PIECE SYSTEMS
ACUTE KITS
<div>PEDIATRIC PRODUCTS</div> <ul style="list-style-type: none"> • Extra-Small Systems • Two-Piece Systems • One-Piece Systems
ACCESSORIES
WOUND & FISTULA POUCHES
RESOURCES
Request Product Sample

LITTLE ONES® EXTRA SMALL TWO-PIECE OSTOMY SYSTEM

For babies smaller than 10lbs/4.5kg

LITTLE ONES® EXTRA SMALL SKIN BARRIER

Cut-to-Fit

- Stomahesive® (Standard Wear)
- Adhesive Coupling Technology™
- No Starter Hole
- Quantity = 5

HCPCS A4414

Cutting Range	Without Tape Collar
● 0 - 9/10" (0-23mm)	411632

LITTLE ONES® EXTRA SMALL POUCH

Drainable | InvisiClose™ Tail Closure

- Adhesive Coupling Technology™
- Pouch Length = 6"
- Quantity = 10

HCPCS A5063 | No Filter

Stoma Size	Transparent
● 0 - 9/10" (0-23mm)	411631

LITTLE ONES® EXTRA SMALL ONE-PIECE OSTOMY SYSTEM

For babies smaller than 10lbs/4.5kg

LITTLE ONES® EXTRA SMALL POUCH

Cut-to-Fit | Drainable | InvisiClose™ Tail Closure

- Stomahesive® (Standard Wear)
- Without Tape Collar
- No Starter Hole
- Pouch Length = 6"
- Quantity = 10

HCPCS A5061 | No Filter

Stoma Size	Transparent
0 - 9/10" (0-23mm)	411633

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC
PRODUCTS

PEDIATRIC PRODUCTS

- Extra-Small Systems
- Two-Piece Systems
- One-Piece Systems

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

LITTLE ONES® TWO-PIECE OSTOMY SYSTEM

Suggested for children over 10lbs/4.5kg


LITTLE ONES® SKIN BARRIER

Cut-to-Fit

- Stomahesive® (Standard Wear)
- Adhesive Coupling Technology™
- Quantity = 5

HCPCS A4414

Cutting Range	Hydrocolloid Tape Collar WHITE	No Tape Collar
● 1/5 - 1 1/4" (4.8-31mm)	411642	411641

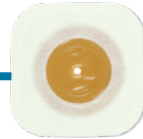

LITTLE ONES® POUCH

Drainable | InvisiClose™ Tail Closure

- Adhesive Coupling Technology™
- Pouch Length = 8"
- Quantity = 10

HCPCS A4425 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque
● 1/5 - 1 1/4" (4.8-31 mm)	411637	411638



LITTLE ONES® POUCH

Closed-End

- Adhesive Coupling Technology™
- Pouch Length = 6"
- Quantity = 10

HCPCS A4419 | Filter (Filter covers included)

Stoma Size	Transparent	Opaque
● 1/5 - 1 1/4" (4.8-31 mm)	411640	411639



LITTLE ONES® SKIN BARRIER

Cut-to-Fit

- Stomahesive® (Standard Wear)
- Flange System
- Quantity = 5

HCPCS A4414

Flange Size	Cutting Range	Acrylic Tape Collar WHITE
● 1 1/4" (32mm)	1/2 - 3/4" (13-19mm)	401925
● 1 3/4" (45mm)	1 1/8 - 1 1/4" (13-32mm)	401926



LITTLE ONES® POUCH

Drainable | Clipped Closure

- Stomahesive® (Standard Wear)
- Flange System
- Pouch Length = 6"
- 10 Soft, Moldable Tail Closures per box (Designed for Pediatric Use)
- Quantity = 10

HCPCS A5063 | No Filter

Flange Size	Transparent
● 1 1/4" (32mm)	401927
● 1 3/4" (45mm)	401928



TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC
PRODUCTS

PEDIATRIC PRODUCTS

- Extra-Small Systems
- Two-Piece Systems
- One-Piece Systems

ACCESSORIES

WOUND & FISTULA POUCHES


RESOURCES

Request Product Sample

LITTLE ONES® POUCH

Closed-End

- Flange System
- Pouch Length = 5"
- Quantity = 20

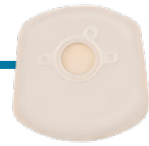**HPCPS A5054 | No Filter**

	Flange Size	Opaque
●	1 1/4" (32mm)	401930
●	1 3/4" (45mm)	401931

LITTLE ONES® POUCH

Urostomy | Fold Over Tap

- Flange System
- Pouch Length = 6"
- 2 Universal Night Drainage Adapters per box
- Quantity = 10

**HPCPS A5073 | No Filter**

	Flange Size	Transparent
●	1 1/4" (32mm)	401929

LITTLE ONES® ONE-PIECE OSTOMY SYSTEM

Suggested for children over 10lbs/4.5kg

LITTLE ONES® POUCH

Cut-to-Fit | Drainable | Clipped Closure

- Stomahesive® (Standard Wear)
- Without Tape Collar
- Pouch Length = 6"
- 15 Soft, Moldable Tail Closures per box (Designed for Pediatric Use)
- Quantity = 15

**HPCPS A5061 | No Filter**

Cutting Range	Transparent
5/16 - 2" (8-50mm)	020922

LITTLE ONES® POUCH

**Urostomy | Cut-to-Fit |
Fold Over Tap**

- Stomahesive® (Standard Wear)
- Without Tape Collar
- Pouch Length = 5"
- 2 Universal Night Drainage Adapters per box
- Quantity = 15

**HPCPS A5071**

Cutting Range	Transparent
5/16 - 1" (8-25mm)	020917

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

**PEDIATRIC
PRODUCTS****PEDIATRIC PRODUCTS**

- Extra-Small Systems
- Two-Piece Systems
- One-Piece Systems

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

ACCESSORIES

CHOOSING THE RIGHT ACCESSORIES

Whether you are looking for better skin protection or added security to your ostomy system, with the right combination of accessories, you can feel even more comfortable and at ease.




Skin Preparation



Pouching Support



Added Security

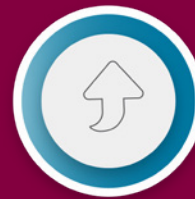

Pouch Removal

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

ACCESSORIES

- Skin Preparation
- Pouching Support
- Added Security
- Pouch Removal
- Irrigation

WOUND
& FISTULA
POUCHES

RESOURCES

[Request Product Sample](#)



SKIN PREPARATION

Solutions that help promote healthy skin by reducing irritation and damage.

ESENTA™



Sting-Free Skin Barrier Wipe

- Creates a breathable barrier on skin
- Protects skin for up to 72 hours¹⁷
- Reduces the risk of skin damage and irritation caused by body wastes, enzymes, and adhesives
- Safe and effective even with frequent use
- Silicone based
- Alcohol free

HCPCS A5120

Liquid per Wipe	Product Code	Box of	Boxes per Case*
3 mL	423392	25	20

*Sold by the case

ESENTA™



Sting-Free Skin Barrier Pump Spray

- Creates a breathable barrier on skin
- Protects skin for up to 72 hours¹⁷
- Reduces the risk of skin damage and irritation caused by body wastes, enzymes, and adhesives
- Safe and effective even with frequent use
- Silicone based
- Alcohol free

HCPCS A4369

Spray Bottle Size	Product Code	Box of	Boxes per Case*
28 mL	423286	1	12
50 mL	423288	1	12

*Sold by the case

ALLKARE®



Protective Barrier Wipe

- Designed to create a protective layer that shields skin from damage
- Helps reduce skin stripping and irritation from body waste
- May help increase the contact or adhesion to skin
- Contains Alcohol

HCPCS A5120

Product Code	Box of
037439	50
037444	100

STOMAHESIVE®



Protective Powder

- Helps protect skin by absorbing excess moisture
- For use on moist skin around a stoma
- Lightly dust affected area with powder, wiping off excess. May be followed up with a skin barrier wipe or spray

HCPCS A4371

Size	Product Code	Box of
1 oz.	025510	1 bottle

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

ACCESSORIES

- Skin Preparation
- Pouching Support
- Added Security
- Pouch Removal
- Irrigation

WOUND
& FISTULA
POUCHES

RESOURCES

Request Product Sample



POUCHING SUPPORT

Products that help create the best fit around the stoma, prevent leakage and help level out the application surface area for better pouch performance.

EAKIN COHESIVE®

Seals

- Moldable ring that helps seal the area around hard-to-fit stomas from any leakage under the skin barrier
- Acts as a filler on uneven skin surfaces to help create a flatter area for the skin barrier to adhere to
- Designed to absorb moisture and prevent it coming in contact with, and causing harm to the skin—forms a gel that protects the skin



HCPCS A4385

Size	Product Code	Depth	Box of
Small 2"	839002	4.2 mm	20
Slims 2"	839005	3 mm	10
Large 4"	839001	2.1 mm	10
StomaWrap 3 1/3"	839006	3 mm	10



WATCH NOW

to learn more about
eakin Cohesive® Seals

STOMAHESIVE®

Paste

- Helps to protect exposed skin between the base of the stoma and the opening of the skin barrier. Can also be used to fill crevices, skin folds and scars to help create a flat pouching surface for optimal appliance attachment
- Helps to prevent leakage under the skin barrier that may result in skin breakdown
- Pectin-based
- Contains alcohol



HCPCS A4406

Size	Product Code	Box of
2 oz.	183910	1 tube

EAKIN COHESIVE®

Paste

- Use to fill crevices, skin folds and scars to help create an even pouching surface
- Leakage prevention by helping form a waterproof barrier around the stoma
- Alcohol free



HCPCS A4405

Size	Product Code	Box of
2.1 oz.	839010	1 tube

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

ACCESSORIES

- Skin Preparation
- **Pouching Support**
- Added Security
- Pouch Removal
- Irrigation

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

STOMAHESIVE®

Skin Barrier

- Helps fill uneven skin surfaces to create a flat pouching surface for optimal appliance attachment
- Designed to help protect the skin from contact with bodily fluids
- Large cutting surface to accommodate large and swollen stomas
- Sterile Stomahesive® skin barrier can be used post-operatively

HCPCS A4362

Size	Product Code	Box of
4 x 4"	021712	5
4 x 4 (Sterile)	021901	5

HCPCS A5122

Size	Product Code	Box of
8 x 8"	021715	3

EAKIN COHESIVE®

Skin Barrier

- Designed to help protect the skin from contact with body fluids
- Acts as a filler on uneven skin surfaces to help create a flat pouching surface for optimal appliance attachment
- Moldable and moisture absorbing

HCPCS A4385

Size	Product Code	Box of
4 x 4"	839004	5
4 x 8"	839003	5

STOMAHESIVE®

Strips

- Helps to create a secure seal between the stoma and the appliance, to protect the skin against stoma output
- Can be cut, torn and shaped into a variety of shapes to help prevent leakage around the stoma by acting as a filler
- Designed to absorb moisture from the stoma, helping to create a snug fit around the stoma.

HCPCS A4362

Size	Product Code	Box of
3 cards of 5 strips each	025542	15 strips

SUR-FIT® NATURA™

Disposable Convex Inserts

- Convex inserts that can be added to a flat skin barrier to create a shallow convexity—insert into the inner diameter of the plastic flange on the pouch facing side of a SUR-FIT® Natura™ Durahesive® and Stomahesive® skin barrier
- Select an insert that is at least ¼" larger than the maximum diameter of the stoma
- For SUR-FIT® Natura™ skin barriers only. Not for use with ConvaTec Moldable Technology™ skin barriers, or any other product
- Quantity = 5

HCPCS A5093

	For Use With Flange Size	Stoma Size no Wider Than	Convex Insert Opening Size	Product Code
●	1 1/2" (38mm)	1/2"	3/4" (19mm)	404006
●	1 1/2" (38mm)	3/4"	1" (25mm)	404008
●	1 3/4" (45mm)	7/8"	1 1/8" (29mm)	404009
●	1 3/4" (45mm)	1"	1 1/4" (32mm)	404010
●	1 3/4" (45mm)	1 1/8"	1 3/8" (35mm)	404011
●	2 1/4" (57mm)	1 1/4"	1 1/2" (38mm)	404012
●	2 1/4" (57mm)	1 3/8"	1 5/8" (41mm)	404013

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

ACCESSORIES

- Skin Preparation
- Pouching Support
- Added Security
- Pouch Removal
- Irrigation

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample



ADDED SECURITY

An extra layer of protection to keep your ostomy system in place.

EASE™



Strips

- Helps keep your skin barrier in place
- Can be applied to cover the edges of the skin barrier
- Designed to be thin and flexible
- Water-resistant

HCPCS A4362 (4 strips = 1 billable unit of A4362)

Size	Product Code	Quantity
1 1/5 x 3 1/2" (3 x 9cm)	422163	80 strips

CONVATEC®

Ostomy Appliance Belt

- Used to help keep your ostomy system in place



HCPCS A4367

Size	Product Code	Box of
Adjustable up to 42"	175507	1

DIAMONDS™

Gelling and Odor Control Sachets

- Gels liquid contents
- Contains activated charcoal which helps to eliminate odors



HCPCS A4422

Product Code	Box of
TR105	100 Sachets

DUOLOCK™

Curved Tail Closure

- Helps secure the end of ConvaTec® Drainable pouch



HCPCS A4363

Product Code	Box of
175652	10

CONVATEC®

Tail Closure

- Helps secure the end of ConvaTec® Drainable pouch



HCPCS A4363

Product Code	Box of
175635	10

CONVATEC®

Night Drainage Container Set*

- Plastic vessel for use for overnight urine collection
- Container capacity of 2000 cc/ml
- Reusable, can be washed out between uses



HCPCS A5102

Product Code	Quantity
027060	1 set

*Contains: 1 ConvaTec® Night Drainage Container with tubing,
1 ConvaTec® Night Drainage Container Cover and 1 Universal Adapter

CONVATEC®

Night Drainage Container Tubing

- Replacement tubing for the ConvaTec® Night Drainage Container



HCPCS A4331

Product Code	Box of
027062	1

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

ACCESSORIES

- Skin Preparation
- Pouching Support
- Added Security
- Pouch Removal
- Irrigation

WOUND
& FISTULA
POUCHES

RESOURCES

Request Product Sample



POUCH REMOVAL

Solutions that help promote healthy skin by protecting it from irritation and damage caused by adhesives.

ESENTA™

Sting-Free Adhesive Remover Wipe

- For fast and easy adhesive release
- Reduces painful skin stripping and skin damage
- Quickly removes adhesive residue
- Feels silky and smooth, no residue left behind
- Silicone based
- Alcohol free



HCPCS A4456

Liquid per Wipe	Product Code	Box of	Boxes per Case*
3 mL	423391	25	20

*Sold by the case

ESENTA™

Sting-Free Adhesive Remover Spray

- For fast and easy adhesive release
- Reduces painful skin stripping and skin damage
- Quickly removes adhesive residue
- Feels silky and smooth, no residue left behind
- Silicone based
- Alcohol free



HCPCS A4455

Spray Bottle Size	Product Code	Box of	Boxes per Case*
50 mL	423289	1	12

*Sold by the case

ALLKARE®

Adhesive Remover Wipe

- Helps gently remove adhesives, including skin barriers, tapes and hydrocolloid dressings
- Contains alcohol



HCPCS A4456

Product Code	Box of
037436	50
037443	100

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

ACCESSORIES

- Skin Preparation
- Pouching Support
- Added Security
- **Pouch Removal**
- Irrigation

WOUND
& FISTULA
POUCHES

RESOURCES

Request Product Sample

IRRIGATION

VISI-FLOW®

Irrigator with Stoma Cone

HCPCS A4398 & A4399

Product Code	Box of
401989	1



VISI-FLOW®

Irrigation Sleeve

- Compatible with Natura™ and SUR-FIT® Natura™ skin barriers
- Color: Transparent
- 1 Tail Closure per box
- Quantity = 5



HCPCS A4436

	Flange Size	Product Code
●	1 1/2" (38mm)	401911
●	1 3/4" (45mm)	401912
●	2 1/4" (57mm)	401913
●	2 3/4" (70mm)	401914

VISI-FLOW®

Stoma Cone

HCPCS A4399

Product Code	Box of
022736	1



VISI-FLOW®

Irrigation Sleeve Tail Closure

HCPCS A4363

Product Code	Box of
175650	10



VISI-FLOW®

Irrigation Adaptor Faceplate

- Compatible with Visi-Flow® Irrigation Sleeves
- Quantity = 1



HCPCS A4361

	Flange Size	Product Code
●	1 3/4" (45mm)	401918
●	2 3/4" (70mm)	401919

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

ACCESSORIES

- Skin Preparation
- Pouching Support
- Added Security
- Pouch Removal
- Irrigation


WOUND
& FISTULA
POUCHES

RESOURCES

Request Product Sample

eakin® WOUND & FISTULA Pouches

eakin® wound pouches can help contain fistula and handle exuding wounds. Pouch options include the ability to connect tubing for remote drainage, or a fold & tuck option.

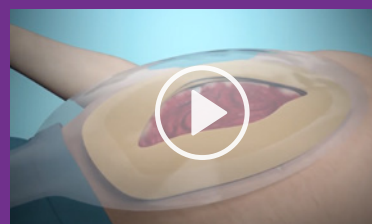

WATCH NOW

to learn how to
apply an eakin®
Wound Pouch

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

EAKIN®

Wound Pouch (1.2 x 1.8" Cutting Area)

- Quantity = 10

HCPCS A6154

Cutting Range	Fold and Tuck Closure
1.2 x 1.8" (3 x 4.5cm)	839250



EAKIN®

Wound Pouch (3.14 x 2.4" Cutting Area)

- Quantity = 10

HCPCS A6154

Cutting Range	Remote Drainage Attachment and Tap Closure
3.14 x 2.4" (8 x 6.1cm)	839267



EAKIN®

Wound Pouch (3.0 x 4.3" Cutting Area)

- Quantity = 10

HCPCS A6154

Cutting Range	Remote Drainage Attachment and Tap Closure	Fold and Tuck Closure
3.0 x 4.3" (7.6 x 10.9cm)	839261	839251



EAKIN®

Wound Pouch (4.3 x 6.9" Cutting Area)

- Quantity = 10

HCPCS A6154

Cutting Range	Remote Drainage Attachment and Tap Closure	Fold and Tuck Closure
4.3 x 6.9" (10.9 x 17.5cm)	839262	839252



EAKIN®

Horizontal Cutting Wound Pouch (9.7 x 6.3" Cutting Area)

- Quantity = 5

HCPCS A6154

Cutting Range	Remote Drainage Attachment and Tap Closure with Access Window	Fold and Tuck Closure
9.7 x 6.3" (24.6 x 16cm)	839264	839253



EAKIN®

Vertical Cutting Wound Pouch (6.3 x 9.7" Cutting Area)

- Quantity = 5

HCPCS A6154

Cutting Range	Remote Drainage Attachment and Tap Closure	Remote Drainage Attachment and Tap Closure with Access Window
6.3 x 9.7" (16 x 24.6cm)	839265	839268



EAKIN®

Vertical Cutting Wound Pouch (5.1 x 11.4" Cutting Area)

- Quantity = 5

HCPCS A6154

Cutting Range	Remote Drainage Attachment and Tap Closure
5.1 x 11.4" (13 x 29cm)	839266



EAKIN®

Eakin Access Window

- Add to Eakin wound and fistula pouches to allow instant access to the wound (*may be used with 839252, 839253, 839262, 839265, 839266*)
- Quantity = 5

HCPCS N/A

Size	Remote Drainage Attachment and Tap Closure
N/A	839280



TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

Request Product Sample

RESOURCES

Our resources will help you support your patients living with an ostomy.



- ✓ Discontinued Products
- ✓ Patient Ostomy Support & Resources
- ✓ Quick Links Guide
- ✓ Medicare Utilization Charts
- ✓ Product Code Index
- ✓ References

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND
& FISTULA
POUCHES

RESOURCES

RESOURCES

- Discontinued Products
- Patient Ostomy Support & Resources
- Quick Links Guide
- Medicare Utilization Chart
- Product Code Index
- References

[Request Product Sample](#)


DISCONTINUED PRODUCTS

SPECIFIC PRODUCT NUMBERS HAVE CHANGED[†]

We've made changes to our Product Portfolio

- Some legacy pouches have been upgraded to our newest technology.
- Sensi-Care™ Sting Free is now ESENTA™, available in skin barrier wipes and spray and adhesive remover wipes and spray.
- We have provided a cross reference (*see the following tables*) to map the discontinued products with its best replacement.





Use this guide to identify the product numbers being upgraded to + pouches or being discontinued.

A  sign indicates these products feature + pouch technology. Other modifications noted in the tables include: ConvaTec filter technology, package quantities, pouch lengths, color, adhesive, larger cutting areas and pre-cut holes.

Natura™ + Color Key & Sizing Guide

COLOR MATCH						
FLANGE SIZE	1 1/4" (32mm)	1 1/2" (38mm)	1 3/4" (45mm)	2 1/4" (57mm)	2 3/4" (70mm)	4" (100mm)

Esteem Synergy®+ Color Key & Sizing Guide

COLOR MATCH				
COUPLING SIZE	Small Stoma Opening	Medium Stoma Opening	Large Stoma Opening	Extra-Large Stoma Opening

[†] Only product numbers listed are being converted to + Pouches.
If an item is not listed, it is not changing.

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

RESOURCES

- Discontinued Products
- Patient Ostomy Support & Resources
- Quick Links Guide
- Medicare Utilization Chart
- Product Code Index
- References

[Request Product Sample](#)

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND
& FISTULA
POUCHES

RESOURCES

RESOURCES

• Discontinued Products

• Patient Ostomy
Support & Resources

• Quick Links Guide

• Medicare Utilization
Chart

• Product Code Index

• References

Request Product Sample

Previous SKU #	Previous QTY	New SKU #	New QTY	Flange or Stoma Size	Differences	New SKU Page #
125352	10	175777	5	19mm	Pack size change	66
125353	10	175778	5	22mm	Pack size change	66
125354	10	175779	5	25mm	Pack size change	66
125355	10	175780	5	28mm	Pack size change	66
125356	10	175781	5	32mm	Pack size change	66
125357	10	175782	5	35mm	Pack size change	66
125358	10	175783	5	35mm	Pack size change	66
125359	10	175784	5	45mm	Pack size change	66
175761	15	416700	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	61
175762	15	416703	30	25mm	⊕ Modified Stomahesive®	62
175763	15	416709	30	35mm	⊕ 3mm Larger Pre-cut Hole Modified Stomahesive®	62
175764	15	416712	30	40mm	⊕ 2mm Larger Pre-cut Hole Modified Stomahesive®	62
175765	15	421823	30	50mm	⊕ 5mm Larger Pre-cut Hole Opaque Modified Stomahesive®	62
175766	15	421823	30	50mm	⊕ Opaque Modified Stomahesive®	62
175768	15	421821	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	61
175769	15	421687	30	25mm	⊕ Modified Stomahesive®	62
175770	15	421689	30	35mm	⊕ 3mm Larger Pre-cut Hole Modified Stomahesive®	62
175771	15	421822	30	40mm	⊕ 2mm Larger Pre-cut Hole Modified Stomahesive®	62
175772	15	421823	30	50mm	⊕ 5mm Larger Pre-cut Hole Modified Stomahesive®	62
175773	15	421823	30	50mm	⊕ Modified Stomahesive®	62
401521	30	421679	30	38mm flange	⊕	45
401522	30	421681	30	45mm flange	⊕	45
401523	30	421683	30	57mm flange	⊕	45
401524	30	421893	30	70mm flange	⊕	45
401525	30	421678	30	38mm flange	⊕	45
401526	30	416406	30	45mm flange	⊕	45
401527	30	416409	30	57mm flange	⊕	45
401528	30	416412	30	70mm flange	⊕	45
401529	30	421677	30	32mm flange	⊕ Length = 8"	45
401530	20	421679	30	38mm flange	⊕ Length = 8"	45
401531	20	421796	30	45mm flange	⊕ Length = 6" Filter HCPCS now A4419	46
401532	20	421797	30	57mm flange	⊕ Length = 6" Filter HCPCS now A4419	46
401916	1	See Notes		45mm flange	Irrigation components available separately	
401917	1	See Notes		70mm flange	Irrigation components available separately	
401920	5	416919	5	45mm flange	Upgraded to our improved InvisiClose® closure system	71
401921	5	416927	5	70mm flange	Upgraded to our improved InvisiClose® closure system	71
401922	5	416929	5	100mm flange	Upgraded to our improved InvisiClose® closure system	71
401993	10	None		45mm flange	Refer to our line of low-pressure Accordion products	
401994	10	None		57mm flange	Refer to our line of low-pressure Accordion products	
401995	10	None		70mm flange	Refer to our line of low-pressure Accordion products	

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND
& FISTULA
POUCHES

RESOURCES

RESOURCES

• Discontinued Products

• Patient Ostomy
Support & Resources

• Quick Links Guide

• Medicare Utilization
Chart

• Product Code Index

• References

Request Product Sample

Previous SKU #	Previous QTY	New SKU #	New QTY		Flange or Stoma Size	Differences	New SKU Page #
402239	5	416923	5	●	57mm flange	Upgraded to our improved InvisiClose® closure system	71
404007	5	413180	10	●	38mm flange	Look at our extended convex line	38
405406	10	See Notes				Discontinued; see Natura®+ Pouch options	
409275	5	416929	5	●	100mm flange	Upgraded to our improved InvisiClose® closure system	71
411309	20	416415	10	●	38mm flange	⊕	42
411310	20	416417	10	●	45mm flange	⊕ Pack size change	42
411311	20	416420	10	●	57mm flange	⊕ Pack size change	42
411312	20	416423	10	●	70mm flange	⊕	42
411360	10	416416	10	●	45mm flange	⊕	42
411361	10	416419	10	●	57mm flange	⊕	42
411362	10	416422	10	●	70mm flange	⊕	42
413141	60	421821	30		20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	61
413142	60	421687	30		25mm	⊕ Modified Stomahesive®	62
413143	60	421689	30		35mm	⊕ 3mm Larger Pre-cut Hole Modified Stomahesive®	62
413144	60	421822	30		40mm	⊕ 2mm Larger Pre-cut Hole Modified Stomahesive®	62
413145	60	421823	30		50mm	⊕ 5mm Larger Pre-cut Hole Modified Stomahesive®	62
413146	60	421823	30		50mm	⊕ Modified Stomahesive®	62
413147	60	416700	30		20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	61
413148	60	416703	30		25mm	⊕ Modified Stomahesive®	62
413149	60	416709	30		35mm	⊕ 3mm Larger Pre-cut Hole Modified Stomahesive®	62
413150	60	416712	30		40mm	⊕ 2mm Larger Pre-cut Hole Modified Stomahesive®	62
413151	60	421823	30		50mm	⊕ 5mm Larger Pre-cut Hole Opaque Modified Stomahesive®	62
413152	60	421823	30		50mm	⊕ Opaque Modified Stomahesive®	62
413169	60	421679	30	●	38mm flange	⊕	45
413170	60	421681	30	●	45mm flange	⊕ Pack size change	45
413171	60	421683	30	●	57mm flange	⊕ Pack size change	45
413172	60	421893	30	●	70mm flange	⊕	45
413173	60	421678	30	●	38mm flange	⊕	45
413174	60	416406	30	●	45mm flange	⊕ Pack size change	45
413175	60	416409	30	●	57mm flange	⊕ Pack size change	42
413176	60	421798	30	●	70mm flange	⊕	42
413312	10	416418	10	●	45mm flange	⊕	42
413313	10	416421	10	●	57mm flange	⊕	42
413314	10	416424	10	●	70mm flange	⊕	42
413334	20	421738	10	●	45mm flange	⊕ Pack size change	42
413335	20	421739	10	●	57mm flange	⊕	42
413336	20	421740	10	●	70mm flange	⊕	42
413435	10	421074	10	●	32mm flange	Pouch length = 9"	47
413436	10	None		●	38mm flange	Discontinued	
413499	1	423289	1			Now branded ESENTA™	99

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND
& FISTULA
POUCHES

RESOURCES

RESOURCES

• Discontinued Products

• Patient Ostomy
Support & Resources

• Quick Links Guide

• Medicare Utilization
Chart

• Product Code Index

• References

Request Product Sample

Previous SKU #	Previous QTY	New SKU #	New QTY	Flange or Stoma Size	Differences	New SKU Page #
413500	30	423391	25		Now branded ESENTA™	99
413501	30	423392	25		Now branded ESENTA™	91
413502	1	423288	1		Now branded ESENTA™	91
413511	30	None		40-50mm	Discontinued	
416402	30	421677	30	● 32mm flange	Length = 8" Filter HCPCS now A5054	45
416404	30	421679	30	● 38mm flange	Enhanced with see-through window	45
416405	30	421679	30	● 38mm flange	Length = 8" Filter HCPCS now A5054	45
416407	30	421681	30	● 45mm flange	Enhanced with see-through window	45
416408	30	421796	30	● 45mm flange	Filter added - HCPCS now A5054	46
416410	30	421683	30	● 57mm flange	Enhanced with see-through window	45
416411	30	421797	30	● 57mm flange	Filter added - HCPCS now A5054	46
416413	30	421893	30	● 70mm flange	Enhanced with see-through window	45
416702	30	421816	30	20-70mm	⊕	61
416780	30	416781	30	● Medium	Opaque	52
416797	10	416795	10	● Small	Filter HCPCS now A4425	51
416900	10	416798	10	● Medium	Filter HCPCS now A4425	51
420793	5	Coming Soon			Discontinued	
420794	25	Coming Soon			Discontinued	
420795	5	Coming Soon			Discontinued	
420796	25	Coming Soon			Discontinued	
420797	1	423286	1		Now branded ESENTA™	91
420798	1	Coming Soon			Discontinued	
421075	10	None		● 38mm flange	Discontinued	
421076	10	413437	10	● 45mm flange	Pouch length = 10"	47
421612	30	None		20-25mm	Discontinued	
421613	30	None		20-35mm	Discontinued	
421614	30	None		20-43mm	Discontinued	
650408	30	421821	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	61
650422	30	416700	30	20-70mm	⊕ Larger Cutting Area Modified Stomahesive®	61
839260	10	839267	10	4.3 x 3"	New size	105
839263	5	839265	5	9.7 x 6.3"	New size	106

OSTOMY RESOURCES & SUPPORT

for You and Your Patient

QUICK LINKS GUIDE

Your guide for quick access to ConvaTec's library of resources in ostomy care, products and support.

ConvaTec

ConvaTec is a global medical products and technologies company developing ostomy products and solutions for almost 40 years.

Toll Free: 1-800-422-8811
www.convatec.com/ostomy

Ostomysecrets®

A ConvaTec company, offers a full line of stylish and functional undergarments, swim wear and accessories designed to support your ostomy solution.

Toll Free: 1-877-613-6246
www.ostomysecrets.com

OTHER OSTOMY RESOURCES

Wound Ostomy Continence Nurses Society™ (WOCN®)

Professional nursing society dedicated to providing expert health care to individuals with wound, ostomy and continence concerns.

Toll Free: 1-888-224-9626
www.wocn.org

American Cancer Society

Nationwide community-based, voluntary healthcare organization.

Toll Free: 1-800-ACS-2345
www.cancer.org

me+™

Support, insights & product information for those living with an ostomy. Enrolling in me+ is fast, easy and free.

Toll Free: 1-800-422-8811
e-mail: CIC@convatec.com
www.mepluscare.com

Ostomy Nurse Solutions App

A fast and easy way to help you support your patient without leaving the bedside.

www.ostomynursesolutions.com

The United Ostomy Associations of America (UOAA)

National organization providing support, information and advocacy.

Toll Free: 1-800-826-0826
www.ostomy.org

Crohn's & Colitis Foundation of America

Non-profit, volunteer-driven organization dedicated to finding a cure for Crohn's Disease and ulcerative colitis.

Toll Free: 1-800-932-2423
www.ccfa.org



Find Your
Category



Open
the Link



Explore the
Content

Ostomy Nurse Solutions App

Application Videos

ConvaTec Moldable Technology™
Product Information

me+™ Patient Enrollment

Accessories

Convexity Product Information

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND
& FISTULA
POUCHES

RESOURCES

RESOURCES

- Discontinued Products
- Patient Ostomy Support & Resources
- Quick Links Guide
- Medicare Utilization Chart
- Product Code Index
- References

Request Product Sample

HCPCS	HCPCS Descriptors	Medicare Usual Maximum Qty.
A4331	Extension Drainage Tube	1 per 6 months
A4361	Ostomy Faceplate	3 per 6 months
A4362	Skin Barrier, Solid, 4x4 or Equivalent	20 per month
A4363	Ostomy Clamp, Any Type, Replacement Only	As needed
A4367	Ostomy Belt	1 per month
A4369	Ostomy Skin Barrier Liquid (spray or brush)	2 oz per month
A4371	Ostomy Skin Barrier Powder	10 oz per 6 months
A4385	Skin Barrier, Solid 4x4 or Equivalent, Extended Wear	As needed
A4388	Pouch, Drainable, w/Extended Wear Barrier Attached	As needed
A4389	Ostomy Pouch, Drainable, w/ Barrier Attached, w/ Built-In Convexity, w/ Filter (1-piece)	20 per month
A4390	Pouch, Drainable, w/Extended Wear Barrier, Attached, w/convexity	As needed
A4392	Ostomy Pouch, Urinary, w/ Barrier Attached, w/ Built-In Convexity (1-piece)	20 per month
A4393	Ostomy Pouch, Urinary, w/ Extended Wear Barrier Attached, w/ Built-In Convexity (1-piece)	20 per month
A4397	Irrigation Sleeve	4 per month
A4398	Ostomy Irrigation Bag	2 per 6 months
A4399	Ostomy Irrigation Cone/Catheter, and Brush	2 per 6 months
A4402	Lubricant	4 oz per month
A4405	Ostomy Skin Barrier, Non-Pectin Based, Paste	4 oz per month
A4406	Ostomy Skin Barrier, Pectin-Based, Paste	4 oz per month
A4407	Skin Barrier, w/Flange, Extended Wear w/ Convexity	As needed
A4409	Skin Barrier, w/Flange, Extended Wear <= 4x4	As needed
A4410	Skin Barrier, w/Flange, Extended Wear > 4x4	As needed
A4413	Pouch, Drainable, High Output, for Barrier w/Flange	As needed
A4414	Skin Barrier, w/Flange, <= 4x4	20 per month
A4415	Skin Barrier, w/Flange, > 4x4	20 per month
A4416	Pouch, closed, w/Barrier Attached, w/Filter	60 per month
A4417	Ostomy Pouch, Closed, w/ Barrier Attached, w/ Built-In Convexity, w/ Filter (1-piece)	60 per month
A4419	Pouch, closed for use on Barrier w/Flange, w/Filter	60 per month
A4421	Ostomy Supply, Miscellaneous	As needed
A4422	Ostomy Absorbent Material (Sheet/Pad/Crystal Packet)	No Fee Schedule
A4424	Ostomy Pouch, Drainable, w/ Barrier Attached, w/ Filter (1-piece)	20 per month
A4425	Pouch, Drainable, with filter for use w/Flange	20 per month
A4428	Pouch, Urinary, w/ Extended Wear Barrier w/Faucet-type Tap	As needed
A4432	Pouch, Urinary, for use on Barrier w/Flange, w/Faucet-type Tap	20 per month
A4455	Adhesive Remover or Solvent	2 oz per month
A4456	Adhesive Remover Wipes, Any Type	50 per month
A5054	Pouch Closed, for Barrier w/Flange	60 per month
A5055	Stoma Cap, each	31 per month
A5056	Pouch, Drainable, w/Extended Wear Barrier Attached, w/filter	20 per month
A5057	Pouch, Drainable, w/Extended Wear Barrier, w/convexity, w/filter	20 per month
A5061	Pouch, Drainable, w/Barrier Attached	20 per month
A5063	Pouch, Drainable, for use w/Flange	20 per month
A5071	Pouch, Urinary, w/Barrier Attached	20 per month
A5073	Pouch, Urinary, for Barrier w/Flange	20 per month
A5093	Ostomy Accessory, Convex Insert	10 per month
A5102	Bedside Drainage Bottle w/ or wo/ tubing, rigid or expandable	2 per 6 months
A5120	Skin Barrier, Wipes, each	150 per 6 months
A5122	Skin Barrier, Solid, 8 X 8 or Equivalent	20 per month
A6154	Wound Pouch	As needed

125257.....	39	175635.....	98
125258.....	39	175650.....	102
125259.....	39	175652.....	98
125260.....	39	175777.....	66
125261.....	39	175778.....	66
125262.....	39	175779.....	66
125263.....	39	175780.....	66
125264.....	39	175781.....	66
125265.....	39	175782.....	66
125266.....	39	175783.....	66
125267.....	40	175784.....	66
125268.....	40	175785.....	66
125269.....	40	175790.....	68
125270.....	40	175791.....	68
125271.....	40	175792.....	68
125272.....	40	175793.....	68
125273.....	40	175794.....	68
125274.....	40	175795.....	68
125275.....	40	175796.....	68
125276.....	40	175797.....	68
125277.....	40	175798.....	68
125330.....	65	183910.....	94
125331.....	66	400598.....	65
125333.....	66	400599.....	65
125334.....	66	401500.....	43
125335.....	66	401501.....	43
125336.....	66	401502.....	43
125337.....	66	401503.....	43
125338.....	66	401504.....	43
125339.....	66	401505.....	43
125340.....	66	401506.....	43
125341.....	66	401507.....	43
125342.....	66	401508.....	43
125343.....	66	401509.....	43
125344.....	66	401510.....	43
125361.....	68	401511.....	43
125362.....	68	401512.....	43
125363.....	68	401513.....	43
125364.....	68	401514.....	43
125365.....	68	401515.....	43
125366.....	68	401533.....	48
125367.....	68	401534.....	48
125368.....	68	401535.....	48
125369.....	68	401536.....	48
175507.....	97	401537.....	48
175611.....	67	401538.....	48

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND
& FISTULA
POUCHES

RESOURCES

RESOURCES

- Discontinued Products
- Patient Ostomy Support & Resources
- Quick Links Guide
- Medicare Utilization Chart
- Product Code Index
- References

Request Product Sample

401539.....	48	401936.....	43	405480.....	50	413160.....	35
401540.....	48	401989.....	101	405481.....	50	413161.....	35
401541.....	48	404006.....	96	405482.....	50	413162.....	35
401542.....	47	404008.....	96	403948.....	49	413163.....	35
401543.....	47	404009.....	96	409268.....	50	413164.....	35
401544.....	47	404010.....	96	409269.....	50	413165.....	35
401545.....	47	404011.....	96	409270.....	50	413166.....	35
401546.....	47	404012.....	96	411264.....	44	413167.....	35
401547.....	47	404013.....	96	411265.....	44	413168.....	35
401548.....	47	404014.....	43	411266.....	44	413177.....	38
401549.....	47	404015.....	43	411267.....	44	413178.....	38
401550.....	47	404016.....	43	411288.....	44	413179.....	38
401551.....	47	404017.....	43	411289.....	44	413180.....	38
401552.....	47	404018.....	43	411290.....	44	413181.....	38
401553.....	47	404025.....	43	411291.....	44	413182.....	38
401554.....	47	404026.....	43	411292.....	44	413183.....	38
401555.....	47	404027.....	43	411489.....	44	413184.....	38
401557.....	44	404028.....	43	411490.....	44	413185.....	38
401558.....	44	404029.....	43	411491.....	44	413186.....	38
401559.....	44	404030.....	43	411492.....	44	413187.....	38
401573.....	39	404031.....	43	411493.....	44	413188.....	38
401574.....	39	404032.....	43	411631.....	83	413417.....	37
401575.....	39	404033.....	43	411632.....	83	413418.....	37
401576.....	39	404034.....	43	411633.....	84	413419.....	37
401577.....	39	404592.....	38	411637.....	85	413420.....	37
401905.....	35	404593.....	38	411638.....	85	413421.....	41
401906.....	39	404594.....	38	411639.....	85	413422.....	41
401909.....	46	405446.....	52	411640.....	85	413423.....	41
401911.....	101	405447.....	52	411641.....	85	413424.....	41
401912.....	101	405448.....	52	411642.....	85	413437.....	47
401913.....	101	405449.....	52	411657.....	49	413438.....	47
401914.....	101	405450.....	52	411658.....	49	413439.....	47
401918.....	102	405451.....	52	411659.....	49	413500.....	62
401919.....	102	405452.....	52	411800.....	37	413510.....	62
401923.....	80	405453.....	52	411801.....	41	413515.....	58
401924.....	80	405456.....	49	411802.....	37	413516.....	58
401925.....	86	405457.....	49	411803.....	41	413517.....	58
401926.....	86	405458.....	49	411804.....	37	413518.....	58
401927.....	86	405467.....	49	411805.....	41	413519.....	58
401928.....	86	405468.....	49	411806.....	37	413520.....	58
401929.....	87	405473.....	50	411807.....	41	413521.....	58
401930.....	87	405474.....	50	413153.....	35	413522.....	58
401931.....	87	405475.....	50	413154.....	35	413523.....	58
401932.....	43	405476.....	50	413155.....	35	416403.....	45
401933.....	43	405477.....	50	413156.....	35	416406.....	45
401934.....	43	405478.....	50	413157.....	35	416409.....	45
401935.....	43	405479.....	50	413159.....	35	416412.....	45

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND
& FISTULA
POUCHES

RESOURCES

RESOURCES

- Discontinued Products
- Patient Ostomy Support & Resources
- Quick Links Guide
- Medicare Utilization Chart
- Product Code Index
- References

Request Product Sample

416414.....	42	416751.....	57	421033.....	41	421682.....	45
416415.....	42	416774.....	51	421034.....	41	421683.....	45
416416.....	42	416775.....	51	421035.....	41	421687.....	62
416417.....	42	416776.....	52	421039.....	36	421688.....	62
416418.....	42	416777.....	52	421040.....	36	421689.....	62
416419.....	42	416778.....	51	421041.....	36	421738.....	42
416420.....	42	416779.....	51	421042.....	75	421739.....	42
416421.....	42	416781.....	52	421043.....	75	421740.....	42
416422.....	42	416782.....	51	421044.....	75	421746.....	42
416423.....	42	416783.....	51	421045.....	74	421747.....	42
416424.....	42	416795.....	51	421046.....	74	421748.....	42
416472.....	42	416796.....	51	421047.....	74	421749.....	42
416700.....	61	416798.....	51	421074.....	47	421796.....	46
416701.....	61	416799.....	51	421454.....	36	421797.....	46
416703.....	62	416901.....	51	421456.....	40	421798.....	45
416704.....	62	416902.....	51	421458.....	36	421816.....	61
416706.....	62	416908.....	56	421460.....	40	421821.....	61
416707.....	62	416915.....	79	421462.....	36	421822.....	62
416709.....	62	416916.....	79	421464.....	40	421823.....	62
416710.....	62	416917.....	72	421615.....	59	421848.....	56
416712.....	62	416918.....	73	421616.....	59	421853.....	55
416713.....	62	416919.....	71	421617.....	59	421854.....	57
416715.....	62	416920.....	71	421618.....	59	421855.....	57
416718.....	55	416921.....	72	421619.....	59	421856.....	57
416719.....	55	416922.....	73	421620.....	59	421857.....	57
416721.....	55	416923.....	71	421621.....	59	421858.....	55
416724.....	57	416924.....	71	421622.....	59	421859.....	57
416725.....	57	416925.....	72	421623.....	63	421860.....	57
416728.....	57	416926.....	73	421624.....	63	421861.....	57
416729.....	57	416927.....	71	421625.....	63	421862.....	57
416732.....	57	416928.....	71	421626.....	63	421893.....	45
416733.....	57	416929.....	71	421627.....	63	422159.....	59
416736.....	57	416930.....	71	421628.....	63	422163.....	97
416737.....	57	416932.....	76	421629.....	63	422363.....	60
416738.....	57	416934.....	77	421631.....	63	422364.....	60
416739.....	57	416938.....	78	421638.....	37	422365.....	60
416740.....	57	416940.....	76	421639.....	37	422539.....	60
416741.....	57	416942.....	77	421640.....	37	422540.....	60
416742.....	57	416946.....	78	421641.....	37	422541.....	60
416743.....	57	416948.....	76	421642.....	37	422545.....	60
416744.....	57	416950.....	77	421643.....	37	422546.....	60
416745.....	57	416954.....	78	421676.....	45	422551.....	64
416746.....	57	416975.....	55	421677.....	45	422552.....	64
416747.....	57	416976.....	55	421678.....	45	422553.....	64
416748.....	57	420695.....	44	421679.....	45	422557.....	64
416749.....	57	420696.....	44	421680.....	45	422558.....	64
416750.....	57	420697.....	44	421681.....	45	422559.....	64

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND
& FISTULA
POUCHES

RESOURCES

RESOURCES

- Discontinued Products
- Patient Ostomy Support & Resources
- Quick Links Guide
- Medicare Utilization Chart
- Product Code Index
- References

Request Product Sample

422560.....	64	022756.....	66
422562.....	64	022757.....	66
422564.....	64	022758.....	66
422565.....	64	022759.....	66
423286.....	91	022760.....	66
423288.....	91	022761.....	66
423289.....	99	022762.....	66
423391.....	99	022763.....	66
423392.....	91	022764.....	66
650828.....	67	022765.....	66
650829.....	67	022766.....	66
650830.....	67	022767.....	66
650831.....	67	022768.....	66
650832.....	67	022769.....	66
650833.....	67	022770.....	66
839001.....	93	022771.....	65
839002.....	93	025433.....	65
839003.....	95	025510.....	92
839004.....	95	025542.....	96
839005.....	93	027060.....	98
839006.....	93	027062.....	98
839010.....	94	037436.....	100
839250.....	105	037439.....	92
839251.....	105	037443.....	100
839252.....	105	037444.....	92
839253.....	106	064927.....	68
839261.....	105	TR105.....	97
839262.....	105		
839264.....	106		
839265.....	106		
839266.....	106		
839267.....	105		
839268.....	106		
839280.....	106		
020917.....	88		
020922.....	88		
021712.....	95		
021715.....	95		
021901.....	95		
022355.....	78		
022356.....	78		
022750.....	66		
022751.....	66		
022752.....	66		
022753.....	66		
022754.....	66		
022755.....	66		

1. Bulkeley, J.E., McMullen, C.K., Grant, M., Wendel, C., Hornbrook, M.C. Krouse, R.S. (2018). Ongoing Ostomy Self-Care Challenges of Long-Term Rectal Cancer Survivors. Supportive Care in Cancer, 26 (11), 3933-3939.
2. Consumer survey following me+ nurse consultation, N=166, ConvaTec, Data on File.
3. Consumer survey following telehealth nurse consultation, ConvaTec & Corstrata, 2020, data on file.
4. ConvaTec Internal in vitro lab testing. OCT 0112
5. Maria Teresa Szewczyk, MD, PhD; Grazyna Majewska, RN, ETN; Mary V. Cabral, MS, FNP-BC, CWOCN-AP; and Karin Hölzel-Piontek, RN; The Effects of Using a Moldable Skin Barrier on Peristomal Skin Condition in Persons with an Ostomy: Results of a Prospective, Observational, Multinational Study, Ostomy Wound Management 2014;60(12):16–26.
6. ConvaTec Q2 2021 NPC US Marketrack Report.
7. Consumer sampling survey, N=34, ConvaTec, 2021, data on file.
8. Bosio G., Pisani F., Lucibello, L. A Proposal for Classifying Peristomal Skin Disorders: Results of a Multicenter Observational Study. OWM. 2007; 53(9): 38-43.
9. Clinical Evaluation of Sur Fit Natura Cut to Fit Skin Wafers Gillian Seaton. Technical data on file, reference document WHRI3940 TA238 2014.
10. Rolstad, BS. Principles and techniques in the use of convexity. Ostomy Wound Manage. 1996. 42(1); 24-32.
11. Pontieri, Lewis V. Basics of ostomy care. Medsurg Nurs. 2006; 15(4): 199-202.
12. Hanley, J. Convex stoma appliances: are we getting it right? British Journal of Nursing (Stoma Supplement). 2013; 22(16).
13. Turnbull, GB. The convexity controversy. Ostomy Wound Manage. 2003; 49(1): 16-www.o-wm.com/article/1203. Accessed 29 Feb, 2016.
14. DRG data 2015, in T5 markets: Japan, UK, USA, Germany, China .
15. Boyd, K. Clinical Protocols for Stoma Care – Use of convex appliances. Nursing Standard. 2004; 18(20): 37-38.
16. Burch, 2011: Management of Stoma Complications. NursingTimes.net. 10 Nov 2011.
17. Voegeli D. LBF® “No-Sting” Barrier Wipes: Skin Care Using Advanced Silicone Technology. British Journal of Nursing, 2008; Vol. 17: No. 7

Disclaimer:

The reimbursement information provided by ConvaTec is intended to provide general information relevant to coding and reimbursement of ConvaTec's products only. Coverage and payment policies for the same insurer can vary from one region to another and may change from time to time because of ongoing changes in government and insurance industry rules and regulations. Therefore, please confirm HCPCS codes with your local DME-MAC, private insurer, or Medicaid agency before processing claims. ConvaTec does not guarantee coverage or payment of its products listed herein.

TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

RESOURCES

- Discontinued Products
- Patient Ostomy Support & Resources
- Quick Links Guide
- Medicare Utilization Chart
- Product Code Index
- References

Request Product Sample

OSTOMY CARE *has met its* APP

Discover a fast and easy way to help you select the right products for patients, order and ship samples quickly, and sign patients up for the me+™ support program, all without leaving the bedside.



SCAN NOW

Scan the code or use the buttons to download the ConvaTec Ostomy Nurse Solutions App



[OstomyNurseSolutions.com](https://www.convatec.com/ostomy)



For general product information,
call 1-800-422-8811 or visit convatec.com/ostomy

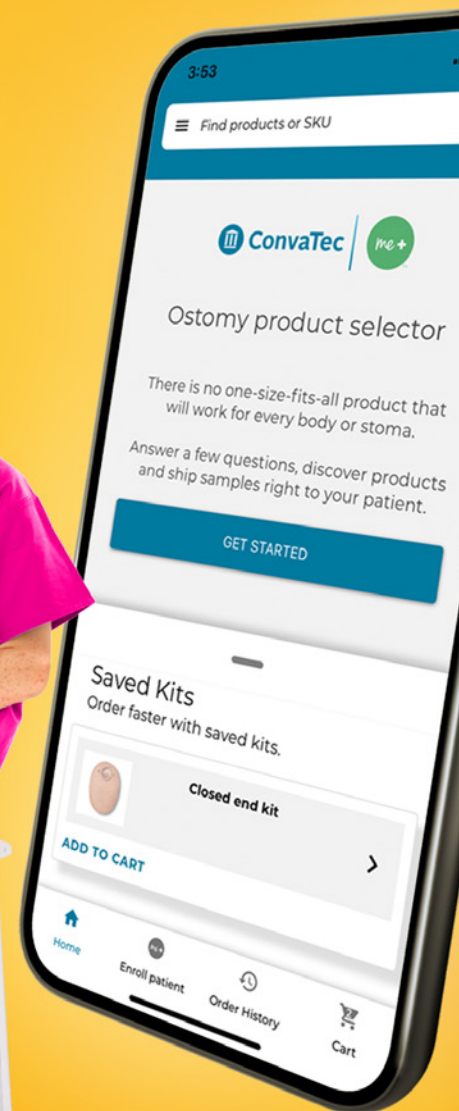


TABLE OF CONTENTS

me+™ SUPPORT

EDUCATION

TWO-PIECE SYSTEMS

ONE-PIECE SYSTEMS

ACUTE KITS

PEDIATRIC PRODUCTS

ACCESSORIES

WOUND & FISTULA POUCHES

RESOURCES

[Request Product Sample](#)