

Tracheostomy Tubes

Standard

Intended use:

This product is intended to be placed into a surgical opening (stoma) of the trachea to facilitate ventilation of the lungs. The cuff, when inflated, is intended to establish a seal between the tracheal wall and the tracheostomy tube in order to prevent patient's aspiration of substances (such as blood or vomit) and provide a means for positive pressure ventilation of the patient. This tube is attached to an external mechanical ventilator by means of a connector. The product is intended for single use only.




List of product range:

Tracheostomy Tube, Standard with and without cuff

Size			Length				Bending angle	CRD**, mm
Size	I.D., mm	O.D., mm	CLL*, mm	A, mm	B, mm	C, mm		
3.0	3.0	4.1	42.0	15.0	0	27.0	95°	-
3.5	3.5	4.8	47.0	17.0	0	30.0	95°	-
4.0	4.0	5.4	54.0	18.0	0	36.0	95°	-
4.5	9.0	6.1	54.0	19.0	0	35.0	95°	-
5.0	5.0	6.7	57.0	20.0	0	37.0	95°	-
6.0	6.0	8.0	64.0	25.5	0	38.5	95°	18
7.0	7.0	9.3	70.0	29.2	0	40.8	95°	23
8.0	8.0	10.7	81.0	35.5	0	45.5	95°	26
9.0	9.0	12.0	87.0	39.9	0	47.1	95°	29
10.0	10.0	13.3	95.5	48.4	0	47.1	95°	29

*CLL: Center line length (from flange to tip)

**CRD: Cuff resting diameter

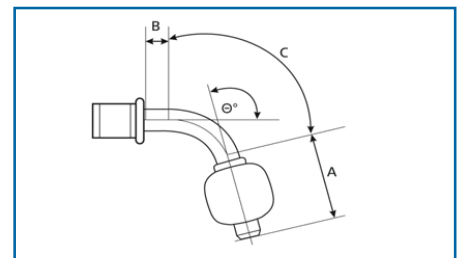

Material specifications:

- Contains latex? No
- Contains DEHP? Yes
- Tube shaft Polyvinyl chloride
- Neck flange Polyvinyl chloride
- Cuff Polyvinyl chloride
- Inflation line Polyvinyl chloride
- Connector Polypropylene
- Valve Polyvinyl chloride, Nitrile, Stainless steel
- Radiopaque line Barium sulphate steel

Certification:

- CE-marking? Yes (CE 0123)
- Product classification Class IIB (EU)
- Quality system certified to ISO 13485
- Global Medical Device Nomenclature (GMDN). 35404

Sterilisation:

- Sterile? Yes
- Sterilized by Ethylene oxide
- Shelf life 5 years

Tracheostomy Tubes

Standard

Packaging:

REF*	Shelf box		Qty., pcs.	Shipping carton		MOQ, pcs.
	Qty., pcs.	Dimensions, mm		Dimensions, mm	Gross Wt., kg	
MM62525XXX**	5	282 x 50 x 113	50	610 x 330 x 310	3.6	300
MM62545XXX**	5	282 x 50 x 113	50	610 x 330 x 310	3.6	300

* Products may not be available in all countries. Please contact your regional ConvaTec sales representative for further information.

** XXX = size, e.g 060 = 6,0

Models:

MM62525XXX**	Tracheostomy Tube with HVLP cuff, obturator and neckstrap
MM62545XXX**	Tracheostomy Tube Plain without cuff, with obturator and neckstrap

** XXX = size, e.g 060 = 6.0

Storage and handling:



With reservations for misprints and alterations in product specifications. Printed version not controlled.

Manufacturer: **Unomedical Sdn Bhd**
(15290-P)
Address: Bakar Arang Industrial Estate
08000 Sungai Petani, Kedah, Malaysia

European Representative: **Unomedical a/s**
Address: Kongevejen 2
3460 Birkerød, Denmark

For further information on these products, please contact your local ConvaTec sales representative.

TM indicates a trademark of Unomedical a/s.
The ConvaTec logo is a trademark of ConvaTec Inc.

