

Ostomy Care • Wound and Skin Care • Continence and Critical Care

2017
2018

Product Catalog

Effective **August 1, 2017**

1-800-465-6302 | www.convatec.ca

ConvaTec

Customer Relations Center:

What is it?

The Customer Relations Center (CRC) is ConvaTec's living communication hub of ostomy, wound and skin care as well as continence and critical care expertise.

For healthcare professionals:

- Clinical information and technical data on our products
- New product information
- Samples for your patients
- Patient literature on ostomy, wound and skin care, continence and critical care
- Advanced wound care expertise

For retailers:

- Technical assistance & product questions
- Samples for your customers
- Contact with Territory Managers

For people living with an ostomy:

- Free product samples
- Helpful literature
- Expert advice on how to care for their stoma

The CRC team:

Registered Nurses



1-800-465-6302

Monday to Friday from 8:00 AM to 6:00 PM (EST)

Toll-free fax: 1-877-437-1777

You can also reach us by email at convatec.canada@convatec.com
or visit our website at www.convatec.ca

TABLE OF CONTENTS

	Page
GENERAL INFORMATION (Terms and conditions of sale)	4
ORDERING INFORMATION	5
OSTOMY CARE PRODUCTS	
Natura®/ Natura®  Ostomy System (Snap Flange).....	6
Natura® Moldable Accordion Skin Barrier	13
Natura® Cut-to-fit Accordion Skin Barrier	14
Esteem <i>synergy</i> ®/ Esteem <i>synergy</i> ®  Ostomy System	20
Active Life® Ostomy System.....	25
Esteem®  Ostomy System	28
PEDIATRIC CARE PRODUCTS	
Little Ones® One-Piece Ostomy System	33
Little Ones® Two-Piece Ostomy System	34
Little Ones® Adhesive Coupling Technology™ (ACT) Ostomy System.....	35
OSTOMY SKIN CARE PRODUCTS & ACCESSORIES	
Stomahesive® Skin Barrier, Powder and Paste	36
AllKare® Wipes	36
Sillesse™ and Niltac™ Skin Protection Products	37
Diamonds® Gelling and Odor Control Sachets	37
Stomahesive® Strip.....	37
Orabase® Oral Protective Paste	37
Orahesive® Denture Adhesive (Powder)	37
Eakin Cohesive® Seal and Skin Barrier	38
Eakin Cohesive® Paste.....	NEW 38
Visi-Flow® Irrigation System	39
Night Drainage Container Set & Tube	40
Accuseal® Adapter.....	40
Belt and Tail Closures for Ostomy Pouches	40
POSTOPERATIVE SYSTEMS & KITS	
Natura®/ Natura®  Postoperative Kits	41
Esteem <i>synergy</i> ® Professional & Postoperative Kits	43
Esteem®  Postoperative Kits.....	43
ConvaTec Loop Rod.....	43
FISTULA & WOUND POUCHES	
Eakin Cohesive® Fistula & Wound Pouch.....	44
CONTINENCE AND CRITICAL CARE PRODUCTS	
Flexi-Seal® Fecal Management System (FMS).....	45
SECUREMENT & MONITORING DEVICES	
FLEXI-TRAK® Anchoring Device.....	45
Drain-Fix™ Securement Devices	45
GentleCath™ Intermittent Urinary Cathether™	46

EXUDATE MANAGEMENT

Absorbent Cover Dressings

AQUACEL® Foam Hydrofiber® Dressing	47
AQUACEL® Ag Foam Hydrofiber® Dressing	48
Foam Lite™ ConvaTec Dressing	48
AQUACEL® Surgical Hydrofiber® Wound Dressing	49
AQUACEL® Ag Surgical Hydrofiber® Wound Dressing	49
CarboFlex® Odor Control Dressing	49

Absorbent Filler & Antimicrobial Dressings

AQUACEL® Ag+ Extra Dressing (for Biofilm Management).....	50
AQUACEL® Ag+ Ribbon Dressing (for Biofilm Management)	50
AQUACEL® Hydrofiber® Wound Dressing	51
AQUACEL® Ag Hydrofiber® Wound Dressing	51
AQUACEL® Extra™ Hydrofiber® Wound Dressing	52
AQUACEL® Ag Extra™ Hydrofiber® Wound Dressing	52
AQUACEL® Ag Burn Hydrofiber® Dressing	53
KALTOSTAT® Calcium Alginate Wound Dressing	53

MOISTURE RETENTIVE DRESSINGS

DuoDERM Signal® Hydrocolloid Dressing	54
DuoDERM® Extra Thin CGF® Hydrocolloid Dressing	54
DuoDERM® CGF® Hydrocolloid Dressing	55
DuoDERM® CGF® Border Hydrocolloid Dressing	55
DuoDERM® Hydroactive® Hydrocolloid Dressing	55

WOUND FILLERS AND WOUND HYDRATION

DuoDERM® Hydroactive® Gel	56
---------------------------------	----

COMPRESSION

SurePress® High Compression Bandage	56
SurePress® Absorbent Padding.....	56

SKIN CARE PRODUCTS

Aloe Vesta® & Sensi-Care® Bathing and Perineal Cleansing.....	57
Aloe Vesta® & Sensi-Care® Moisturizing and Protective Skin Care.....	58

NEGATIVE PRESSURE WOUND THERAPY

Avelle™ Negative Pressure Wound Therapy System	NEW 59
--	---------------

GLOSSARY	60
----------------	----

ALPHABETICAL INDEX	63
--------------------------	----

NUMERICAL INDEX	66
-----------------------	----

USEFUL RESOURCES	72
------------------------	----

AUTHORIZATION TO RETURN FORM	73
------------------------------------	----

PACKAGE CODE NUMBERING SYSTEM	74
-------------------------------------	----

® Registered trademark of ConvaTec Inc.

™ Trademark of ConvaTec Inc.

Licensed for use in Canada by ConvaTec Inc.

Eakin Cohesive® products are trademarks of T.G. Eakin Limited.

GENERAL INFORMATION

Note: In order to avoid misunderstanding and to protect the interest of our customers, we respectfully request careful reading of the "Terms and Conditions of Sale", as outlined. Non-compliance with these terms may result in suspension of your account.

TERMS AND CONDITIONS OF SALE

PRICES AND TERMS OF SALE

Prices and terms listed herein supersede those set forth on all previous price lists and all contracts and agreements that are expired. Due to market fluctuation, prices shown are subject to change without notice. **RESTRICTION OF PURCHASE & SALES:** Purchase of ConvaTec products must be done through the Canadian order desk at 1-866-987-1471 or indirectly through a distributor who purchases **directly from ConvaTec Canada**. These products are **not for distribution or resale outside of Canada** or to anyone who - to the customers' knowledge - is susceptible to this practice.

PAYMENT TERMS

TERMS ARE NET 30 DAYS FROM DATE OF INVOICE.

Electronic Fund Transfers (EFT) are possible, please contact the Order Desk.

CHANGES IN PRICES

In all cases, merchandise will be invoiced at prices in effect on the date the order is received. Any pricing discrepancy **MUST** be reported within 60 days.

MINIMUM SHIPMENTS AND TRANSPORTATION COSTS

All orders amounting to a minimum of \$250 before tax for retail accounts and for health care facilities will be shipped prepaid via seller's choice of carrier. Any orders of value less than the required minimum will be subject to a \$50 service charge. Extra shipping costs incurred because of special services requested by a customer will be added to the invoice.

SHIPMENTS

Title to merchandise passes to the customer upon delivery of the goods to the customer.

SHORTAGE, DAMAGE OR LOSS IN SHIPMENT

Merchandise released for shipment leaves the depot in good order. Claims for damage, defacement or loss occurring in transit should be filed promptly with the transportation company. Customers should contact the Order Desk for processing of such claims and in supplying proof of delivery, providing these requests are received within 21 DAYS following the date of shipment. **ANY COSTS INCURRED WILL BE PASSED ALONG TO THE CUSTOMER IF PROOF OF DELIVERY IS ESTABLISHED.**

RETURNED GOODS POLICY

All returns **MUST** have written authorization from the Order Desk (telephone: 1-866-987-1471). This authorization will be confirmed by the representative and the credit will only be issued after inspection of the returned goods by ConvaTec.

RETURNED GOODS WILL BE CREDITED AT ZERO IF NOT ACCOMPANIED BY A RETURN AUTHORIZATION NUMBER OR DOCUMENT. Order desk may ask for copy of purchase order for compliance purposes.

- All merchandise returned for credit must be prepaid and received in the original unopened containers. Partially full containers, those from which the seal, carton or label have been defaced or removed, will not be credited under any circumstances. Goods for return should be packaged carefully before sending to ConvaTec; products and containers received with any signs of damage will not be credited. **ConvaTec will not ship back products that have been sent to us, as such, please ensure you respect the terms and conditions before sending products to ConvaTec.**
- Products shipped in error by ConvaTec Canada will be fully credited if returned within 21 days of invoice, freight collect.
- Products ordered in error will be fully credited **if returned within 21 days of invoice, freight prepaid.**
- There will be a restocking service charge equal to 15% of the invoice price or the lowest price indicated for that product in the current price list.
- All returned goods must have been purchased in the last 12 months.
- All returned goods must have at least 6 months of shelf life remaining.
- All returned goods must be received by ConvaTec within 21 days after the written authorization was issued.
- Products that have been deleted from the line for a period of 6 months or more after date of deletion notification from ConvaTec, including special offers, will not be accepted for credit or exchange.
- Merchandise that has been involved in a fire, flood, sacrifice or bankruptcy sale or returned by other than original purchaser will not be accepted for credit or exchange.

CONSUMER SATISFACTION

We are quality driven and we want the consumers to be completely satisfied with our products. In the event of a product complaint from a consumer, the retailer should refer this consumer to the ConvaTec Customer Relations Center at 1-800-465-6302.

TAXES

Provincial sales tax and/or GST, HST, QST will be charged, if applicable.

THE ABOVE TERMS CANNOT BE MODIFIED BY OUR REPRESENTATIVES.

ConvaTec

Ordering Information

PRODUCT ORDER DESK - TEL: 1-866-987-1471

- FAX: 1-866-987-1472

TO REACH CONVATEC CUSTOMER RELATIONS CENTER:

TEL: 1-800-465-6302

FAX: 1-877-437-1777

E-MAIL: convatec.canada@convatec.com

WEBSITE: www.convatec.ca

MERCHANDISE RETURNS

Returns of merchandise must be sent to:

ConvaTec Canada

RETURNS DEPARTMENT

1375 Newton

Boucherville, QC

J4B 5H2

ACCOUNT PAYMENT:

In order to expedite and facilitate the payment of your account, we maintain a Lock Box System. Please forward your cheque to the Post Office Box below. This is for cheques only.

TORONTO

ConvaTec Canada Ltd.

PO Box 15000

Station A

Toronto, ON

M5W 1C1



1425, TransCanada, Suite 100 - Dorval , QC H9P 2V3 Tel: (514)-822-5900 Fax : (514)-822-5930

www.convatec.ca



Natura® Two-Piece Ostomy System

Stoma Sizing Usage Chart

Skin Barriers with Flange														
Stoma Size (Diameter)	13 mm ½"	16 mm ⅝"	19 mm ¾"	22 mm ⅞"	25 mm 1"	28 mm 1⅛"	32 mm 1¼"	35 mm 1⅜"	38 mm 1½"	45 mm 1¾"	50 mm 2"	57 mm 2¼"	64 mm 2½"	70-90 mm 2¾"-3½"
Flange Size	32 mm 1¼"			38 mm 1½"		45 mm 1¾"		57 mm 2¼"			70 mm 2¾"		100 mm 4"	

NOTE: Always choose a skin barrier with with flange ring at least 13 mm (1/2 inch) **larger** than the diameter of the stoma.

Moldable Skin Barriers with Flange											
Stoma Size (Diameter)	13 mm 1/2"	16 mm 5/8"	19 mm 3/4"	22 mm 7/8"	25 mm 1"	28 mm 1 1/8"	33 mm 1 5/16"	38 mm 1 1/2"	45 mm 1 3/4"	50 mm 2"	56 mm 2 1/8"
Flange Size	45 mm 1 3/4" Small				45 mm 1 3/4" Medium			57 mm 2 1/4" Large		70 mm 2 3/4" Extra-Large	
										Oval	

Moldable Skin Barriers with Accordion Flange				
Stoma Size (Diameter)	13 mm 1/2"	22 mm 7/8"	33 mm 1 5/16"	45 mm 1 3/4"
Flange Size	57 mm 2 1/4" Small		57 mm 2 1/4" Medium	70 mm 2 3/4" Large



Cut-to-Fit Skin Barriers with Accordion Flange									
Stoma Size (Diameter)	13 mm 1/2"	16 mm 5/8"	19 mm 3/4"	21 mm 13/16"	25 mm 1"	28 mm 1 1/8"	33 mm 1 5/16"	38 mm 1 1/2"	45 mm 1 3/4"
Flange Size	45 mm 1 3/4"			57 mm 2 1/4"			70 mm 2 3/4"		

NOTE: Always choose a skin barrier with flange ring at least 13 mm (1/2 inch) **larger** than the diameter of the stoma.

Cut-to-fit Convex Skin Barriers with Accordion Flange									
Stoma Size (Diameter)	13 mm 1/2"	16 mm 5/8"	19 mm 3/4"	22 mm 7/8"	25 mm 1"	28 mm 1 1/8"	35 mm 1 3/8"	38 mm 1 1/2"	48 mm 1 7/8"
Flange Size	45 mm 1 3/4"			57 mm 2 1/4"			70 mm 2 3/4"		



Natura® Low-Pressure Adaptor



The Natura® Low-Pressure Adaptor converts the Natura® ostomy system into a low-pressure solution post-surgery. The product is an insert between the Natura® pouch and skin barrier to create enough clearance for a nurse to insert fingers and snap the pouch coupling ring into place without applying force directly to the patient's abdomen, which may cause pain or otherwise harm the post-op area.

Use of the Natura® Low-Pressure Adaptor can be discontinued when abdominal tenderness subsides. The pouch and skin barrier will connect together without using the Natura® Low-Pressure Adaptor.

The Natura® Low-Pressure Adaptor fits with the Natura® skin barrier and pouch 45 mm, 57 mm, 70 mm and 100 mm coupling sizes. It is a disposable one time use product to be changed with each skin barrier change.

TIP!—
Ensure that the flange ring sizes of the Natura® skin barrier, Adaptor and Pouch are the same.

Reminder —
Make sure to attach the Adaptor to the skin barrier before applying it to the skin.



401993 - 45 mm (1¾") flange
Box of 10
401994 - 57 mm (2¼") flange
Box of 10
401995 - 70 mm (2¾") flange
Box of 10
401996 - 100 mm (4") flange
Box of 10'

When Using Natura® Disposable Convex Inserts with Natura® Stomahesive® or Durahesive® Skin Barriers

Stoma Size (Diameter)	13 mm ½"	16 mm ⅝"	19 mm ¾"	22 mm 7⁄8"	25 mm 1"	28 mm 1⅛"	32 mm 1¼"	35 mm 1⅜"
Flange Size	38 mm 1½"			45 mm 1¾"			57 mm 2¼"	
Use Convex Insert	404006	404007	404008	404009	404010	404011	404012	404013


Natura® Disposable Convex Inserts*



Natura® Disposable Convex Inserts are designed to snap into the flange ring adding convexity to the standard skin barrier for improved management of retracted stomas or other conditions which call for additional barrier-to-skin contact.


	FLANGE SIZE	INSERT INNER DIAMETER	STOMA SIZE
404006	38 mm (1½")	19 mm (¾")	13 mm (½")
404007		22 mm (⅞")	16 mm (⅝")
404008		25 mm (1")	19 mm (¾")
404009	45 mm (1¾")	28 mm (1⅛")	22 mm (⅞")
404010		32 mm (1¼")	25 mm (1")
404011		35 mm (1⅜")	28 mm (1⅛")
404012	57 mm (2¼")	38 mm (1½")	32 mm (1¼")
404013		41 mm (1⅝")	35 mm (1⅜")

Box of 5 convex inserts

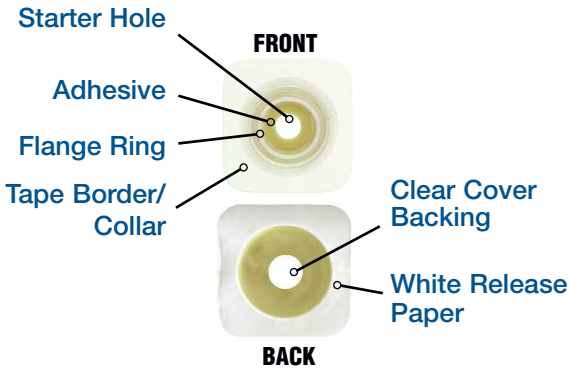

* Natura® inserts should only be used with Natura® Skin Barriers.

Natura® Stomahesive® Skin Barrier (with low profile flange)

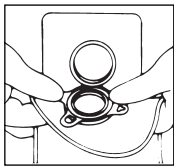

Stomahesive® Skin Barriers are ideal for people whose stoma output is semi-formed to formed stool. This reliable skin barrier with low-profile flange provides excellent security because of a unique formula that enables the skin barrier to adhere to both dry and moist skin. The outer protective coating helps keep the skin barrier in place to provide added security when showering, bathing or swimming. Natura® Stomahesive® Skin Barriers are specifically designed for use with Natura® pouches.

Did you know?

1972: The Stomahesive® Skin Barrier was the 1st advanced skin barrier to improve ostomy care.

Natura® Stomahesive® Flexible Skin Barrier (with low profile flange and white collar), Cut-to-Fit

Natura® Flexible Skin Barriers with flange and tape collar have rounded edges for a more comfortable fit, low-profile flange and split release paper backing for easy application.



10 cm x 10 cm (approx. 4" x 4")

125257 – 32 mm (1¼")

125258 – 38 mm (1½")

125259 – 45 mm (1¾")

Box of 10 skin barriers*

12.5 cm x 12.5 cm (approx. 5" x 5")

125260 – 57 mm (2¼")

125261 – 70 mm (2¾")

Box of 10 skin barriers*

Natura® Stomahesive® Flexible Skin Barrier (with low profile flange and tan collar), Cut-to-Fit



10 cm x 10 cm (approx. 4" x 4")

125262 – 32 mm (1¼")

125263 – 38 mm (1½")

125264 – 45 mm (1¾")

Box of 10 skin barriers*

12.5 cm x 12.5 cm (approx. 5" x 5")

125265 – 57 mm (2¼")

125266 – 70 mm (2¾")

Box of 10 skin barriers*

Natura® Stomahesive® Solid Skin Barrier (with low profile flange), Cut-to-Fit

Stomahesive® Solid Skin Barrier with no tape collar.



10 cm x 10 cm (approx. 4" x 4")

401573 – 32 mm (1¼")

401574 – 38 mm (1½")

401575 – 45 mm (1¾")

401576 – 57 mm (2¼")

Box of 10 skin barriers*

12.5 cm x 12.5 cm (approx. 5" x 5")

401577 – 70 mm (2¾")

Box of 10 skin barriers*

15 cm x 15 cm (approx. 6" x 6")

401906 – 100 mm (4")

Box of 5 skin barriers*

* Stoma measuring guides are provided in each box of skin barriers.
Additional guides are available upon request from ConvaTec Customer Relations Center.



Natura® Stomahesive® Flexible Skin Barrier
(with flange and tan collar), Pre-Cut



Protective Skin Barrier with pre-cut stoma sizes

The pre-cut opening eliminates the need for cutting the skin barrier, accommodating those people with limited manual dexterity.

10 cm x 10 cm (approx. 4" x 4")

FLANGE SIZE	STOMA SIZE
125267 - 45 mm (1¾")	13 mm (½")
125268 - 45 mm (1¾")	16 mm (⅝")
125269 - 45 mm (1¾")	19 mm (¾")
125270 - 45 mm (1¾")	22 mm (7⁄8")
125271 - 45 mm (1¾")	25 mm (1")
125272 - 45 mm (1¾")	28 mm (1⅛")
125273 - 45 mm (1¾")	32 mm (1¼")
125274 - 45 mm (1¾")	35 mm (1⅜")

Box of 10 skin barriers



12.5 cm x 12.5 cm (approx. 5" x 5")

FLANGE SIZE	STOMA SIZE
125275 - 57 mm (2¼")	38 mm (1½")
125276 - 57 mm (2¼")	41 mm (1⅝")
125277 - 57 mm (2¼")	45 mm (1¾")

Box of 10 skin barriers

Did you know? _____
If the product description of a skin barrier contains the word "flexible", it means the skin barrier has a tape edge.

Durahesive® Skin Barrier maintains an effective protective seal longer than conventional skin barriers. A “turtlenecking” effect helps to prevent leakage around the stoma and skin breakdown, for increased patient comfort and wear time. This solid, flexible or convex skin barrier, with tan or white collar, is designed with a low-profile flange for wearing comfort.

Did you know?

The solid Durahesive® Skin Barrier is more flexible than the solid Stomahesive® Skin Barrier.

Natura® Durahesive® Flexible Skin Barrier
(with low profile flange and white collar), Cut-to-Fit


10 cm x 10 cm (approx. 4" x 4")

413159 – 32 mm (1 1/4")

413160 – 38 mm (1 1/2")

413161 – 45 mm (1 3/4")

Box of 10 skin barriers*

12.5 cm x 12.5 cm (approx. 5" x 5")

413162 – 57 mm (2 1/4")

Box of 10 skin barriers*

12.5 cm x 12.5 cm (approx. 5" x 5")

413163 – 70 mm (2 3/4")

Box of 10 skin barriers*

Natura® Durahesive® Flexible Skin Barrier
(with low profile flange and tan collar), Cut-to-Fit


10 cm x 10 cm (approx. 4" x 4")

413164 – 32 mm (1 1/4")

413165 – 38 mm (1 1/2")

413166 – 45 mm (1 3/4")

Box of 10 skin barriers*

12.5 cm x 12.5 cm (approx. 5" x 5")

413167 – 57 mm (2 1/4")

Box of 10 skin barriers*

12.5 cm x 12.5 cm (approx. 5" x 5")

413168 – 70 mm (2 3/4")

Box of 10 skin barriers*

Natura® Durahesive® Solid Skin Barrier
(with low profile flange), Cut-to-Fit


Durahesive® Solid Skin Barrier with no tape border.

10 cm x 10 cm (approx. 4" x 4")

413153 – 32 mm (1 1/4")

413154 – 38 mm (1 1/2")

413155 – 45 mm (1 3/4")

413156 – 57 mm (2 1/4")

Box of 10 skin barriers*

12.5 cm x 12.5 cm (approx. 5" x 5")

413157 – 70 mm (2 3/4")

Box of 10 skin barriers*

15 cm x 15 cm (approx. 6" x 6")

401905 – 100 mm (4")

Box of 5 skin barriers*

* Stoma measuring guides are provided in each box of skin barriers.
Additional guides are available upon request from ConvaTec Customer Relations Center.



Natura® Durahesive® Flexible Skin Barrier with Convex-It®
(with low profile flange and white collar), Pre-Cut



The Durahesive® Flexible Skin Barrier with Convex-It® maintains an effective protective seal longer than conventional skin barriers. It provides deep convexity for extra security for retracted or difficult-to-manage stomas. A turtlenecking effect helps to prevent leakage around the stoma and skin breakdown, for increased patient comfort and wear-time.

11.5 cm x 11.5 cm (approx. 4" x 4")

FLANGE SIZE	STOMA SIZE
413177 - 45 mm (1¾")	13 mm (½")
413178 - 45 mm (1¾")	16 mm (⅝")
413179 - 45 mm (1¾")	19 mm (¾")
413180 - 45 mm (1¾")	22 mm (⅞")
413181 - 45 mm (1¾")	25 mm (1")
413182 - 45 mm (1¾")	28 mm (1⅛")
413183 - 45 mm (1¾")	32 mm (1¼")
413184 - 45 mm (1¾")	35 mm (1⅜")

Box of 10 skin barriers

12.5 cm x 12.5 cm (approx. 5" x 5")

FLANGE SIZE	STOMA SIZE
413185 - 57 mm (2¼")	38 mm (1½")
413186 - 57 mm (2¼")	41 mm (1⅝")
413187 - 57 mm (2¼")	45 mm (1¾")
413188 - 57 mm (2¼")	50 mm (2")

Box of 10 skin barriers




Natura® Skin Barrier with ConvaTec Moldable Technology™

The Natura® Moldable skin barriers allow the stoma opening to be molded to the shape and size of the stoma without the need for scissors.

Natura® Skin Barrier with ConvaTec Moldable Technology™ with white Hydrocolloid Collar Adhesive (HC) or white Acrylic Collar Adhesive (AC)



The Natura® skin barrier with ConvaTec Moldable Technology™ provides a worry-free fit because the molded adhesive with Rebounding Memory Technology™ gently hugs individual stoma contours, creating a secure, personalized fit. The Durahesive® technology top-side layer swells or “turtlenecks” up the stoma improving fit and providing a gapless seal.

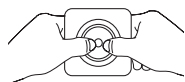

Durahesive®			
Collar		10 cm x 10 cm (approx. 4.5" x 4.5")	
HC	AC	FLANGE SIZE	STOMA SIZE
411800	413417 - Small	45 mm (1 3/4")	13 mm - 22 mm (1/2" - 7/8")
411802	413418 - Medium	45 mm (1 3/4")	22 mm - 33 mm (7/8" - 1 1/4")
		12.5 cm x 12.5 cm (approx. 5" x 5")	
		FLANGE SIZE	STOMA SIZE
411804	413419 - Large	57 mm (2 1/4")	33 mm - 45 mm (1 1/4" - 1 3/4")
411806	413420 - Extra Large	70 mm (2 3/4")	45 mm - 56 mm (1 3/4" - 2 1/8") - oval
Box of 10 skin barriers			


Stomahesive®			
Collar		10 cm x 10 cm (approx. 4" x 4")	
HC	AC	FLANGE SIZE	STOMA SIZE
411801	413421 - Small	45 mm (1 3/4")	13 mm - 22 mm (1/2" - 7/8")
411803	413422 - Medium	45 mm (1 3/4")	22 mm - 33 mm (7/8" - 1 1/4")
		12.5 cm x 12.5 cm (approx. 5" x 5")	
		FLANGE SIZE	STOMA SIZE
411805	413423 - Large	57 mm (2 1/4")	33 mm - 45 mm (1 1/4" - 1 3/4")
411807	413424 - Extra Large	70 mm (2 3/4")	45 mm - 56 mm (1 3/4" - 2 1/8") - oval
Box of 10 skin barriers			

Natura® Durahesive® Convex Flexible Skin Barrier with ConvaTec Moldable Technology™ (with low profile flange and white collar) Acrylic Collar Adhesive (AC)



Creates a customized stomal opening for retracted, flush, round, oval or irregular stomas. The Natura® Durahesive® Convex Flexible Skin Barrier with ConvaTec Moldable Technology™ must be used together with the Natura® pouches.

**MOLDABLE
CONVEXITY**

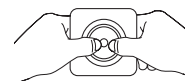

11.5 cm x 11.5 cm (approx. 4" x 4")			
		FLANGE SIZE	STOMA SIZE
404592 - Small		45 mm (1 3/4")	13 mm - 22 mm (1/2" - 7/8")
404593 - Medium		45 mm (1 3/4")	22 mm - 33 mm (7/8" - 1 1/4")
Box of 10 skin barriers			
12.5 cm x 12.5 cm (approx. 5" x 5")			
		FLANGE SIZE	STOMA SIZE
404594 - Large		57 mm (2 1/4")	33 mm - 45 mm (1 1/4" - 1 3/4")
Box of 10 skin barriers			



Natura® Skin Barriers with ACCORDION flange


The Accordion flange makes pouch application much easier and more comfortable - The pressure is off! With ample room for coupling, it's ideal for tender post-op abdomens. In the "down" position, it is low profile, extremely flexible and conforms to body contours. Compatible with our state-of-the-art pouches.



Natura® Stomahesive® Moldable Skin Barrier with ACCORDION Flange, Hydrocolloid Collar Adhesive (HC), tan



Stomahesive® Skin Barriers are ideal for people whose stoma output is semi-formed to formed stool.

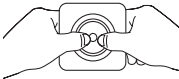

11.5 cm x 11.5 cm (approx. 4.5" x 4.5")

	FLANGE SIZE	STOMA SIZE
421033 - Small	57 mm (2¼")	13 mm - 22 mm (½" - 7⁄8")
421034 - Medium	57 mm (2¼")	22 mm - 33 mm (7⁄8" - 1¼")

12 cm x 12 cm (approx. 5" x 5")


	FLANGE SIZE	STOMA SIZE
421035 - Large	70 mm (2¾")	33 mm - 45 mm (1¼" - 1¾")

Box of 10 skin barriers

Natura® Durahesive® Moldable Skin Barrier with ACCORDION Flange, Hydrocolloid Collar Adhesive (HC), tan



Durahesive® Skin Barrier maintains an effective protective seal longer than conventional skin barriers.

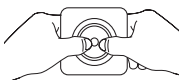

11.5 cm x 11.5 cm (approx. 4.5" x 4.5")

	FLANGE SIZE	STOMA SIZE
421039 - Small	57 mm (2¼")	13 mm - 22 mm (½" - 7⁄8")
421040 - Medium	57 mm (2¼")	22 mm - 33 mm (7⁄8" - 1¼")

12 cm x 12 cm (approx. 5" x 5")

	FLANGE SIZE	STOMA SIZE
421041 - Large	70 mm (2¾")	33 mm - 45 mm (1¼" - 1¾")

Box of 10 skin barriers



Take the pressure off of pouch application



Natura® Stomahesive® Cut-to-Fit Skin Barrier with ACCORDION flange, Acrylic Collar Adhesive (AC), white



10 cm x 10 cm (approx. 4" x 4")

421456 - Small 45 mm (1³/₄") for stomas up to 21mm (1³/₁₆")

421460 - Medium 57 mm (2¹/₄") for stomas up to 33mm (1¹/₄")

11.5 cm x 11.5 cm (approx. 4.5" x 4.5")

421464 - Large 70 mm (2³/₄") for stomas up to 45mm (1³/₄")

Box of 10 skin barriers*



421468 - X Large 100 mm (4") for stomas up to 70mm (2³/₄")

Box of 10 skin barriers

Natura® Durahesive® Cut-to-Fit Skin Barrier with ACCORDION flange, Acrylic Collar Adhesive (AC), white



10 cm x 10 cm (approx. 4" x 4")

421454 - Small 45 mm (1³/₄") for stomas up to 21mm (3/₁₆")

421458 - Medium 57 mm (2¹/₄") for stomas up to 33mm (1¹/₄")

11.5 cm x 11.5 cm (approx. 4.5" x 4.5")

421462 - Large 70 mm (2³/₄") for stomas up to 45mm (1³/₄")

Box of 10 skin barriers*



421466 - X Large 100 mm (4") for stomas up to 70mm (2³/₄")

Box of 10 skin barriers

Natura® Durahesive® Convex Cut-to-Fit Skin Barrier with ACCORDION flange, Acrylic Collar Adhesive (AC) or Hydrocolloid Adhesive (HC)



Collar

AC	HC
----	----

10 cm x 10 cm (approx. 4" x 4")

421638 **421639 - Small**

45 mm (1³/₄") for stomas up to 22mm (7/₁₆")

421640 **421641 - Medium**

57 mm (2¹/₄") for stomas up to 35mm (1³/₈")

11.5 cm x 11.5 cm (approx. 4.5" x 4.5")

421642 **421643 - Large**

70 mm (2³/₄") for stomas up to 48mm (1⁷/₈")

Box of 10 skin barriers



"UP" position



"DOWN" position



The two-piece system provides versatility and can be removed or interchanged with other Natura® pouch types without removal of the skin barrier. Available with one or two-sided comfort panel.

Natura® Drainable Pouch, without Filter

Natura® drainable pouches are composed of a lightweight plastic material that is odor-proof and quiet under clothing.



Small, 25.4 cm (10"), Opaque, 2-sided comfort panel, without Filter

Box of 10

- 401505 – 32 mm (1¼")
- 401506 – 38 mm (1½")
- 401507 – 45 mm (1¾")
- 401508 – 57 mm (2¼")
- 401509 – 70 mm (2¾")

Box of 20

- 404030 – 32 mm (1¼")
- 404031 – 38 mm (1½")
- 404032 – 45 mm (1¾")
- 404033 – 57 mm (2¼")
- 404034 – 70 mm (2¾")

Box of 10 pouches (2 tail closures per box)



Standard, 30.5 cm (12"), Opaque, 2-sided comfort panel, without Filter

Box of 10

- 401500 – 32 mm (1¼")
- 401501 – 38 mm (1½")
- 401502 – 45 mm (1¾")
- 401503 – 57 mm (2¼")
- 401504 – 70 mm (2¾")

Box of 20

- 404025 – 32 mm (1¼")
- 404026 – 38 mm (1½")
- 404027 – 45 mm (1¾")
- 404028 – 57 mm (2¼")
- 404029 – 70 mm (2¾")

Box of 10 pouches (2 tail closures per box)

Box of 20 pouches (2 tail closures per box)



Standard, 30.5 cm (12"), Opaque, 1-sided comfort panel, without Filter

- 404014 – 32 mm (1¼")
- 404015 – 38 mm (1½")
- 404016 – 45 mm (1¾")
- 404017 – 57 mm (2¼")
- 404018 – 70 mm (2¾")

Box of 20 pouches (2 tail closures per box)



Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel, without Filter

- 401510 – 32 mm (1¼")
- 401511 – 38 mm (1½")
- 401512 – 45 mm (1¾")
- 401513 – 57 mm (2¼")
- 401514 – 70 mm (2¾")

Box of 10 pouches (2 tail closures per box)



Large, 35.5 cm (14"), Transparent, 1-sided comfort panel, without Filter

For loop ostomy, the management of abdominal fistulas and draining wounds.

- 401515 – 100 mm (4")

Box of 10 pouches (2 tail closures per box)

Did you know?

The comfort panel of our drainable pouches is made of a fine plastic mesh that feels silky and comfortable on the skin.

Natura® Drainable Pouch with Filter

Natura® Drainable Pouch with filter is designed to stay flatter, help prevent manual air release and control odor.

Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel, with Filter

411264 – 38 mm (1½")

411265 – 45 mm (1¾")

411266 – 57 mm (2¼")

411267 – 70 mm (2¾")

Box of 10 pouches (2 tail closures per box)
(Waterproof filter tabs provided in each box)

Standard, 30.5 cm (12"), Opaque, 2-sided comfort panel, with Filter

411490 – 38 mm (1½")

411491 – 45 mm (1¾")

411492 – 57 mm (2¼")

411493 – 70 mm (2¾")

Box of 20 pouches (2 tail closures per box)
(Waterproof filter tabs provided in each box)

Small, 25.4 cm (10"), Opaque, 2-sided comfort panel, with Filter

411289 – 38 mm (1½")

411290 – 45 mm (1¾")

411291 – 57 mm (2¼")

Box of 20 pouches (2 tail closures per box)
(Waterproof filter tabs provided in each box)

Natura® High Output Drainable Pouch with Filter

Used for the management of high stomal output due to short bowel syndrome or other causes.

Standard, 30.5 cm (12"), Opaque, 1-sided comfort panel with replaceable filter

401557 – 45 mm (1¾")

401558 – 57 mm (2¼")

401559 – 70 mm (2¾")

Box of 5 pouches
(1 tail closure and 1 secondary drainage adapter per box)

Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel with replaceable filter

420695 – 45 mm (1¾")

420696 – 57 mm (2¼")

420697 – 70 mm (2¾")

Box of 5 pouches
(1 tail closure and 1 secondary drainage adapter per box)

TIP!

The filters in both drainable and closed-end pouches need to be protected by a filter tab when the client is bathing.

OSTOMY CARE PRODUCTS

Natura®/Natura® Ostomy System

Natura®/ Natura®

LATEX

Natura® Drainable Pouch with InvisiClose® (Clipless) Tail Closure, without Filter

The InvisiClose® clipless tail closure offers double security, ease of drainage, comfort and cleanliness. The easy-to-open outlet end strips and tail design help control flow and reduce splashing. The locking mechanism is enhanced by an additional security flap which is designed to provide comfort and protect skin and clothing from locking material.



Regular

Standard, 30.5 cm (12"), Opaque, 1-sided comfort panel, without Filter

413332 – 32 mm (1¼")

413334 – 45 mm (1¾")

413335 – 57 mm (2¼")

413336 – 70 mm (2¾")

Box of 20 pouches

Regular

Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel, without Filter

413312 – 45 mm (1¾")

413313 – 57 mm (2¼")

413314 – 70 mm (2¾")

Box of 10 pouches

Natura® Drainable pouches with InvisiClose® and Lock-It Pocket™, without Filter



Regular

Standard, 30.5 cm (12") Transparent, 1-sided comfort panel, without Filter

416418 – 45 mm (1¾")

416421 – 57 mm (2¼")

416424 – 70 mm (2¾")

416472 – 100 mm (4")

Box of 10 pouches

Natura® Drainable Pouch and InvisiClose® (Clipless) Tail Closure, with Filter

The Natura® Drainable Pouch with Filter and InvisiClose® clipless tail closure offers all the benefits of the InvisiClose® clipless tail closure plus the added comfort, convenience and confidence of ConvaTec's filter.



Regular

Standard, 30.5 cm (12"), Opaque, 2-sided comfort panel with Filter

411309 – 38 mm (1½")

411310 – 45 mm (1¾")

411311 – 57 mm (2¼")

411312 – 70 mm (2¾")

Box of 20 pouches
(Waterproof filter tabs provided in each box)

Regular

Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel with Filter

411360 – 45 mm (1¾")

411361 – 57 mm (2¼")

411362 – 70 mm (2¾")

Box of 10 pouches
(Waterproof filter tabs provided in each box)

Natura® Drainable pouches InvisiClose® and Lock-It Pocket™, with Filter



Regular

Standard, 30.5 cm (12"), Opaque, 2-sided comfort panel with Filter

416415 – 38 mm (1½")

416417 – 45 mm (1¾")

416420 – 57 mm (2¼")

416423 – 70 mm (2¾")

Box of 10 pouches
(Waterproof filter tabs provided in each box)

Regular

Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel with Filter

416416 – 45 mm (1¾")

416419 – 57 mm (2¼")

416422 – 70 mm (2¾")

Box of 10 pouches
(Waterproof filter tabs provided in each box)

Natura®/Natura® Closed-End Pouch

Closed-end pouches are designed for one-time usage, which makes them ideal for everyday use, special occasions or for intimate moments. Eliminates the need to manipulate a tail closure and emptying a pouch.

Natura® Closed-End Pouch, without Filter



Regular

Standard, opaque, 20.3 cm (8"), 2-sided comfort panel, without Filter

401521 – 38 mm (1½")

401522 – 45 mm (1¾")

401523 – 57 mm (2¼")

401524 – 70 mm (2¾")

Box of 30 pouches

Natura®

Standard, opaque, 20.3 cm (8"), 2-sided comfort panel, without Filter

416401 – 32 mm (1¼")

416404 – 38 mm (1½")

416407 – 45 mm (1¾")

416410 – 57 mm (2¼")

416413 – 70 mm (2¾")

Box of 30 pouches



Small, opaque, 15.2 cm (6"), 2-sided comfort panel, without Filter

416402 – 32 mm (1¼")

416405 – 38 mm (1½")

416408 – 45 mm (1¾")

416411 – 57 mm (2¼")

Box of 30 pouches



Natura® Closed-End Pouch with Filter



Regular

Standard, opaque, 20.3 cm (8"), 2-sided comfort panel with Filter

401525 – 38 mm (1½")

401526 – 45 mm (1¾")

401527 – 57 mm (2¼")

401528 – 70 mm (2¾")

Box of 30 pouches

(Waterproof filter tabs provided in each box)

Natura®

Standard, opaque, 20.3 cm (8"), 2-sided comfort panel with Filter

416400 – 32 mm (1¼")

416403 – 38 mm (1½")

416406 – 45 mm (1¾")

416409 – 57 mm (2¼")

416412 – 70 mm (2¾")

Box of 30 pouches

(Waterproof filter tabs provided in each box)



Natura® Mini-Pouch without Filter



The Natura® Mini-Pouch is a small, closed-end pouch designed for individuals desiring a very low profile or children requiring a smaller pouch. Ideal design to be used during social settings, travelling, swimming or during intimate moments.

Opaque, 12.5 cm (5"), 1-sided comfort panel, without Filter

401529 – 32 mm (1¼")

401530 – 38 mm (1½")

401531 – 45 mm (1¾")

401532 – 57 mm (2¼")

Box of 20 mini-pouches

Natura® Flange Cap with Filter



Made with carbon cloth filter, built-in gas release, a soft inner comfort panel and absorbent pad. Designed specifically for individuals who irrigate and do not wear a pouch between irrigations, or for short-term use during times of stomal inactivity.

Opaque, 1-sided comfort panel with Filter

401909 – 45 mm (1¾")

Box of 25 flange caps

(Waterproof filter tabs provided in each box)



Natura® Urostomy Pouch

Designed for use with any Natura® Stomahesive® or Durahesive® Skin Barrier.
 All Natura® Urostomy Pouches have a **1-sided comfort panel**.

Natura® Urostomy Pouch with Accuseal® Tap

Natura® Urostomy Pouch with Accuseal® Tap features a swivel drainage valve and a stem plug for added security. An Accuseal® adapter is provided to facilitate connection to a night drainage system. Natura® Urostomy Pouch with Accuseal® Tap is constructed with an effective anti-reflux valve and its shape is contoured for quick and easy emptying. Natura® Urostomy Pouches with low profile Accuseal® Tap provide handling ease, comfort and a more discreet appearance.

TIP!
 The Natura® urostomy pouches with an Accuseal® tap can be connected directly to the night drainage system. No connector is required.



Standard, 25.4 cm (10"), Transparent

- 401542 – 32 mm (1¼")
- 401543 – 38 mm (1½")
- 401544 – 45 mm (1¾")
- 401545 – 57 mm (2¼")
- 401546 – 70 mm (2¾")

Box of 10 pouches (2 Accuseal® adapters per box)



Small, 23 cm (9"), Opaque

- 401547 – 32 mm (1¼")
- 401548 – 38 mm (1½")
- 401549 – 45 mm (1¾")
- 401550 – 57 mm (2¼")

Box of 10 pouches (2 Accuseal® adapters per box)



Standard, 25.4 cm (10"), Opaque

- 401551 – 32 mm (1¼")
- 401552 – 38 mm (1½")
- 401553 – 45 mm (1¾")
- 401554 – 57 mm (2¼")
- 401555 – 70 mm (2¾")

Box of 10 pouches (2 Accuseal® adapters per box)

Natura® Urostomy Pouch with Bendable Tap

Natura® Urostomy Pouches with Bendable Tap provide handling ease, comfort and a discreet appearance.



Small, 20.3 cm (8"), Transparent

- 401538 – 32 mm (1¼")
- 401539 – 38 mm (1½")
- 401540 – 45 mm (1¾")
- 401541 – 57 mm (2¼")

Box of 10 pouches (2 clear plastic adapters per box)



Standard, 23 cm (9"), Transparent

- 401533 – 32 mm (1¼")
- 401534 – 38 mm (1½")
- 401535 – 45 mm (1¾")
- 401536 – 57 mm (2¼")
- 401537 – 70 mm (2¾")

Box of 10 pouches (2 clear plastic adapters per box)

Natura® Urostomy Pouch with Soft Tap



Natura® Urostomy Pouch with Soft Tap incorporates a new pouch shape that has been redesigned with a soft body-side comfort panel. Pouch baffling evenly distributes urine for a slimmer profile and less tugging. An upgraded Flexible Soft Tap is designed for security and added comfort.

Standard, 10" Natura® Urostomy Pouch w/ Soft Tap, Transparent

- 413435 – 32 mm (1¼")
- 413436 – 38 mm (1½")
- 413437 – 45 mm (1¾")
- 413438 – 57 mm (2¼")
- 413439 – 70 mm (2¾")

Box of 10 pouches (2 clear plastic adapters per box)



Esteem synergy® Two-Piece Ostomy System

Skin Barrier - To - Pouch Cross Reference

The table below shows the stoma sizes that each barrier can accommodate and shows the pouch collar diameter to be used with that skin barrier.

Flat Skin Barriers

Size	Label Colour	Reference Number to Ensure Proper Sizing	Maximum Stoma Size
Small	Blue	35	35 mm (1 ³ / ₈ ")
Medium	Green	48	48 mm (1 ⁷ / ₈ ")
Large	Yellow	61	61 mm (2 ³ / ₈ ")
Extra Large	Peach	89	89 mm (3 ¹ / ₂ ")

Durahesive® Moldable Convex Flexible Skin Barrier

Size	Label Colour	Reference Number to Ensure Proper Sizing	Maximum Stoma Size
Small	Blue	35	22 mm (7/ ₈ ")
Medium	Green	48	33 mm (1 ¹ / ₄ ")
Large	Yellow	61	45 mm (1 ³ / ₄ ")

Pouches

Size	Label Colour	Reference Number to Ensure Proper Sizing	Pouch Collar Diameter
Small	Blue	35	35 mm (1 ³ / ₈ ")
Medium	Green	48	48 mm (1 ⁷ / ₈ ")
Large	Yellow	61	61 mm (2 ³ / ₈ ")
Extra Large	Peach	89	89 mm (3 ¹ / ₂ ")



Esteem synergy® Skin Barrier with Transparent Landing Zone

The Esteem synergy® ostomy system utilizes an Adhesive Coupling Technology™ that connects the pouch to the skin barrier and eliminates the need for flange rings that have to be snapped or locked together. Esteem synergy® ostomy system combines the benefits of one and two-piece systems. The clear, flat, flexible plastic “Landing Zone” on the skin barrier is square and is slightly larger in size compared to the round adhesive ring on the pouch. This makes the application process easy – precise alignment is unnecessary.

Esteem synergy® Stomahesive® Flexible Skin Barrier (White Collar), Cut-to-Fit



- 405456 – Small**, for stomas up to 35 mm (1³/₈")
- 405457 – Medium**, for stomas up to 48 mm (1⁷/₈")
- Box of 10 skin barriers*
- 405458 – Large**, for stomas up to 61 mm (2³/₈")
- Box of 10 skin barriers*
- Note: starter hole 13 mm (1/2")*

Esteem synergy® Stomahesive® Solid Skin Barrier, Cut-to-Fit



- 405467 – Small**, for stomas up to 35 mm (1³/₈")
- 405468 – Medium**, for stomas up to 48 mm (1⁷/₈")
- Box of 10 skin barriers*
- 403948 – Large**, for stomas up to 61 mm (2³/₈")
- Box of 10 skin barriers*
- 405483 – Extra-Large**, for stomas up to 89 mm (3¹/₂")
- Box of 10 skin barriers*
- Note: starter hole 13 mm (1/2")*

Esteem synergy® Stomahesive® Flexible Skin Barrier (White Collar), Pre-cut



- 405474 – For stoma size 19 mm (3/4") - pouch 35 mm (1³/₈")**
- 405475 – For stoma size 22 mm (7/8") - pouch 35 mm (1³/₈")**
- 405476 – For stoma size 25 mm (1") - pouch 35 mm (1³/₈")**
- 405477 – For stoma size 28 mm (1¹/₈") - pouch 35 mm (1³/₈")**
- 405478 – For stoma size 32 mm (1¹/₄") - pouch 35 mm (1³/₈")**
- 405479 – For stoma size 35 mm (1³/₈") - pouch 35 mm (1³/₈")**
- Box of 10 skin barriers*
- 405480 – For stoma size 38 mm (1¹/₂") - pouch 48 mm (1⁷/₈")**
- 405482 – For stoma size 45 mm (1³/₄") - pouch 48 mm (1⁷/₈")**
- Box of 10 skin barriers*

* Stoma measuring guides are provided in each box of skin barriers. Additional guides are available upon request from ConvaTec Customer Relations Center.

Esteem synergy® Skin Barrier with ConvaTec Moldable Technology™

The Esteem synergy® System also offers skin barriers with ConvaTec Moldable Technology™ that hug the individual contours of your stoma to allow for a secure, personalised fit without the need for scissors.



Esteem synergy® Stomahesive® Flat Skin Barrier with ConvaTec Moldable Technology™, white hydrocolloid tape collar



411657 – Small

For stomas size 13mm to 22mm (1/2" to 7/8") - pouch 48mm (1 7/8")

411658 – Medium

For stomas size 22mm to 33mm (7/8" to 1 1/4") - pouch 48mm (1 7/8")

411659 – Large


For stomas size 33mm to 45mm (1 1/4" to 1 3/4") - pouch 61mm (2 3/8")

Box of 10 skin barriers

Esteem synergy® Durahesive® Convex Flexible Skin Barrier with ConvaTec Moldable Technology™, white acrylic tape collar



**MOLDABLE
CONVEXITY**

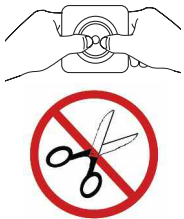

This innovative mold-to-fit convex skin barrier easily shapes to create a customized stoma opening to conform to irregularly-shaped stomas. It provides moderate to deep convexity and hugs the base of the stoma for optimum security and support. It eliminates the need for scissors and is perfect for people living with an ostomy who have dexterity and/or vision problems.

409268 – Small

For stomas size 13 mm to 22 mm (1/2" to 7/8") - pouch 35 mm (1 3/8")

409269 – Medium

For stomas size 22 mm to 33 mm (7/8" to 1 1/4") - pouch 48 mm (1 7/8")

409270 – Large

For stomas size 33 mm to 45 mm (1 1/4" to 1 3/4") - pouch 61 mm (2 3/8")

Box of 10 skin barriers*

Esteem synergy® Drainable Pouch with Filter (with Cotton-Like Comfort Panel)

Combines the benefits of one and two-piece systems. Esteem synergy® drainable pouches offer the patient simplicity, discretion and comfort. For added security, the transparent landing zone helps determine when the pouch should be changed. Easy alignment allows for a consistently secure seal. The pouches feature an advanced filter system designed to absorb and neutralize a broad range of odor-causing agents. A durable soft-touch cotton-like backing is designed to dry easily after showering and be comfortable to wear throughout the day. The Esteem synergy® pouch is easy to apply and easy to remove.

Esteem synergy® Straight Drainable Pouch, with Filter




Esteem synergy® combines the low profile, comfort and flexibility of a one-piece system with the lifestyle versatility and skin friendliness of a two-piece system.

Extra Large, 35.5 cm (14"), Transparent, with Filter

405406 – Straight - 89 mm (3 1/2")

Box of 10 pouches (10 straight tail closures per box)

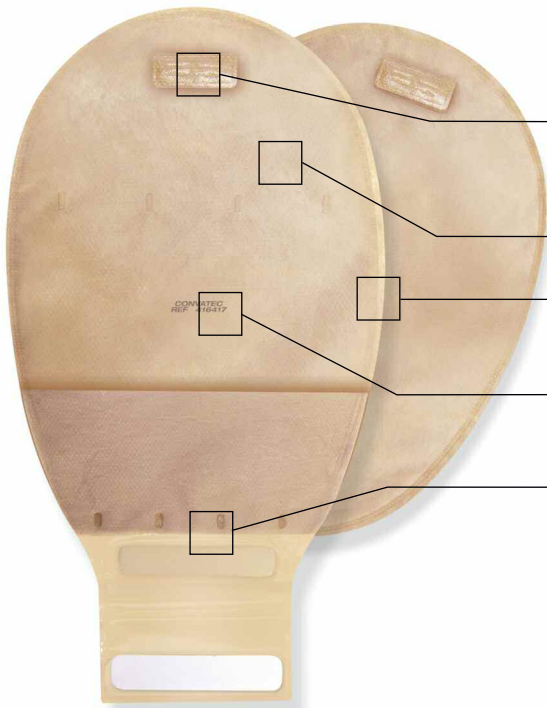

Esteem synergy® Features:

- Smart Filter**
 - Reduces clogging allowing it to work longer, keeping the pouch flat – while substantially reducing odor. Filter covers available.
- Comfortable Material**
 - Soft and smooth Comfort Panel that resists water.
- Streamlined Shape**
 - Innovative streamlined shape is rounded for added comfort and discretion.
- Easy to Read Code**
 - Product code number is printed right on the pouch for quick identification.
- Discreet Closure**
 - Lock-it pocket™ on InvisiClose® pouches hides away the flapping tail to provide the security and comfort of a closed-end pouch. Closure is even easier to use with less snaps and fewer steps.

Esteem synergy® Drainable pouches with InvisiClose® and Lock-It Pocket™



Esteem synergy®  drainable pouches offer simplicity, discretion and comfort. The InvisiClose® clipless tail closure offers double security, ease of drainage, comfort and cleanliness. The easy-to-open outlet and widened tail design help control flow and reduce splashing. The locking mechanism is enhanced by a Lock-It Pocket™ that hides away the tail to provide the security and comfort of a closed-end pouch.



Standard, 30.5 cm (12”), Opaque, 2-sided comfort panel with Filter

- 416796** – 35 mm (1³/₈”)
- 416799** – 48 mm (1⁷/₈”)
- 416902** – 61 mm (2³/₈”)

Box of 10 pouches

Standard, 30.5 cm (12”), Transparent, 1-sided comfort panel with Filter

- 416795** – 35 mm (1³/₈”)
- 416798** – 48 mm (1⁷/₈”)
- 416901** – 61 mm (2³/₈”)

Box of 10 pouches

Standard, 30.5 cm (12”),Transparent, 1-sided comfort panel without Filter

- 416797** – 35 mm (1³/₈”)
- 416900** – 48 mm (1⁷/₈”)
- 416903** – 61 mm (2³/₈”)

Box of 10 pouches



Esteem synergy® Closed-end pouches with Filter

Closed-end pouches are designed for one-time usage, which makes them ideal for everyday use, special occasions or for intimate moments. Eliminate the need to manipulate a tail closure and emptying a pouch.



Standard, 20.3 cm (8''), Opaque, 2-sided comfort panel with Filter

416775 – 35 mm (1³/₈'')

416779 – 48 mm (1⁷/₈'')

416783 – 61 mm (2³/₈'')

Box of 30 pouches



Small, 15.2 cm (6''), Opaque, 2-sided comfort panel with Filter,

416777 – 35 mm (1³/₈'')

416781 – 48 mm (1⁷/₈'')

Box of 30 pouches

Standard, 20.3 cm (8''), Transparent, 1-sided comfort panel with Filter

416774 – 35 mm (1³/₈'')

416778 – 48 mm (1⁷/₈'')

416782 – 61 mm (2³/₈'')

Box of 30 pouches

Small, 15.2 cm (6''), Transparent, 1-sided comfort panel with Filter

416776 – 35 mm (1³/₈'')

416780 – 48 mm (1⁷/₈'')

Box of 30 pouches

Esteem synergy® Urostomy Pouch with Accuseal® Tap

The Esteem synergy® Urostomy Drainable pouch is comfortable, thin, and easy to use. Its anti-reflux valve keeps urine away from the stoma. Its rounded corners are gentle against the skin.



Standard, 26.2 cm (10.3''), Transparent, 1-sided comfort panel

405451 – Small - 35 mm (1³/₈'')

Box of 10 pouches (2 Accuseal® adapters per box)

Active Life® One-Piece System

Active Life® appliances are made of quiet, soft odor-proof film. A microporous, built-in hypo-allergenic adhesive support collar is for utmost security and comfort.

Active Life® Drainable Pouch with Stomahesive® Flexible Skin Barrier, Pre-Cut, without Filter

Active Life® Drainable Pouch features pre-cut stoma openings, a non-clinging “fabric” comfort panel and a wide tail.

Small, opaque, 25.4 cm (10”), without Filter

(Pre-cut stoma openings in 1 size)

22751 – 25 mm (1”)

Box of 10 pouches (1 tail closure per box)

Small, opaque, 25.4 cm (10”), without Filter

(Pre-cut stoma openings in 5 sizes)

125331 – 19 mm ($\frac{3}{4}$ ”)

125333 – 32 mm (1 $\frac{1}{4}$ ”)

125334 – 38 mm (1 $\frac{1}{2}$ ”)

125335 – 45 mm (1 $\frac{3}{4}$ ”)

125336 – 50 mm (2”)

Box of 20 pouches (1 tail closure per box)

Standard, opaque, 30.5 cm (12”), without Filter

(Pre-cut stoma openings in 6 sizes)

125338 – 19 mm ($\frac{3}{4}$ ”)

125339 – 25 mm (1”)

125340 – 32 mm (1 $\frac{1}{4}$ ”)

125341 – 38 mm (1 $\frac{1}{2}$ ”)

125342 – 45 mm (1 $\frac{3}{4}$ ”)

125343 – 50 mm (2”)

Box of 20 pouches (1 tail closure per box)

Standard, transparent, 30.5 cm (12”), without Filter

(Pre-cut stoma openings in 6 sizes)

22764 – 19 mm ($\frac{3}{4}$ ”)

22765 – 25 mm (1”)

22766 – 32 mm (1 $\frac{1}{4}$ ”)

22767 – 38 mm (1 $\frac{1}{2}$ ”)

22768 – 45 mm (1 $\frac{3}{4}$ ”)

22769 – 50 mm (2”)

Box of 10 pouches (1 tail closure per box)



Active Life® Drainable Pouch with Stomahesive® Flexible Skin Barrier, Cut-to-Fit, without Filter



The cut-to-fit style provides the ability to cut the opening to fit the shape of the stoma; ideal for irregular stomas or during the postoperative period. Stoma opening can be enlarged to any shape and size up to 64 mm (2 $\frac{1}{2}$ ”).

22771 – Standard, transparent, 30.5 cm (12”), without Filter

Box of 10 pouches (1 tail closure per box)

125330 – Standard, transparent, 30.5 cm (12”), without Filter

Box of 20 pouches (1 tail closure per box)

Active Life® Drainable Pouch with Flexible Durahesive Plus® Skin Barrier, Cut-to-Fit, without Filter



Ideal for irregular stomas or during postoperative period. Provides the benefits of Active Life® Drainable Pouch with Flexible Durahesive Plus® Skin Barrier for extended wear time. Stoma opening can be enlarged to any shape and size up to 64 mm (2 $\frac{1}{2}$ ”).

400598 – Standard, transparent, 30.5 cm (12”), without Filter

400599 – Standard, opaque, 30.5 cm (12”), without Filter

Box of 10 pouches (1 tail closure per box)



Active Life® Drainable Pouch with Stomahesive® Solid Skin Barrier, Cut-to-Fit, without Filter



Stoma opening can be enlarged to any shape and size up to 100 mm (4").

650803 – Grand, transparent, 35.5 cm (14"), without Filter

Box of 10 pouches
(10 DuoLock® tail closures per box)

The benefits of the Active Life® *One-Piece* Closed-End system are offered with a specially formulated Stomahesive® Skin Barrier, designed for those who change pouches frequently. A filter is built in for odor control and a non-clinging, soft "fabric" comfort panel is incorporated.

Active Life® Closed-End Pouch with Filter and Stomahesive® Flexible Skin Barrier, Cut-to-Fit



Stoma opening can be enlarged up to 50 mm (2")

650408 – Opaque, 20.3 cm (8") with Filter

650422 – Transparent, 20.3 cm (8") with Filter

Box of 30 pouches
(Waterproof filter tabs provided in each box)

Active Life® Closed-End Pouch with Filter and Stomahesive® Solid Skin Barrier, Pre-Cut



The pre-cut style utilizes the large oval skin barrier that eliminates the need for both a support collar and a tape edge. Pre-cut stoma openings also eliminate the need for cutting, trimming or shaping the skin barrier.

Opaque, 20.3 cm (8") with Filter

(Pre-cut stoma openings in 5 sizes)

175769 – 25 mm (1")

175770 – 32 mm (1¼")

175771 – 38 mm (1½")

175772 – 45 mm (1¾")

175773 – 50 mm (2")

Box of 15 pouches
(Waterproof filter tabs provided in each box)

Opaque, 20.3 cm (8") with Filter

(Pre-cut stoma openings in 5 sizes)

413142 – 25 mm (1")

413143 – 32 mm (1¼")

413144 – 38 mm (1½")

413145 – 45 mm (1¾")

413146 – 50 mm (2")

Box of 60 pouches
(Waterproof filter tabs provided in each box)

Active Life® Stoma Cap with Filter and Stomahesive® Solid Skin Barrier, Cut-to-Fit



Discreet stoma coverage for individuals who irrigate and do not wear a pouch between irrigations, or for short-term use during times of stomal inactivity. Includes a soft inner comfort panel and filter. The cut-to-fit style provides the ability to cut the opening to fit the shape of the stoma.

Stoma opening can be enlarged up to 50 mm (2")

175611 – Opaque, diameter of 10 cm (approx. 4") with Filter

Box of 30 stoma caps
Waterproof filter tabs provided in each box)

Active Life® One-Piece Convex System

The Active Life® One-Piece Convex system features the long lasting protection of the Durahesive® skin barrier in an integral convex faceplate. The pre-cut stoma opening eliminates the need to cut, trim, or shape.

Active Life® Urostomy Pouch with Bendable Tap and Convex Durahesive® Flexible Skin Barrier, Pre-Cut

CONVEX



Standard, 23 cm (9"), Transparent

- 125361 – 13 mm ($\frac{1}{2}$ ") stoma opening size
- 125362 – 16 mm ($\frac{5}{8}$ ") stoma opening size
- 125363 – 19 mm ($\frac{3}{4}$ ") stoma opening size
- 125364 – 22 mm ($\frac{7}{8}$ ") stoma opening size
- 125365 – 25 mm (1") stoma opening size
- 125366 – 28 mm (1 $\frac{1}{8}$ ") stoma opening size
- 125367 – 32 mm (1 $\frac{1}{4}$ ") stoma opening size
- 125368 – 35 mm (1 $\frac{3}{8}$ ") stoma opening size
- 125369 – 38 mm (1 $\frac{1}{2}$ ") stoma opening size

Box of 10 pouches (1 plastic adapter per box)

Did you know?

All the pouches on this page have built-in belt tabs.

Active Life® Drainable Pouch with Convex Durahesive® Flexible Skin Barrier, Pre-Cut

CONVEX



Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel

- 175777 – 19 mm ($\frac{3}{4}$ ") stoma opening size
- 175778 – 22 mm ($\frac{7}{8}$ ") stoma opening size
- 175779 – 25 mm (1") stoma opening size
- 175780 – 28 mm (1 $\frac{1}{8}$ ") stoma opening size

Box of 5 pouches (1 tail closure per box)

Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel

- 125352 – 19 mm ($\frac{3}{4}$ ") stoma opening size
- 125353 – 22 mm ($\frac{7}{8}$ ") stoma opening size
- 125354 – 25 mm (1") stoma opening size
- 125355 – 28 mm (1 $\frac{1}{8}$ ") stoma opening size
- 125356 – 32 mm (1 $\frac{1}{4}$ ") stoma opening size
- 125357 – 35 mm (1 $\frac{3}{8}$ ") stoma opening size
- 125358 – 38 mm (1 $\frac{1}{2}$ ") stoma opening size
- 125359 – 45 mm (1 $\frac{3}{4}$ ") stoma opening size

Box of 10 pouches (1 tail closure per box)

Active Life® Urostomy Pouch with Accuseal® Tap and Durahesive® Flexible Skin Barrier, Pre-Cut



The Active Life® Urostomy Pouch with Durahesive® Skin Barrier is designed for people who want the security and skin protection found in Durahesive® but who do not require convexity. It is constructed of odor-proof, transparent film and has a soft comfort panel on the body side of the pouch. The built-in anti-reflux valve resists the back flow of urine. The Accuseal® tap is designed for easy manipulation and rapid emptying. An Accuseal® adapter is supplied for connection to night drainage.

Standard, 24.8 cm (9 $\frac{3}{4}$ "), Transparent, 1-sided comfort panel

- 650828 – 19 mm ($\frac{3}{4}$ ") stoma opening size
- 650829 – 22 mm ($\frac{7}{8}$ ") stoma opening size
- 650830 – 25 mm (1") stoma opening size
- 650831 – 28 mm (1 $\frac{1}{8}$ ") stoma opening size
- 650832 – 32 mm (1 $\frac{1}{4}$ ") stoma opening size
- 650833 – 38 mm (1 $\frac{1}{2}$ ") stoma opening size

Box of 10 pouches (1 Accuseal® adapter per box)

Active Life® Urostomy Pouch with Bendable Tap and Stomahesive® Skin Barrier, Cut-to-Fit



Designed for people who need the security and skin protection found in Stomahesive®. It is constructed of odor-proof transparent film with a built-in anti-reflux valve to resist back flow of urine. The Bendable Tap provides a low-profile design.

Stoma opening can be enlarged to any shape and size up to 45 mm (1 $\frac{3}{4}$ ").

Standard, 23 cm (9"), Transparent

- 64927 – 19 mm ($\frac{3}{4}$ ") starter hole stoma opening size

Box of 10 pouches (1 plastic adapter per box)

Esteem® ⊕ One-piece System



Esteem® ⊕ system with Stomahesive® Technology is a one-piece system which is secure, flexible and most of all, kind to the skin.

SECURITY through durable adhesion

A skin-friendly skin barrier based on Stomahesive® technology that adheres on first contact with the skin and can be safely and easily removed, helping damaged and irritated skin to recover.

COMFORT both in fitting and during wear through unique design and choice of material.

ConvaTec worked with nurses and patients when designing the Esteem® ⊕ system using a computerized pouch design process. The ergonomic shape which contours with the body's natural shape makes Esteem® ⊕ comfortable and discreet to wear.



Esteem® ⊕ One-piece System

- **Easy** to tailor to the stoma
- Starter opening **for both round and oval** stomas
- Measuring options **up to 70mm**
- **Optimal choice** of material and design
- **Soft**, durable backing, rounded, flowing form



Esteem® ⊕ Convex system


Offers all the benefits of one-piece Esteem® system even when convexity is required

- choice of sizes

No cutting is required on Esteem® ⊕ Convex which is available in a wide range of pre-cut pouches. The rigid convex skin barrier, made from Durahesive® Technology, is designed to provide extended wear time.

Esteem® ⊕ Moldable One-piece System

The convenience of a one piece plus proven skin protection – it's the perfect combination. Elastic-like seal “rebounds” to fit the size and shape of stoma, providing a secure, snug fit. Just roll, mold and hold the adhesive to accommodate any stoma shape and size. No cutting, no stretching, no guesswork. Thin flexible design is made to move with your patients. The state-of-the-art pouch offers a soft comfort panel, easy-view window and the InvisiClose® clip-less tail closure with Lock-It Pocket™.

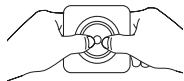**Esteem® ⊕ Drainable System, Moldable, Stomahesive®****Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel, with Filter****413515 – Medium** - 20-30 mm**413516 – Large** - 30-40 mm**413517 – Extra Large** - 40-50 mm

Box of 10 pouches

Standard, 30.5 cm (12"), Transparent, 1-sided comfort panel, without Filter**413518 – Medium** - 20-30 mm**413519 – Large** - 30-40 mm


Box of 10 pouches

Standard, 30.5 cm (12"), Opaque, 2-sided comfort panel, with Filter**413521 – Medium** - 20-30 mm**413522 – Large** - 30-40 mm**413523 – Extra Large** - 40-50 mm

Box of 10 pouches

Esteem® ⊕ Closed System, Moldable, Stomahesive®**Standard, 20.3 cm (8"), Opaque, 2-sided comfort panel, with Filter****413509 – Medium** - 20-30 mm**413510 – Large** - 30-40 mm**413511 – Extra Large** - 40-50 mm

Box of 30 pouches

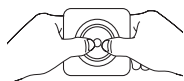



Esteem® Solid Stomahesive® Closed-end Pouches with Filter, Cut-to-Fit



Standard, 20.3 cm (8''), Transparent, 1-sided comfort panel, with Filter

416700 – 60/70 mm (2²/₅'' by 2³/₄'')

Box of 30 pouches

Standard, 20.3 cm (8''), Opaque, 2-sided comfort panel, with Filter

416701 – 60/70 mm (2²/₅'' by 2³/₄'')

Box of 30 pouches



Small, 15.2 cm (6''), Opaque, 2-sided comfort panel, with Filter

416702 – 60/70 mm (2²/₅'' by 2³/₄'')

Box of 30 pouches

Esteem® Solid Stomahesive® Closed-end Pouches with Filter, Pre-Cut



Standard, 20.3 cm (8''), Transparent, 1-sided comfort panel, with Filter

416703 – 25 mm (1'')

416706 – 30 mm (1³/₁₆'')

416709 – 35 mm (1³/₈'')

416712 – 40 mm (1⁹/₁₆'')

Box of 30 pouches

Standard, 20.3 cm (8''), Opaque, 2-sided comfort panel, with Filter

416704 – 25 mm (1'')

416707 – 30 mm (1³/₁₆'')

416710 – 35 mm (1³/₈'')

416713 – 40 mm (1⁹/₁₆'')

416715 – 50 mm (2'')

Box of 30 pouches



Small, 15.2 cm (6''), Opaque, 2-sided comfort panel, with Filter

416705 – 25 mm (1'')

416711 – 35 mm (1³/₈'')

416714 – 40 mm (1⁹/₁₆'')

Box of 30 pouches



Esteem® Solid Stomahesive®, Drainable Pouches with InvisiClose®, Lock-It Pocket™, Cut-to-Fit



Standard, 30.5 cm (12''), Transparent, 1-sided comfort panel, with Filter

416718 – 60/70 mm (2²/₅" by 2³/₄")

Box of 10 pouches



Standard, 30.5 cm (12''), Opaque, 2-sided comfort panel, with Filter

416719 – 60/70 mm (2²/₅" by 2³/₄")

Box of 10 pouches

Standard, 30.5 cm (12''), Transparent, 1-sided comfort panel, without Filter

416721 – 60/70 mm (2²/₅" by 2³/₄")

Box of 10 pouches

Standard, 30.5 cm (14''), Transparent, 1-sided comfort panel, without Filter

416908 – 100 mm (4")

Box of 10 pouches

Esteem® Solid Stomahesive®, Drainable Pouches with InvisiClose®, Lock-It Pocket™, Filter, Pre-Cut



Standard, 30.5 cm (12''), Opaque, 2-sided comfort panel, with Filter

416725 – 25 mm (1")

416729 – 30 mm (1³/₁₆")

416733 – 35 mm (1³/₈")

416737 – 40 mm (1⁹/₁₆")

Box of 10 pouches



Standard, 30.5 cm (12''), Transparent, 1-sided comfort panel, with Filter

416724 – 25 mm (1")

416728 – 30 mm (1³/₁₆")

416732 – 35 mm (1³/₈")

416736 – 40 mm (1⁹/₁₆")

Box of 10 pouches

Esteem® flexible Durahesive®, Drainable Pouches with InvisiClose®, Lock-It Pocket™, Filter, Cut-to-Fit



Standard, 30.5 cm (12''), Opaque, 1-sided comfort panel, with Filter

416975 – 19/64mm (3/4" by 2 1/2")

Box of 10 pouches

Standard, 30.5 cm (12''), Transparent, 1-sided comfort panel, with Filter

416976 – 19/64mm (3/4" by 2 1/2")

Box of 10 pouches

Esteem® Durahesive® Flexible skin barrier with Accuseal Tap, Cut-to-Fit

NEW

421631– 13/45mm (1/2" - 1 3/4")

Box of 10 pouches



Esteem® Convex Flexible Durahesive® Drainable Pouches with InvisiClose®, Lock-It Pocket™, Filter, Pre-Cut



Standard, 30.5 cm (12''), Opaque, 2-sided comfort panel with Filter

- 416739 – 19 mm ($\frac{3}{4}$ ")
- 416741 – 22 mm ($\frac{7}{8}$ ")
- 416743 – 25 mm (1")
- 416745 – 28 mm ($1\frac{1}{8}$ ")
- 416747 – 32 mm ($1\frac{1}{4}$ ")
- 416749 – 35 mm ($1\frac{3}{8}$ ")
- 416751 – 38 mm ($1\frac{1}{2}$ ")

Box of 10 pouches

Standard, 30.5 cm (12''), Transparent, 1-sided comfort panel with Filter

- 416738 – 19 mm ($\frac{3}{4}$ ")
- 416740 – 22 mm ($\frac{7}{8}$ ")
- 416742 – 25 mm (1")
- 416744 – 28 mm ($1\frac{1}{8}$ ")
- 416746 – 32 mm ($1\frac{1}{4}$ ")
- 416748 – 35 mm ($1\frac{3}{8}$ ")
- 416750 – 38 mm ($1\frac{1}{2}$ ")

Box of 10 pouches

Esteem® Flex Solid Convex

As the latest addition to our range of one-piece ostomy solutions, the Esteem® Flex Convex system combines the comfort and freedom of flexibility with the firmness of convexity.

Key benefits include:

- **Moves with you** – flexible convexity fits to body contours for comfort and easy positioning
- **Designed for skin health and security** – hydrocolloid skin barrier is made to stay firmly in place, helping to minimise leakage. Soft tapered edges also help resist snagging on clothing
- **Comfort and discretion** – soft quiet fabric
- **Trust and reassurance** – filter on closed and drainable options helps minimise odours and ballooning, while the anti-reflux valve on our urostomy pouches helps prevent back flow
- **Peace of mind** – split-fabric, easy-view window for accurate positioning and observation
- **Simple secure fastening** – the Safe Seal clipless closure on drainable options helps make cleaning and emptying easy
- **Support and security** – soft flexible belt loops available for added security



Drainable with filter

- 421615 - V1 - Cut-to-fit, 20mm - 43mm
- 421619 - V2 - Cut-to-fit, 20mm - 35mm
- 421622 - V3 - Cut-to-fit, 20mm - 25mm

Box of 10 pouches

Urostomy

- 421623 - V1 - Cut-to-fit, 10mm - 43mm
- 421626 - V2 - Cut-to-fit, 10mm - 35mm
- 421628 - V3 - Cut-to-fit, 10mm - 25mm

Box of 10 pouches

Closed End with filter

- 421612 - V1 - Cut-to-fit, 20mm - 43mm
- 421613 - V2 - Cut-to-fit, 20mm - 35mm
- 421614 - V3 - Cut-to-fit, 20mm - 25mm

Box of 30 pouches

THREE PLATEAU SIZES

The Esteem® Flex Convex system is available in 3 plateau sizes to accommodate a range of stomas





Little Ones® One-Piece Ostomy System*

Little Ones® One-Piece Urostomy Pouch with Solid Stomahesive®
Skin Barrier, Cut-to-Fit



The Little Ones® One-Piece ostomy system provides easy, one-step application and removal with a Bendable drainage Tap for easy emptying. An anti-reflux valve helps prevent urine from pooling at the base of the stoma which extends wear time.

Stoma opening can be enlarged to any shape up to 25 mm (1").

Transparent, approx. 15 cm (6")

20917 – For stoma sizes 8 mm to 25 mm ($\frac{5}{16}$ " to 1")

Box of 15 pouches (2 night drainage adapters per box)

Little Ones® One-Piece Drainable Pouch with Solid Stomahesive®
Skin Barrier, Cut-to-Fit



The skin barrier on the Little Ones® ostomy system provides similar skin protective qualities as those offered by Stomahesive®, but is modified to be more gentle. This skin barrier is thinner, more flexible and easy to remove, minimizing potential irritation to a young child's skin.

Stoma opening can be enlarged to any shape and size, up to 50 mm (2").

Transparent, approx. 15 cm (6")

20922 – For stoma sizes 8 mm to 50 mm ($\frac{5}{16}$ " to 2")

Box of 15 pouches (16 individual tail closures per box)

Little Ones® One-Piece Drainable Pouch with InvisiClose®
(Clipless) Tail Closure, Cut-to-Fit

To fit babies smaller than 4.5kg/10lbs



Trusted Stomahesive® technology assists in providing a safe, comfortable fit. This system offers an integrated cutting guide and no starter hole to allow for customization. The InvisiClose® outlet facilitates easy drainage, yet provides two levels of closure security you can trust.

Extra Small, approx. 11 cm (4"), Transparent, 1-sided printed comfort panel

411633 – For stomas up to 23 mm ($\frac{7}{8}$ ")

Box of 10 pouches

* Little Ones® measuring guides are available in each box of skin barriers.



Little Ones® Two-Piece Ostomy System*

Little Ones® Two-Piece Stomahesive® Flexible Skin Barrier (with Flange and White Collar), Cut-to-Fit



Smaller dimensions to fit a young child's body comfortably. Soft, pliable, adhesive collar applies easily, anchors securely, dries quickly after bathing. The Stomahesive® Skin Barrier offers excellent protection to the skin. The skin barrier's flexibility contours nicely to the child's body, allowing flexibility of movement.

Approx. 8 cm x 8 cm (3¼" x 3¼")

401925 – 32 mm (1¼")

Box of 5 skin barriers

Little Ones® Two-Piece Drainable Pouch



Two-piece system permits easy pouch removal and easy repositioning to ensure unobstructed draining even while the child is asleep. Wide-tail opening for faster, easier emptying.

Transparent, approx. 15 cm (6")

401927 – 32 mm (1¼")

401928 – 45 mm (1¾")

Box of 10 pouches (10 individual tail closures per box)

Natura® Flexible Skin Barrier 125259 45 mm (1¾") can be used with Little Ones® pouch 401928.

Little Ones® Two-Piece Urostomy Pouch



The anti-reflux valve keeps urine from pooling at base of stoma. Secure drainage tap helps provide leakproof dependability. Pouch can be easily repositioned to allow for unobstructed drainage, even if the child is asleep. Drainage adapter provides easy attachment to night drainage container.

Transparent, approx. 13 cm (5")

401929 – 32 mm (1¼")

Box of 10 pouches (2 night drainage adapters per box)

Little Ones® Two-Piece Mini-Pouch



Designed for children and/or caregivers who require the convenience and security of a two-piece closed-end system. Eliminates the need for a tail clip and allows for disposal after one-time use.

Opaque, approx. 13 cm (5"), 1-sided comfort panel

401930 – 32 mm (1¼")

Box of 20 mini-pouches

* Little Ones® measuring guides are available in each box of skin barriers.



Adhesive Coupling Technology™ (ACT) Little Ones® Ostomy System*

Adhesive Coupling Technology™ lets you remove the pouch without removing the skin barrier. This makes pouch changing easy, is gentler on a baby's tender skin, and allows you to view the stoma without removing the skin barrier. The soft, colorful flocking and rounded corners of the pouch are comfortable on a baby's sensitive skin.

Adhesive Coupling Technology™ (ACT) Little Ones® Stomahesive® Solid Skin Barrier



The thin, flexible Stomahesive® skin barrier moves with the baby's belly, minimizing discomfort.

Standard, 8.5 cm x 8.5 cm (3³/₈" x 3³/₈")

411641 – For stomas size 4.8 mm to 31 mm (3¹/₁₆" to 1¹/₄")

Box of 5 skin barriers

Adhesive Coupling Technology™ (ACT) Little Ones® Stomahesive® Flexible Skin Barrier with Hydrocolloid Collar Adhesive (HC)



The thin, flexible Stomahesive® skin barrier moves with the baby's belly, minimizing discomfort. The hydrocolloid flexible fabric collar provides an extra level of security and improves flexibility and comfort by gently adhering to sensitive peristomal skin.

Standard, 8.5 cm x 8.5 cm (3³/₈" x 3³/₈")

411642 – For stomas size 4.8 mm to 31 mm (3¹/₁₆" to 1¹/₄")

Box of 5 skin barriers

Adhesive Coupling Technology™ (ACT) Little Ones® Closed-End Pouch with Filter



Standard, approx. 15 cm (6"), Opaque, 2-sided printed comfort panel

411639 – For stomas size 4.8 mm to 31 mm (3¹/₁₆" to 1¹/₄")

Box of 10 pouches

Standard, approx. 15 cm (6"), Transparent, 1-sided printed comfort panel

411640 – For stomas size 4.8 mm to 31 mm (3¹/₁₆" to 1¹/₄")

Box of 10 pouches

Adhesive Coupling Technology™ (ACT) Little Ones® Drainable Pouch with Filter and InvisiClose® (Clipless) Tail Closure



The InvisiClose® outlet facilitates easy drainage, yet provides two levels of closure security you can trust.

Standard, approx. 21 cm (8"), Transparent, 1-sided printed comfort panel

411637 – For stomas up to 31 mm (1¹/₄")

Box of 10 pouches

* Little Ones® measuring guides are available in each box of skin barriers.

OSTOMY SKIN CARE PRODUCTS & ACCESSORIES

Skin Care Products

**Stomahesive® Skin Barrier**

Stomahesive® Skin Barrier is a multi-purpose hydrocolloid skin barrier composed of gelatin, pectin, sodium carboxymethyl-cellulose and polyisobutylene laminated between polyethylene and an opaque paper backing for use as a peristomal covering with ostomies and around fistulas.

21712 – approx. 10 cm x 10 cm (4" x 4")

Box of 5 skin barriers

21715 – approx. 20 cm x 20 cm (8" x 8")

Box of 3 skin barriers

Stomahesive® Protective Powder

Stomahesive® Protective Powder is composed of gelatin, pectin and sodium carboxymethylcellulose. Stomahesive® Protective Powder helps form a protective barrier against excoriating discharge on weeping areas. Increases wear time by absorbing moisture while helping to prevent skin irritation.

025510 – Bottle of 28.3 g (1 oz)

Stomahesive® Paste**Pectin-based Protective Skin Barrier + Filler**

Stomahesive® Paste is composed of ingredients similar to Stomahesive® Protective Powder. Stomahesive® Paste is a highly effective filler, sealant and protective skin barrier for use with colostomy and ileostomy. The paste is also used to protect exposed skin around fistula sites, between the base of the stoma and the opening in the appliance and as a filler for skin folds, uneven skin surfaces and scars.

183910 – Tube of 56.7 g (2 oz)

TIP!

Immediately close the tube after each use to prevent the paste from hardening over time.

AllKare® Adhesive Remover Wipe*

Adhesive Remover Wipe is a pad saturated with a unique iso-alcohol-free adhesive remover. Its formula effectively removes tape, skin barrier, and adhesive residues from the skin that can reduce wear times, and it contains emollients that leave the skin feeling soft and supple.

037436 – Box of 50 wipes

37443 – Box of 100 wipes

AllKare® Protective Barrier Wipe*

Non-water soluble formula for skin protection. Provides a barrier film on the skin under tapes, skin barriers, adhesives, and hydro-colloid dressings to help protect against irritation that may be caused by skin stripping, tension blisters, folliculitis, maceration or adhesive residue.

037439 – Box of 50 wipes

37444 – Box of 100 wipes

* AllKare® wipes do not adversely affect Stomahesive® or Durahesive® skin barrier adhesion or wearing time.

OSTOMY SKIN CARE PRODUCTS & ACCESSORIES

Skin Care & Oral Products



Niltac™ & Silesse™ Sting Free Remover and Skin Barrier



A versatile, full line of adhesive releasers, removers and skin barrier products designed to protect vulnerable skin from the damaging effects of adhesives, body waste and friction.

420788 – Niltac™ Adhesive Remover Wipes

box of 30 wipes

420789 – Silesse™ Protective Skin Barrier Wipes

box of 30 wipes

420787 – Niltac™ Sting Free Adhesive Releaser Spray (50ml)

420790 – Silesse™ Sting Free Protective Skin Barrier Spray (50ml)

Diamonds® Gelling and Odor Control Sachets

Diamonds®



Solidifies contents which helps prevent leaks, provides a flatter more discreet pouch and reduces noise

420791

100 Sachets per jar

Stomahesive® Strip



Double sided adhesive for use as a filler to help protect uneven skin surfaces due to scarring, skin folds or depressions which are adjacent to an ostomy, fistula, or draining wound. Individually packed and sealed.

25542 – 15 mm x 120 mm (approx. 1" x 5")

Box of 15 strips

Orabase® Oral Protective Paste



For temporary relief of minor discomfort and protection of irritated areas of the mouth and gums.

Note: Sold only in minimum quantities of one dozen.

062436 – 7.5 g tube

12-tube package

Orahesive® Denture Adhesive



Powder Form

Holds dentures firmly.

Note: Sold only in minimum quantities of one dozen.

025545 – 25 g bottle

12-bottle package



Eakin Cohesive® Seal*



Eakin Cohesive® seal is available as a donut-shaped ostomy accessory that can be molded to a variety of shapes. It acts as a gasket or barrier and provides enhanced protection for hard-to-fit stomas. The unique moldable hydrocolloid barrier forms a gel that helps prevent skin damage by sealing the area between the skin and ostomy appliance.



Cohesive SLIMS®

Leakage prevention, comfort and security for all types, shapes and sizes of stoma.

Cohesive® SMALL

A THICKER SEAL

Small Cohesive® Seals are the same shape and size as Cohesive SLIMS® but are thicker. They are ideal for extra skin protection, problem skin or people with significant leakage problems.

Cohesive® LARGE

DOUBLE SIZE

Large Cohesive® Seals cover a larger area of skin around the stoma providing protection over a wide area. They are ideal for people with urostomies or for people who have sensitive skin.

Cohesive® StomaWrap™

These seals are ideal for large or oval shaped stomas, for those with limited dexterity and useful in protecting skin around wounds (pictureframing).

839005 – SLIMS™, 5 cm (2")

Box of 10 slim seals

839002 – Small, 5 cm (2")

Box of 20 small seals

839001 – Large, 10 cm (4")

Box of 10 large seals

839006 – StomaWrap™, 8.5 cm (3")

Box of 10 large oval seals

Eakin Cohesive® Skin Barrier*



Eakin Cohesive® Skin Barrier is moldable and moisture-absorbing. It is designed to help protect the skin from contact with body fluids. It is used wherever the skin is damaged, excoriated or vulnerable and needs protection from body wastes and fluids.

839004 – Small, 10 cm x 10 cm (4" x 4")

Box of 5 skin barriers

839003 – Large, 10 cm x 20 cm (4" x 8")

Box of 5 skin barriers

Eakin Cohesive® Paste

Eakin Cohesive® Paste is an alcohol-free, non-sting adhesive paste, used to fill all crevices, skin folds and scars to create an even pouching surface.

NEW

839010 – Tube of 60 g (2.1oz)

1 tube

* Please see page 44 for other Eakin Cohesive® products.



Visi-Flow® Irrigation System

A complete irrigation system offering the advantages of flow visibility and control, easy front-filling and a removable stoma cone.

The Visi-Flow® Irrigation System is only compatible with the Natura® Irrigation Sleeves.

Visi-Flow® Irrigator



Visi-Flow® Irrigator with stoma cone is designed for persons with a colostomy for whom irrigation is indicated.

The Visi-Flow® Irrigator features a front-fill reusable bag with an exclusive “paddle wheel” flow indicator, thumb-control flow regulator, bold-print volume graduations and a flexible stoma cone attachment. Special fitting allows easy removal for cleaning.

2000 mL capacity (67 oz)

401989 – 1 irrigator with stoma cone per boxTaxable

Visi-Flow® Stoma Cone



Stoma cone replacement for Visi-Flow® Irrigator. Visi-Flow® stoma cone is soft on the outside, tough on the inside to prevent premature collapse.

22736 – 1 stoma cone per box

Natura® Visi-Flow® Irrigation Sleeve



Natura® Visi-Flow® Irrigation Sleeves are made of clear film with a slide seal closure for easy entry and security and are designed for use with the Natura® system.

401911 – 38 mm (1½")

401912 – 45 mm (1¾")

401913 – 57 mm (2¼")

401914 – 70 mm (2¾")

Box of 5 irrigation sleeves (1 tail closure per box)

Natura® Visi-Flow® Irrigation Adapter Faceplate



For people who irrigate and prefer a one-piece system. This is to be used with a Visi-Flow® Irrigator and a Natura® Visi-Flow® Irrigation Sleeve. A belt is needed to hold the faceplate in place.

401918 – 45 mm (1¾")

401919 – 70 mm (2¾")

1 irrigation adapter faceplate per box

OSTOMY SKIN CARE PRODUCTS & ACCESSORIES

Accessories

**Night Drainage Container Set**

The urine collection system with the Accuseal® connection. Designed for convenient collection of urine over an extended period.

Set includes:

- 1 Night Drainage Container with tubing, capacity of 2000 mL (67 oz)
- 1 Night Drainage Container cover
- 1 Universal adapter

27060 – 1 set per packageTaxable

Night Drainage Container Tube

To be used with night drainage container.

27062 – 1 tube, length of 147 cm (approx. 58") Taxable

Accuseal® Adapter*

Adapter that attaches to the Accuseal® tap on the urostomy pouch, if needed to connect to a drainage system.

21915 – Package of 2 adapters

Ostomy Appliance Belt


Adjustable belt designed for use with the Natura® and Active Life® Systems. Ostomy Appliance Belt made of soft, high quality, white washable material; provides added sense of security.

175507 – Suitable for sizes up to approximately 107 cm (42")

1 belt per boxTaxable

TIP!

The belt should be hand-washed using a mild soap and laid flat to dry.

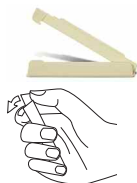ConvaTec Tail Closure

The ConvaTec Tail Closure has been designed with a soft comfortable feel and greater flexibility. For improved security, a clicking noise can be heard when closing the tail closure. Designed for use with drainable pouches. Prior to initial use, open and close a few times.

175635 – Box of 10 tail closuresTaxable

Did you know?

Closing tail closures by pressing only in the middle may not lock the closure properly. You must also press down on the "head".

DuoLock® Curved Tail Closure

DuoLock® Curved Tail Closure for use with drainable pouches. This unique curved tail closure has a dual locking mechanism to prevent accidental opening and conforms to body contours.

175652 – Box of 10 tail closuresTaxable

Natura® Visi-Flow® Irrigation Sleeve Tail Closure

Designed with a hook which allows Natura® Irrigation Sleeve to be folded upward and hung from flange for more convenient wearing, while waiting for second evacuation to occur.

175650 – Box of 10 tail closuresTaxable

* Please contact the ConvaTec Customer Relations Center at 1 800 465-6302 to order the Accuseal® adapter.

POSTOPERATIVE SYSTEMS & KITS

Postoperative Kits



Natura® Postoperative Kit, Non-Sterile*



Each kit contains:

- 1 Standard, Drainable, transparent pouch
- 1 Stomahesive® Flexible skin barrier
- 1 Tail closure

Box of 5 kits

401920 – For colostomy / ileostomy, 45 mm (1³/₄")

401921 – For colostomy / ileostomy, 70 mm (2³/₄")

Each kit contains:

- 1 Standard, Drainable, transparent pouch
- 1 Stomahesive® Solid skin barrier
- 1 Tail closure

Box of 5 kits

401922 – For colostomy / ileostomy, 100 mm (4")



Each kit contains:

- 1 Urostomy pouch with Accuseal® tap
- 1 Duraheasive® Flexible skin barrier
- 1 Accuseal® adapter

Box of 5 kits

401923 – For urostomy, 45 mm (1³/₄")

Natura® Postoperative Kit, Non-Sterile*

Each kit contains:

- 1 standard size, drainable, transparent pouch with comfort panels on the bodyside and InvisiClose® tail closure with Lock-It Pocket™
- 1 skin barrier
- 1 Low Pressure Adaptor

Box of 5 kits

Duraheasive® with ConvaTec Moldable Technology™

416917 – For colostomy / ileostomy, 45 mm (1³/₄"), for stoma sizes 13 mm - 22 mm (1/2" - 7/8")

416921 – For colostomy / ileostomy, 57 mm (2¹/₄"), for stoma sizes 22 mm - 33 mm (7/8" - 1¹/₄")

416925 – For colostomy / ileostomy, 70 mm (2³/₄"), for stoma sizes 33 mm - 45 mm (1¹/₄" - 1³/₄")

Stomahesive® with ConvaTec Moldable Technology™

416918 – For colostomy / ileostomy, 45 mm (1³/₄"), for stoma sizes 13 mm - 22 mm (1/2" - 7/8")

416922 – For colostomy / ileostomy, 57 mm (2¹/₄"), for stoma sizes 22 mm - 33 mm (7/8" - 1¹/₄")

416926 – For colostomy / ileostomy, 70 mm (2³/₄"), for stoma sizes 33 mm - 45 mm (1¹/₄" - 1³/₄")

Stomahesive® Cut-to-Fit

416919 – For colostomy / ileostomy, 45 mm (1³/₄")

416923 – For colostomy / ileostomy, 57 mm (2¹/₄")



* Postoperative Kits are sold to hospitals only.



Natura® Postoperative Kit, Non-Sterile* including Moldable Skin Barrier with ACCORDION flange



Each kit contains:

- 1 standard size, drainable, transparent pouch with comfort panels on the bodyside and InvisiClose® tail closure with Lock-It Pocket™
- 1 Moldable ACCORDION Skin Barrier

Box of 5 kits

Durahesive® with ConvaTec Moldable Technology™

421045 – For colostomy / ileostomy, 57 mm, for stoma sizes 13 mm - 22 mm ($\frac{1}{2}$ " - $\frac{7}{8}$ ")

421046 – For colostomy / ileostomy, 57 mm, for stoma sizes 22 mm - 33 mm ($\frac{7}{8}$ " - $1\frac{1}{4}$ ")

421047 – For colostomy / ileostomy, 70 mm, for stoma sizes 33 mm - 45 mm ($1\frac{1}{4}$ " - $1\frac{3}{4}$ ")

Stomahesive® with ConvaTec Moldable Technology™

421042 – For colostomy / ileostomy, 57 mm, for stoma sizes 13 mm - 22 mm ($\frac{1}{2}$ " - $\frac{7}{8}$ ")

421043 – For colostomy / ileostomy, 57 mm, for stoma sizes 22 mm - 33 mm ($\frac{7}{8}$ " - $1\frac{1}{4}$ ")

421044 – For colostomy / ileostomy, 70 mm, for stoma sizes 33 mm - 45 mm ($1\frac{1}{4}$ " - $1\frac{3}{4}$ ")

Natura® Postoperative Kit, Sterile*

Stomahesive® Cut-to-Fit

Each kit contains:

- 1 skin barrier with cut-to-fit opening
- 1 Drainable Pouch (12" transparent with 1-sided comfort panel) with InvisiClose® Tail Closure System
- 1 Low-Pressure Adaptor



Box of 5 kits

416920 – For colostomy / ileostomy, 45 mm ($1\frac{3}{4}$ ")

416924 – For colostomy / ileostomy, 57 mm ($2\frac{1}{4}$ ")

Stomahesive® Cut-to-Fit for urostomy

Each kit contains:

- 1 skin barrier with cut-to-fit opening
- 1 Urostomy Pouch (10" transparent with 1-sided comfort panel) with Accuseal® Swivel Tap
- 1 Accuseal® night drainage adaptor
- 1 Low-Pressure Adaptor

Box of 5 kits

416938 – For urostomy, 45 mm ($1\frac{3}{4}$ ")

416946 – For urostomy, 57 mm ($2\frac{1}{4}$ ")

416954 – For urostomy, 70 mm ($2\frac{3}{4}$ ")

* Postoperative Kits are sold to hospitals only.



Natura® Loop Ostomy Postoperative Kit*

- Each kit contains:
- 1 standard size, Drainable transparent pouch with comfort panels on the bodyside and InvisiClose® tail closure with Lock-It Pocket™
 - 1 solid skin barrier
 - 1 Low Pressure Adaptor
 - 1 sterile Loop Rod

Non-Sterile Kits

Box of 5 kits



Stomahesive® Cut-to-Fit

- 416927** – 70 mm (2¾")
416929 – 100 mm (4")

Sterile Kits

Box of 5 kits

Stomahesive® Cut-to-Fit

- 416928** – 70 mm (2¾")
416930 – 100 mm (4")

Esteem synergy® Professional Kit, Non-Sterile



- Each kit contains:
- 5 Drainable transparent pouches
 - 5 Stomahesive® skin barriers with flexible collar
 - 1 Tail closure

- 409275** – For colostomy / ileostomy, 89 mm (3½")

Esteem® Stomahesive® Post-Operative Kits with drainable pouches, InvisiClose®, Lock-It Pocket™



- Each kit contains:
- 1 Cut-To-Fit Skin Barrier, standard size pouch with 1-sided comfort panel and InvisiClose® Tail Closure

Box of 5 kits

Transparent, 1-sided comfort panel, Standard Size (12"), non-sterile

- 416915** – 100 mm (4")

Transparent, 1-sided comfort panel, Standard Size (12"), sterile

- 416916** – 100 mm (4")

ConvaTec Loop Rod, Sterile



- T-shape loop rods help keep the rod in place
- Swivel head design allows for easy application and removal
- Designed with holes, gives the option of suturing

22355 – 65mm Loop Ostomy Rod - Sterile

- Box of 10

22356 – 90mm Loop Ostomy Rod - Sterile

- Box of 10

* Postoperative Kits are sold to hospitals only.



Eakin Cohesive® Fistula & Wound Pouch



Eakin Wound Pouches provide a comfortable, skin friendly, cost effective and time saving procedure for managing fistulas, multiple or irregular stomas and high-output wounds.

839250 – Small pouch with “Fold & Tuck” closure for wounds up to 45 x 30mm

Box of 10 pouches

839251 – Medium pouch with “Fold & Tuck” closure for wounds up to 110 x 75mm

Box of 10 pouches

839252 – Large pouch with “Fold & Tuck” closure for wounds up to 175 x 110mm

Box of 10 pouches

839253 – Extra large pouch with “Fold & Tuck” closure for horizontal wounds up to 245 x 160mm

Box of 5 pouches

839260 – Small pouch with remote drainage attachment for wounds up to 45 x 30mm

Box of 10 pouches

839267 – Small Plus pouch with remote drainage attachment for wounds up to 60 x 86mm

Box of 10 pouches

839261 – Medium pouch with remote drainage attachment for wounds up to 110 x 75mm

Box of 10 pouches

839262 – Large pouch with remote drainage attachment for wounds up to 175 x 110mm

Box of 10 pouches

839263 – Extra large pouch with remote drainage attachment for horizontal wounds up to 245 x 160mm

Box of 5 pouches

839264 – Extra large pouch with remote drainage attachment and Access Window for horizontal wounds up to 245 x 160mm

Box of 5 pouches

839265 – Extra large pouch with remote drainage attachment for vertical wounds up to 245 x 160mm

Box of 5 pouches

839268 – Extra large pouch with remote drainage attachment and Access Window for vertical wounds up to 245 x 160mm

Box of 5 pouches

839266 – Extra large pouch with remote drainage attachment for vertical incision wounds up to 290 x 130mm

Box of 5 pouches

839280 – Eakin Access Windows for use with Eakin Wound Pouches

Box of 5 windows



Flexi-Seal® Fecal Management System (FMS)

Dramatically advancing fecal incontinence management.

Flexi-Seal® FMS is a temporary containment device indicated for the fecal management of patients with little or no bowel control and liquid or semi-liquid stool. It is designed to safely and effectively contain and divert fecal matter, protect patients' wounds from fecal contamination and reduce both the risk of skin breakdown and spread of infection.

Flexi-Seal® Fecal Management System Kit contains 1 soft silicone catheter tube assembly, 1 syringe and collection bag(s).



Flexi-Seal[®] PROTECT



Flexi-Seal[®] SIGNAL[™]
The seal of innovation

421630 – Flexi-Seal® PROTECT FMS Kit (1 kit/box, 1bag)

1 kit per box, 1 bag

Flexi-Seal® PROTECT is designed to protect patients and clinicians against over-inflation risk. Flexi-Seal® PROTECT is the only fecal management system that senses when the retention balloon is correctly filled and visually indicates over inflation with color indicators. Odor barrier runs the full length of the catheter, providing an end-to-end odor barrier.

418000 – Flexi-Seal® SIGNAL™ FMS

1 kit per box, with 3 filtered bags

411100 – Flexi-Seal® FMS

1 Kit

411101 – Disposable Replacement Collection Bags, 1000 ml capacity

Box of 10 bags

411102 – Disposable Replacement Collection Bags with Filters, 1000 ml capacity

Box of 10 bags

411108 – Disposable Replacement Collection Bags with APS Filter, 1000 ml capacity

Box of 10 bags

SECUREMENT DEVICES

FLEXI-TRAK® Anchoring Device




For anchoring tubes and catheters to skin, such as: urinary catheters, gastrostomy tubes, I.V. tubes, rectal tubes, chest tubes, central venous lines, pacemaker wires, sump tubes, biliary drainage tubes, nephrostomy tubes, and peritoneal dialysis tubes.

37449 – Large, 38 mm x 100 mm (1.5" x 4")

Box of 50 anchorsTaxable

Drain-Fix™ Securement Device

Securement device for percutaneous catheters and tubes

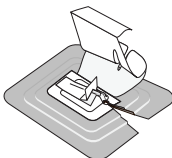

Drain-Fix™ Securement Device is a sterile securement device for percutaneous catheters and tubes from sizes CH 5 to CH 22. Drain-Fix™ Securement Device is applied for fixating for instance Pigtail catheters, Balloon catheters, Suprapubic catheters, PEG feeding tubes and other percutaneous catheters and drains.

680M – Small, 5 CH - 14 CH

Box of 25 pcs

685M – Large, 12 CH - 22 CH

Box of 20 pcs



GentleCath™ Intermittent Urinary Catheter™



GentleCath™ catheter is designed for maximum comfort, safety and ease of use. The tip is soft and rounded and the eyelets are polished and smooth for easy insertion.



Female PVC Urinary Catheter

- 501019 – Size 08 FR
- 501020 – Size 10 FR
- 501021 – Size 12 FR
- 501022 – Size 14 FR
- 501023 – Size 16 FR

Box of 100



Male Straight Tip PVC Urinary Catheter

- 501002 – Size 10 FR
- 501003 – Size 12 FR
- 501004 – Size 14 FR
- 501005 – Size 16 FR
- 501006 – Size 18 FR

Box of 100



Male Coudé Tip PVC Urinary Catheter

- 501011 – Size 08 FR
- 501012 – Size 10 FR
- 501013 – Size 12 FR
- 501014 – Size 14 FR
- 501015 – Size 16 FR
- 501016 – Size 18 FR

Box of 100

Feel the difference!



GentleCath™ is designed to make catheter insertion and removal as smooth and easy as possible for the patient.

- Polished and rounded eyelets minimise friction.
- Soft, rounded tip aids smooth insertion.
- Available in male and female straight tip and male coudé tip, for insertion through complicated passages.



See the difference!

- Color coded connector for easy FR size identification.
- Available in a wide range of sizes to meet your patients' needs.
- Fits drainage bags or other devices.





AQUACEL® Foam Hydrofiber® Dressing, Sterile



Non-Adhesive



Adhesive



Adhesive Heel



Adhesive Sacral



AQUACEL®
Foam

AQUACEL® Foam Dressing offers a range of sterile, adhesive and non-adhesive absorbent foam wound dressings consisting of 4 different components: 1) a waterproof and bacterial/viral barrier polyurethane top layer which protects the wound from external contaminants, reducing the risk of infection. The top layer also helps to manage the moisture vapour transmission of the exudate absorbed by the dressing; 2) a soft absorbent foam pad; 3) an AQUACEL® wound interface made of Hydrofiber® Technology (sodium carboxymethylcellulose) which creates a soft, cohesive gel that intimately micro-contours to the wound bed, maintains a moist environment in the wound to support the body's healing process and aids in the removal of non-viable tissue from the wound (autolytic debridement) without damaging newly formed tissue; 4) and a gentle silicone adhesive border for the adhesive dressing range which provides secure, skin-friendly adhesion to the surrounding skin and not on the wound bed thereby supporting non-traumatic dressing removal.

AQUACEL® Foam dressings may be used as primary or secondary dressing for the management of both chronic and acute wounds, such as leg ulcers, pressure ulcers (Stages II-IV) and diabetic ulcers; surgical wounds (post-operative, donor sites, dermatological); partial thickness (second degree) burns; traumatic or surgical wounds left to heal by secondary intention such as dehiscent surgical incisions; surgical wounds that heal by primary intent, such as dermatological and surgical incisions (e.g. orthopaedic and vascular); etc (see package insert for complete instructions for use).

Non-Adhesive:

420631 – Box of 10 dressings, approx. 5 cm x 5 cm (approx. 2" x 2")	
420633 – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")	
420635 – Box of 5 dressings, approx. 15 cm x 15 cm (approx. 6" x 6")	
420637 – Box of 5 dressings, approx. 15 cm x 20 cm (approx. 6" x 8")	
420636 – Box of 5 dressings, approx. 20 cm x 20 cm (approx. 8" x 8")	
421156 – Box of 5 dressings, approx. 10 cm x 20 cm (approx. 4" x 8")	

Adhesive:

420804 – Box of 10 dressings, approx. 8 cm x 8 cm (approx. 3.2" x 3.2")	
420680 – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")	
420619 – Box of 10 dressings, approx. 12.5 cm x 12.5 cm (approx. 5" x 5")	
420621 – Box of 10 dressings, approx. 17.5 cm x 17.5 cm (approx. 7" x 7")	
420623 – Box of 5 dressings, approx. 21 cm x 21 cm (approx. 8.5" x 8.5")	
420624 – Box of 5 dressings, approx. 25 cm x 30 cm (approx. 10" x 12")	
421150 – Box of 5 dressings, approx. 10 cm X 20 cm (approx. 4" x 8")	
421152 – Box of 5 dressings, approx. 10 cm x 25 cm (approx. 4" x 10")	
421154 – Box of 5 dressings, approx. 10 cm x 30 cm (approx. 4" x 12")	

Adhesive Heel:

420625 – Box of 5 dressings, approx. 19.8 cm x 14 cm (approx. 8" x 5.5")	
---	-------

Adhesive Sacral:

420626 – Box of 5 dressings, approx. 20 cm x 16.9 cm (approx. 8" x 7")	
---	-------

AQUACEL®
Foam
Skin's Guardian



**Protecting.
Defending.
Nurturing.**

*Additional tape, or other means of fixation, will be required to secure the AQUACEL® Foam non-adhesive dressing in place.



AQUACEL® Ag Foam Hydrofiber® Dressing, Sterile



The only silver foam dressing that offers the healing benefits of an AQUACEL® contact layer, the comfort of foam and the bacteria-killing power of ionic silver. Choose from several AQUACEL® Ag Foam dressing shapes and sizes, or cut them to a custom shape or size,* keeping all the benefits of Hydrofiber® Technology.

Non-Adhesive



Non-Adhesive:

- 420639 – Box of 10 dressings, approx. 5 cm x 5 cm (approx. 2" x 2")
- 420642 – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")
- 420645 – Box of 5 dressings, approx. 15 cm x 15 cm (approx. 6" x 6")
- 420646 – Box of 5 dressings, approx. 20 cm x 20 cm (approx. 8" x 8")
- 420806 – Box of 5 dressings, approx. 15 cm x 20 cm (approx. 6" x 8")

Adhesive



Adhesive:

- 420805 – Box of 10 dressings, approx. 8 cm x 8 cm (approx. 3.2" x 3.2")
- 420681 – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")
- 420627 – Box of 10 dressings, approx. 12.5 cm x 12.5 cm (approx. 5" x 5")
- 420628 – Box of 10 dressings, approx. 17.5 cm x 17.5 cm (approx. 7" x 7")
- 420629 – Box of 5 dressings, approx. 21 cm x 21 cm (approx. 8.5" x 8.5")
- 420807 – Box of 5 dressings, approx. 25 cm x 30 cm (approx. 10" x 12")

Adhesive Heel



Adhesive Heel:

- 420647 – Box of 5 dressings, approx. 19.8 cm x 14 cm (approx. 8" x 5.5")

Adhesive Sacral



Adhesive Sacral:

- 420648 – Box of 5 dressings, approx. 20 cm x 16.9 cm (approx. 8" x 7")

FoamLite™ ConvaTec Dressing




Suitable for non-exuding chronic and acute wounds. FoamLite™ ConvaTec is a thin foam dressing with a breathable film and perforated gentle skin-friendly silicone adhesive. Thin and flexible, it conforms to the body even in awkward positions.



- 421557 – Box of 10 dressings, approx. 8 cm x 8 cm (approx. 3.1" x 3.1")
- 421559 – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")
- 421561 – Box of 10 dressings, approx. 15 cm x 15 cm (approx. 5.9" x 5.9")
- 421563 – Box of 10 dressings, approx. 5.5 cm x 12 cm (approx. 2.1" x 4.7")
- 421927 – Box of 10 dressings, approx. 5 cm x 5 cm (approx. 2" x 2")
- 421928 – Box of 10 dressings, approx. 10 cm x 20 cm (approx. 4" x 8")



Dressing Sizes

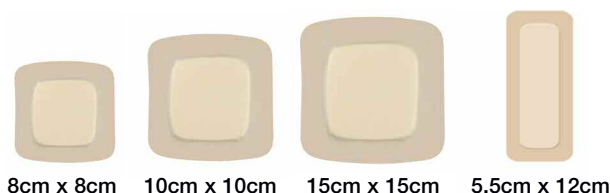

EXUDATE MANAGEMENT

Absorbent Cover Dressings



AQUACEL® Surgical Cover Dressing, Sterile



AQUACEL® Surgical cover dressing is a sterile post operative dressing comprised of an inner (wound contact) non woven pad incorporating Hydrofiber® Technology. The pad is composed of sodium carboxymethylcellulose stitch bonded with nylon and elastane yarns. The pad is held in place between a top layer of hydrocolloid which is bound to an outer polyurethane film and a windowed skin contact layer consisting of polyurethane film sandwiched between two further layers of hydrocolloid.

AQUACEL® Surgical cover dressing is a waterproof, bacterial and viral barrier dressing which is designed to conform to changes in the wound/incision geometry during body movement and post operative edema.

AQUACEL®
Surgical

412017 – Box of 10 dressings, approx. 9 cm x 10 cm (approx. 3.5" x 4")	
412018 – Box of 10 dressings, approx. 9 cm x 15 cm (approx. 3.5" x 6")	
412019 – Box of 10 dressings, approx. 9 cm x 25 cm (approx. 3.5" x 10")	
420669 – Box of 10 dressings, approx. 9 cm x 30 cm (approx. 3.5" x 12")	
412020 – Box of 10 dressings, approx. 9 cm x 35 cm (approx. 3.5" x 14")	

AQUACEL® Ag Surgical Cover Dressing, Sterile



AQUACEL® Ag Surgical cover dressing is a sterile post operative dressing comprised of an inner (wound contact) non woven pad incorporating Hydrofiber® Technology and ionic silver stitch bonded with nylon and elastane yarns. The pad is held in place between a top layer of hydrocolloid which is bound to an outer polyurethane film and a windowed skin contact layer consisting of polyurethane film sandwiched between two further layers of hydrocolloid.

AQUACEL® Ag Surgical cover dressing is a waterproof, bacterial and viral barrier dressing which is designed to conform to changes in the wound/incision geometry during body movement and post operative edema.

AQUACEL® Ag
Surgical

412009 – Box of 10 dressings, approx. 9 cm x 10 cm (approx. 3.5" x 4")	
412010 – Box of 10 dressings, approx. 9 cm x 15 cm (approx. 3.5" x 6")	
412011 – Box of 10 dressings, approx. 9 cm x 25 cm (approx. 3.5" x 10")	
420670 – Box of 10 dressings, approx. 9 cm x 30 cm (approx. 3.5" x 12")	
412012 – Box of 10 dressings, approx. 9 cm x 35 cm (approx. 3.5" x 14")	

CarboFlex® Odor Control Dressing, Sterile

CarboFlex® Odor Control Dressing is a sterile non-adhesive dressing with an absorbent wound contact layer (containing alginate and Hydrofiber® Technology), an activated charcoal central pad and a smooth water resistant top layer. It is indicated for the management of malodorous acute and chronic wounds.

CarboFlex®
Odor Control Dressing

403202 – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")	
403203 – Box of 5 dressings, approx. 7.5 cm x 15 cm (approx. 3" x 6")	
403204 – Box of 5 dressings, approx. 15 cm x 20 cm (approx. 6" x 8")	



AQUACEL® Ag+ Extra™, Sterile



AQUACEL® Ag+ EXTRA™ Enhanced Hydrofiber™ Dressing with Silver and Strengthening Fibre is a soft, sterile, non-woven dressing made from two layers of sodium carboxymethylcellulose impregnated with 1.2% ionic silver (an antimicrobial agent), enhanced by ethylenediaminetetra-acetic acid di-sodium salt (EDTA) and benzethonium chloride (BeCl), and strengthened by regenerated cellulose fibre. This dressing absorbs high amounts of wound fluid and bacteria and creates a soft, cohesive gel that intimately conforms to the wound surface, maintains a moist environment and aids in the removal of non-viable tissue from the wound (autolytic debridement). A moist wound environment and control of wound bacteria supports the body's healing process and helps reduce the risk of wound infection. The ionic silver in the dressing kills pathogenic microorganisms, both planktonic and within bacterial biofilms, including wound bacteria, yeasts and moulds. The dressing also disrupts and absorbs biofilm, prevents biofilm formation/reformation and increases the efficiency of silver transfer to microorganisms. The dressing itself also provides an antimicrobial barrier to protect the wound bed.



- 413566** – Box of 10 dressings, approx. 5 cm x 5 cm (approx. 2" x 2")
- 413567** – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")
- 413568** – Box of 5 dressings, approx. 15 cm x 15 cm (approx. 6" x 6")
- 413569** – Box of 5 dressings, approx. 20 cm x 30 cm (approx. 8" x 12")


AQUACEL® Ag+ Ribbon Dressing, Sterile



AQUACEL® Ag+ Enhanced Hydrofiber™ Ribbon Dressing with Silver and Strengthening Fibre is a soft, sterile, non-woven ribbon dressing made from a layer of sodium carboxymethylcellulose impregnated with 1.2% ionic silver (an antimicrobial agent), enhanced by ethylenediaminetetraacetic acid di-sodium salt (EDTA) and benzethonium chloride (BeCl), and strengthened by regenerated cellulose fibre. This dressing absorbs high amounts of wound fluid and bacteria and creates a soft, cohesive gel that intimately conforms to the wound surface, maintains a moist environment and aids in the removal of non-viable tissue from the wound (autolytic debridement). A moist wound environment and control of wound bacteria supports the body's healing process and helps reduce the risk of wound infection. The ionic silver in the dressing kills pathogenic microorganisms, both planktonic and within bacterial biofilms, including wound bacteria, yeasts and moulds. The dressing also disrupts and absorbs biofilm, prevents biofilm formation/reformation and increases the efficiency of silver transfer to microorganisms. The dressing itself also provides an antimicrobial barrier to protect the wound bed.



- 413570** – Box of 5 dressings, approx. 1 cm x 45 cm (approx. 0.39" x 18")
- 413571** – Box of 5 dressings, approx. 2 cm x 45 cm (approx. 3/4" x 18")

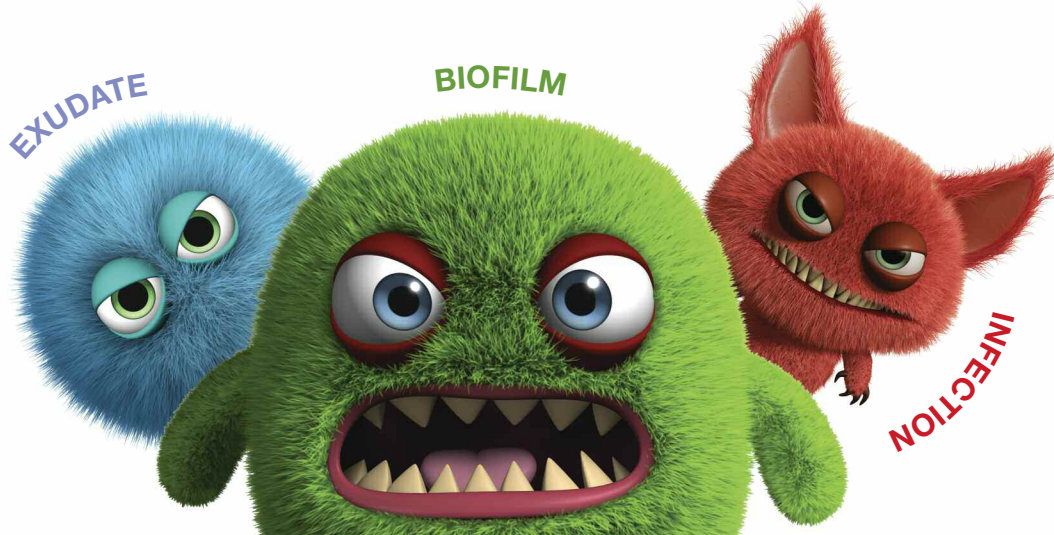



AQUACEL® Hydrofiber® Wound Dressing, Sterile



AQUACEL® Hydrofiber® Dressing is a sterile, white, fibrous dressing derived from 100% pure sodium carboxymethylcellulose (NaCMC). AQUACEL® dressing combines the look, feel and handling properties of gauze and leading alginates, but retains more fluid than alginates and almost three times more fluid than gauze when used under compression. AQUACEL® dressing is indicated for the management of abrasions, lacerations, and under the supervision of a healthcare professional, for pressure ulcers, leg ulcers, incisions, donor sites and second degree burns as well as on infected wounds and for use in managing surgical or traumatic wounds that have been left to heal by secondary intention. Additionally, it can be used on wounds prone to bleeding and to facilitate the control of minor bleeding.

AQUACEL®
Dressings

TIP!

AQUACEL® & AQUACEL® Ag dressings can be used under DuoDERM® Hydrocolloid dressings.

- 177901** – Box of 10 dressings, approx. 5 cm x 5 cm (approx. 2" x 2")
- 177902** – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")
- 177903** – Box of 5 dressings, approx. 15 cm x 15 cm (approx. 6" x 6")

AQUACEL® Hydrofiber® Ribbon Dressing with Strengthening Fiber, Sterile



AQUACEL® Dressings with Strengthening Fiber have all the advantages of the original AQUACEL® Ribbon Dressings and now, even more strength. Indeed, stitch bonding has been added to strengthen the original AQUACEL® Ribbon Dressings. By doing so, the dressings will:

- Be over 20 times stronger*
- Have more than 1.3 times less shrinkage*

- 403770** – Box of 5 dressings, approx. 2 cm x 45 cm (approx. 3/4" x 18")
- 420127** – Box of 5 dressings, approx. 1 cm x 45 cm (approx. 0.39" x 18")

AQUACEL® Ag Hydrofiber® Wound Dressing with Silver, Sterile



AQUACEL® Ag dressing is an absorbent **antimicrobial** dressing that features all of the benefits of AQUACEL® dressing with the important addition of ionic silver. AQUACEL® Ag dressing kills a broad spectrum of wound pathogens (microbicidal) in the dressing that can cause infection including Methicillin Resistant Staphylococcus Aureus (MRSA) and Vancomycin Resistant Enterococcus (VRE)¹. AQUACEL® Ag dressing's antimicrobial is conveniently available within the exudate management dressing eliminating the need for the application of both an antimicrobial product to control bacteria and a dressing to manage exudate.

AQUACEL® Ag
Dressings

Did you know?

You can use AQUACEL® Ag dressing in a wound that has been cleansed with either sterile water or sterile normal saline.

- 403706** – Box of 10 dressings, approx. 5 cm X 5 cm (approx. 2" x 2")
- 403708** – Box of 10 dressings, approx. 10 cm X 10 cm (approx. 4" x 4")
- 403710** – Box of 5 dressings, approx. 15 cm X 15 (approx. 6" x 6")
- 403711** – Box of 5 dressings, approx. 20 cm X 30 cm (approx. 8" x 12")

AQUACEL® Ag Hydrofiber® Ribbon Dressing with Strengthening Fiber, Sterile



AQUACEL® Ag Dressings with Strengthening Fiber have all the advantages of the original AQUACEL® Ag Ribbon Dressings and now, even more strength. Indeed, stitch bonding has been added to strengthen the original AQUACEL® Ag Ribbon Dressings. By doing so, the dressings will:

- Be over 20 times stronger*
- Have more than 1.3 times less shrinkage*

- 403771** – Box of 5 dressings, approx. 2 cm x 45 cm (approx. 3/4" x 18")
- 420128** – Box of 5 dressings, approx. 1 cm x 45 cm (approx. 0.39" x 18")

¹ Jones SA et al., Bowler PG, Walker M, Parsons D. Controlling wound bioburden with a novel silver-containing Hydrofiber® dressing. *Wound Repair Regeneration*. May-June 2004; 12(3):288-294.

* Compared with the original AQUACEL® and AQUACEL® Ag dressings.



AQUACEL® Extra™ Hydrofiber® Dressing with Strengthening Fiber, Sterile



AQUACEL® Extra™ Dressing with Strengthening Fiber is a soft, sterile, non-woven dressing made from two layers of Hydrofiber® (sodium carboxymethylcellulose) dressing stitched together with strengthening fibers. This conformable and highly absorbent dressing absorbs wound fluid and creates a soft cohesive gel that intimately micro-contours to the wound bed, maintains a moist environment which supports the body's healing process and aids in the removal of necrotic material from the wound (autolytic debridement) without damaging newly formed tissue. AQUACEL® EXTRA™ dressing is designed to make it 9X stronger and provides an increased absorbency of 39% higher than current AQUACEL® dressing.¹



- 420671** – Box of 10 dressings, approx. 5 cm x 5 cm (approx. 2" x 2")
- 420672** – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")
- 420673** – Box of 5 dressings, approx. 15 cm x 15 cm (approx. 6" x 6")

AQUACEL® Ag Extra™ Hydrofiber® Dressing with Strengthening Fiber, Sterile



AQUACEL® Ag Extra™ Dressing with Silver and Strengthening Fiber is a soft, sterile, non-woven dressing made from two layers of ionic silver impregnated Hydrofiber® (sodium carboxymethylcellulose) which are stitched together with strengthening fibers. This conformable and highly absorbent dressing absorbs wound fluid and creates a soft cohesive gel that intimately micro-contours to the wound bed, maintains a moist environment which supports the body's healing process and aids in the removal of necrotic material from the wound (autolytic debridement) without damaging newly formed tissue. The silver in the dressing kills wound bacteria held in the dressing, and provides an antimicrobial barrier to protect the wound bed. AQUACEL® Ag Extra™ dressing is designed to make it 9X stronger and provides an increased absorbency of 50% higher than current AQUACEL® Ag dressing.²



- 420675** – Box of 10 dressings, approx. 5 cm x 5 cm (approx. 2" x 2")
- 420676** – Box of 10 dressings, approx. 10 cm x 10 cm (approx. 4" x 4")
- 420678** – Box of 5 dressings, approx. 15 cm x 15 cm (approx. 6" x 6")
- 420679** – Box of 5 dressings, approx. 20 cm x 30 cm (approx. 8" x 12")

1. Preliminary assessment of the physical properties of AQUACEL® EXTRA and AQUACEL® Scientific Background Report. WHRI3461 TA214. 2011, Data on File, ConvaTec Inc.
2. Assessment of the in vitro physical properties of AQUACEL® Ag Extra™ and AQUACEL® Ag dressings. Technical Assessment. WHRI3602 TA235. 2012, Data on File, ConvaTec Inc.

EXUDATE MANAGEMENT

Absorbent Filler & Antimicrobial Dressings



AQUACEL® Ag Burn Dressing, Sterile



The AQUACEL® Ag Burn Hydrofiber® Dressing reinforced with nylon, silver-impregnated, antimicrobial, absorbent dressing is a soft, sterile, non-woven hydroentangled dressing comprised of Hydrofiber® (sodium carboxymethylcellulose) with nylon thread used to stitch bond the product along the length. The dressing contains 1.2% w/w ionic silver. The silver in the dressing kills a broad spectrum of wound bacteria. This dressing absorbs high amounts of wound fluid and bacteria and creates a soft, cohesive gel that intimately conforms to the wound surface, maintains a moist environment and aids in the removal of non-viable tissue from the wound (autolytic debridement).

Under the supervision of a healthcare provider, the AQUACEL® Ag Burn Hydrofiber® Dressing is indicated mainly for the management of partial thickness (second degree) burns and wounds that are prone to bleeding, such as wounds that have been mechanically or surgically debrided and donor sites. Additionally, it can also be used for other indications such as diabetic foot ulcers, leg ulcers (venous stasis ulcers, arterial ulcers and leg ulcers of mixed etiology), pressure ulcers/sores (partial & full thickness), surgical wounds left to heal by secondary intention such as dehiscent surgical incisions, surgical wounds that heal by primary intention such as dermatological and surgical incisions (e.g., orthopedic and vascular), traumatic wounds, oncology wounds with exudate (e.g. fungoides-cutaneous, tumors, fungating carcinoma, cutaneous metastasis, Kaposi's sarcoma, and angiosarcoma), infected or painful wounds.

AQUACEL® Ag.
Burn



403786	– Box of 5 dressings, approx. 13 cm x 10 cm (approx. 5" x 4")	
403787	– Box of 5 dressings, approx. 17 cm x 15 cm (approx. 7" x 6")	
403788	– Box of 5 dressings, approx. 23 cm x 30 cm (approx. 9" x 12")	
403789	– Box of 3 dressings, approx. 23 cm x 100 cm (approx. 9" x 40")	
403791	– Box of 1 dressing, Size 1 (Glove)	
403792	– Box of 1 dressing, Size 2 (Glove)	
403793	– Box of 1 dressing, Size 3 (Glove)	
403794	– Box of 1 dressing, Size 4 (Glove)	
403795	– Box of 1 dressing, Size 5 (Glove)	

KALTOSTAT® Calcium Alginate Wound Dressing, Sterile



KALTOSTAT® wound dressing is a soft, sterile non-woven pad or ribbon dressing of calcium-sodium alginate fibre. The alginate fibres absorb wound fluid and form a moist environment at the wound interface. It may be used for minor abrasions, lacerations, scalds and burns, and under the supervision of a healthcare professional, for stage II to IV pressure ulcers, venous ulcers, arterial and diabetic ulcers, donor sites, surgical incisions, and traumatic wounds. Upon contact with a bleeding wound, KALTOSTAT® dressing promotes hemostasis.

KALTOSTAT®

168117	– Box of 5 dressings, 2g packing	
168210	– Box of 10 dressings, approx. 5 cm x 5 cm (approx. 2" x 2")	
168212	– Box of 10 dressings, approx. 7 cm x 12 cm (approx. 3" x 5")	
168214	– Box of 10 dressings, approx. 10 cm x 19 cm (approx. 4" x 8")	




DuoDERM Signal® A Tapered-Edge CGF® (Control Gel Formula) Hydrocolloid Dressing, Sterile



DuoDERM Signal® dressing is a tapered-edge adhesive hydrocolloid dressing. It is flexible, secure, translucent and has a green indicator line on the dressing which helps to determine when to change the dressing. DuoDERM Signal® dressing is ideal for use on wounds with moderate amounts of exudate.

Multi-Purpose Use:

	403326 – Square - approx. 10 cm x 10 cm (approx. 4" x 4") Box of 5 dressings, individually wrapped
	403327 – Square - approx. 14 cm x 14 cm (approx. 6" x 6") Box of 5 dressings, individually wrapped
	403328 – Square - approx. 20 cm x 20 cm (approx. 8" x 8") Box of 5 dressings, individually wrapped
	403332 – Triangle - approx. 15 cm x 18 cm (approx. 6" x 7") Box of 5 dressings, individually wrapped
	403333 – Triangle - approx. 20 cm x 23 cm (approx. 8" x 9") Box of 5 dressings, individually wrapped
Heel	410500 – Heel shape - approx. 18 cm x 19 cm (approx. 7" x 8") Box of 5 dressings, individually wrapped
Sacrum	410501 – Butterfly shape - approx. 20 cm x 22 cm (approx. 8" x 9") Box of 5 dressings, individually wrapped
Oval	410510 – Oval shape - approx. 11 cm x 19 cm (approx. 4.5" x 7.5") Box of 5 dressings, individually wrapped

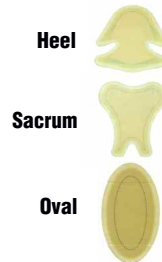

DuoDERM® Extra Thin CGF® (Control Gel Formula) Hydrocolloid Dressing, Sterile



DuoDERM® Extra Thin CGF® dressing is a highly flexible hydrocolloid dressing designed for use on dry to lightly exuding wounds, and for early intervention on skin at risk of breakdown.

	187900 – Rectangle - approx. 5 cm x 10 cm (approx. 2" x 4") Box of 20 dressings, individually wrapped
	187901 – Square - approx. 7 cm x 7 cm (approx. 3" x 3") Box of 20 dressings, individually wrapped
	187902 – Oval - approx. 10 cm x 15 cm (approx. 4" x 6") Box of 10 dressings, individually wrapped
	187903 – Triangle - approx. 15 cm x 18 cm (approx. 6" x 7") Box of 10 dressings, individually wrapped
	187955 – Square - approx. 10 cm x 10 cm (approx. 4" x 4") Box of 10 dressings, individually wrapped
	187957 – Square - approx. 14 cm x 14 cm (approx. 6" x 6") Box of 10 dressings, individually wrapped
	187961 – Rectangle - approx. 5 cm x 20 cm (approx. 2" x 8") Box of 10 dressings, individually wrapped
	413488 – Spots - approx. 4.4 cm X 3.8 cm (approx. 1.75" x 1.5") Box of 100 dressings individually wrapped in 10 packs of 10 spots each.

TIP!

To ensure a longer wearing time when applying any DuoDERM® dressing, always use a large enough piece. Leave a 3.2 cm (1¼") border of DuoDERM® onto the unaffected skin.

MOISTURE RETENTIVE DRESSINGS

Hydrocolloid Dressings

**DuoDERM® CGF® (Control Gel Formula) Hydrocolloid Dressing, Sterile**

DuoDERM® CGF® dressing is an adhesive hydrocolloid wound dressing designed for use on exuding dermal ulcers. It is indicated for use on dermal ulcers including full-thickness wounds, pressure ulcers (stage II-IV), leg ulcers, superficial wounds, partial-thickness burns and donor sites.



- 187658 – Square**, approx. 10 cm x 10 cm (approx. 4" x 4")
Box of 20 dressings, individually wrapped
- 187659 – Square**, approx. 14 cm x 14 cm (approx. 6" x 6")
Box of 20 dressings, individually wrapped
- 187660 – Square**, approx. 10 cm x 10 cm (approx. 4" x 4")
Box of 5 dressings, individually wrapped
- 187661 – Square**, approx. 14 cm x 14 cm (approx. 6" x 6")
Box of 5 dressings, individually wrapped
- 187662 – Square**, approx. 20 cm x 20 cm (approx. 8" x 8")
Box of 5 dressings, individually wrapped
- For larger coverage on leg ulcers, second degree burns or donor sites:**
- 187643 – Rectangle**, approx. 15 cm x 20 cm (approx. 6" x 8")
Box of 5 dressings, individually wrapped
- 187644 – Rectangle**, approx. 20 cm x 29 cm (approx. 8" x 12")
Box of 5 dressings, individually wrapped

DuoDERM® CGF® (Control Gel Formula) Border Hydrocolloid Dressing, Sterile*

DuoDERM® CGF® Border dressing includes an adhesive border so that additional taping (picture framing) is not required. It is designed for use on exuding wounds in hard-to-dress areas.

**Designed for the sacral region:**

- 187973 – Triangle**, Dimension of dressing including border:
Approx. 15 cm x 18 cm (approx. 6" x 7")
Dimension of DuoDERM® Area: approx. 10 cm x 13 cm (approx. 4" x 5")
Box of 5 dressings, individually wrapped
- 187974 – Triangle**, Dimension of dressing including border:
Approx. 20 cm x 22 cm (approx. 8" x 9")
Dimension of DuoDERM® Area: approx. 15 cm x 18 cm (approx. 6" x 7")
Box of 5 dressings, individually wrapped

Multi-purpose use:

- 187970 – Square**, Dimension of dressing including border:
Approx. 10 cm x 10 cm (approx. 4" x 4")
Dimension of DuoDERM® Area: approx. 6 cm x 6 cm (approx. 3" x 3")
Box of 5 dressings, individually wrapped
- 187971 – Square**, Dimension of dressing including border:
Approx. 14 cm x 14 cm (approx. 5" x 5")
Dimension of DuoDERM® Area: approx. 10 cm x 10 cm (approx. 4" x 4")
Box of 5 dressings, individually wrapped
- 187972 – Square**, Dimension of dressing including border:
Approx. 20 cm x 20 cm (approx. 8" x 8")
Dimension of DuoDERM® Area: approx. 15 cm x 15 cm (approx. 6" x 6")
Box of 5 dressings, individually wrapped

DuoDERM® Hydroactive® Hydrocolloid Dressing, Sterile

DuoDERM® Hydroactive® Hydrocolloid dressing is designed for the management of partial and full-thickness pressure ulcers and leg ulcers.



- 187611 – Square**, approx. 10 cm x 10 cm (approx. 4" x 4")
Box of 20 dressings, individually wrapped
(Institutional package)

* To ensure a longer wearing time when applying DuoDERM® CGF® Border dressings, always use a large enough piece. Leave a 3.2 cm (1 1/4") border of the DuoDERM® mass (not calculating the tape border) onto the unaffected skin.



DuoDERM® Hydroactive® Gel, Sterile



DuoDERM® Hydroactive® Gel is a sterile gel composed of natural hydrocolloids in a clear, viscous vehicle. It is designed for the hydration and management of partial and full-thickness wounds. DuoDERM® Gel forms a moist wound environment supportive of the healing process by aiding autolytic debridement.

- 187987 – 3 single use tubes of 30 g each
- 187990 – 10 single use tubes of 15 g each

COMPRESSION



COMPRESSION BANDAGES

SurePress® High Compression Bandage



REUSABLE

SurePress® High Compression Bandage is designed for use over a primary wound dressing for patients with impaired venous return. It is an elasticized, reusable bandage with a yellow center line and two rectangular extension indicators for different size limbs. It delivers 35-45 mmHg sustained compression.

Did you know? _____

The SurePress® High Compression Bandage can be washed up to 20 times. (hand-washed and laid flat to dry)

- 650947 – 10 cm x 3 m (4" x 3.2 yards) unstretched

SurePress® Absorbent Padding

A soft non-woven and absorbent bandage composed of a blend of viscose / bicomponent and super-absorbent fibers. This bandage protects the skin, provides additional absorption and even distribution of pressure under the SurePress® High Compression Bandage.

Did you know? _____

The SurePress® Absorbent Padding helps to distribute graduated compression and protect skin from exudates leaking from the primary dressing.

- 650948 – 10 cm x 3 m (4" x 3.2 yards)
Box of 6 bandages

SKIN CARE PRODUCTS

Bathing and Perineal Cleansing



1
CLEANSE

Cleanse skin of grease, grime, oils, and other bodily fluids.

2
MOISTURIZE

Moisturize to prevent and treat dry skin, improve skin tone and texture.

3
PROTECT

Protects and helps relieve chapped or cracked skin.

Aloe Vesta® Body Wash & Shampoo



Dispensing Pump: #100543



Aloe Vesta® Body Wash & Shampoo is a pH balanced formula designed to condition, soothe and moisturize the skin. Its formula is gentle to skin and hair and may be used for routine tub and shower bathing, or for no-rinse bathing and shampooing.

- 401865 – Case: 48 x 118 ml (4 oz.)
- 401866 – Case: 48 x 236 ml (8 oz.)
- 401867 – Case: 4 x 4 L (1 gal.)
- 324612 – Case: 8 X 1L (34 oz.)

Aloe Vesta® Cleansing Foam



Aloe Vesta® Cleansing Foam is a total body cleansing, no-rinse pH balanced formula, designed to soothe and moisturize. It is easy to use, provides rich foaming action while still being gentle to skin and hair.

- 325204 – Case: 24 x 118 ml (4 oz.)
- 401871 – Case: 12 x 236 ml (8 oz.)

Aloe Vesta® Perineal / Skin Cleanser

Dispensing Pump: #100543



Aloe Vesta® Perineal/Skin Cleanser is a no-rinse pH balanced formula, designed to soothe and moisturize while cleansing. Its lemon fragrance also helps with the management of odours.

- 401868 – Case: 48 x 118 ml (4 oz.)
- 401869 – Case: 48 x 236 ml (8 oz.)
- 401870 – Case: 4 x 3.6 L (approx. 1 gal.)

Sensi-Care® Perineal / Skin Cleanser



Dispensing Pump: #100543



Sensi-Care® Perineal/Skin Cleanser is a no-rinse pH balanced formula that is ideal for fragile, sensitive skin. It is designed to soothe and moisturize while cleansing.

- 401882 – Case: 48 x 118 ml (4 oz.)
- 401883 – Case: 48 x 236 ml (8 oz.)

(T) = Subject to Sales Tax.

SKIN CARE PRODUCTS

Moisturizing and Protective Skin Care



Aloe Vesta® Daily Moisturizer



Aloe Vesta® Daily Moisturizer is a lanolin and fragrance-free product formulated to soften and moisturize the skin, and to maintain that condition by slowing down the drying of the skin.

ALOE VESTA®

Dispensing Pump: #100543

401872 – Case: 48 x 118 ml (4 oz.)
401873 – Case: 48 x 236 ml (8 oz.)

Aloe Vesta® Protective Ointment



Aloe Vesta® Protective Ointment soothes, moisturizes, helps protect skin from irritants and seals out moisture. Routine use of this fragrance-free formula helps to maintain healthy looking skin.

ALOE VESTA®

401898 – Case: 24 x 59 ml (2 oz.)
401899 – Case: 12 x 236 ml (8 oz.)

Sensi-Care® Protective Barrier



Sensi-Care® Protective Barrier is a fragrance-free occlusive formula which seals out moisture and helps protect skin from irritants. It is ideal for sensitive skin and contains petrolatum and 15% zinc oxide.

SENSI-CARE®

325614 – Case: 24 x 102 ml (3.5 oz.)

Did you know?

Cleansing & moisturizing are vital to protect the skin.



Avelle™ Negative Pressure Wound Therapy System



The first system to combine negative pressure wound therapy with Hydrofiber® Technology. Negative Pressure Wound Therapy delivered in a simple, discrete, canister-less system: Easy to use - Small and portable - Quiet, no alarm - Disposable pump with 30 day lifetime† - Designed for flexibility - Suitable for use within and across different care settings - Separate pump and dressing packs - Supports patients resuming normal life.



Feature description of Hydrofiber® dressing:

1. Gentle silicone adhesive border*secures the dressing in place.
2. Stitch-bonded Hydrofiber® wound contact layer gels on contact with wound fluid and is designed to maintain dressing integrity on removal.*
3. Hydrofiber® core gels on contact with wound fluid locking in exudate and the bacteria it contains.* Channels in the Hydrofiber® core are designed to allow fluid to move throughout the dressing ensuring exudate is locked away.*
4. Foam layer is designed to aid distribution of negative pressure across the dressing.*
5. Film layer provides a bacterial and viral barrier, is showerproof, and permits evaporation of exudate aiding overall fluid handling.*

NEW

421551	Pump Unit
---------------	-----------------

NEW

	DRESSING SIZE	INNER PAD DIMENSION
421552	16 cm x 16 cm (approx. 6.3 in. x 6.3 in.)	8 cm x 8 cm (approx. 3.2 in. x 3.2 in.)
	Box of 5 dressings	
421553	16 cm x 21 cm (approx. 6.3 in. x 8.3 in.)	8 cm x 13 cm (approx. 3.2 in. x 5.1 in.)
	Box of 5 dressings	
421554	12 cm x 21 cm (approx. 4.7 in. x 8.3 in.)	4 cm x 13 cm (approx. 1.6 in. x 5.1 in.)
	Box of 5 dressings	
421555	12 cm x 31 cm (approx. 4.7 in. x 12.2 in.)	4 cm x 23 cm (approx. 1.6 in. x 9.1 in.)
	Box of 5 dressings	



1. ASSESSMENT OF THE IN-VITRO PROPERTIES AVELLE™ NEGATIVE PRESSURE WOUND THERAPY DRESSING. WHRI4520 MS128.
 Data on file. 2015. ConvaTec.
 * As demonstrated in-vitro.
 † Battery change may be required during pump lifetime



Glossary
Alphabetical Index
Numerical Index
Useful Resources

Authorization to
return merchandise

GLOSSARY

OSTOMY CARE TERMS

DEFINITION

Adhesive Coupling Technology™ (ACT)	A two-piece ostomy pouching system that sticks together with a tape-like adhesive rather than a plastic ring.
Appliance	The pouch or bag worn over the stoma (often called an ostomy pouching system).
Closed-end (pouch)	A pouch that has no opening at the bottom and is removed and thrown away after each use. This makes them ideal for special occasions, intimate moments, or for everyday use. Some are offered with a filter designed to reduce odor and gas buildup.
Comfort panel	A fine plastic mesh on one or both sides of the pouches, feels silky and comfortable on the skin. Closed-end pouches have a cotton-like finish.
Convexity	The outward curving of a faceplate which begins at the aperture of the faceplate and extends outward. ¹
Cut-to-Fit	The opening of the skin barrier or appliance needs to be enlarged (using scissors) to match the stoma size.
Drainable (pouch)	Has an opening at the bottom which allows stool or urine to be drained and re-closed with a clip or self gripping fastener hook tape, thus can be re-used. Urostomy pouches are drained by the valve or tap.
Durahesive® skin barrier	A durable hydrocolloid skin barrier designed to help manage liquid output. The skin barrier is resistant to caustic effluent.
Filter	Round piece of charcoal built-in to some pouches that helps release gas and/or odor and prevents the pouch from filling with air.
Flange	The plastic connector ring on the ostomy skin barrier or pouch.
InvisiClose® tail closure	A cliplless pouch that is closed via an integrated interlocking system.
Moldable	No scissors are needed to size a moldable skin barrier. The adhesive can be molded or shaped to the size of the stoma using one's fingers.
One-piece system	The skin barrier and pouch are combined into one unit for a one-step application. It is flexible and offers a low profile. The pouch cannot be removed without removing the skin barrier/wafer.
Peristomal skin	The skin surrounding a stoma.
Pre-Cut	The opening of the skin barrier or appliance is already cut to match the stoma size (No need of scissors).
Skin barrier	Wafer portion of an ostomy system that is in direct contact with the skin. It helps protect the skin around the stoma and provides a secure, leak-resistant seal with the pouch. ²
Skin protective barrier	A wipe or spray that forms a protective film on skin to help prevent skin breakdown by shielding it from irritating or caustic substances.
Stomahesive® skin barrier	A hydrocolloid skin barrier designed to adhere to wet/dry areas and is ideal for protection of the skin around the stoma.
Two-piece system	The skin barrier/wafer and pouch are separate and attached together with a round, plastic ring (called a flange) or by adhesive coupling. The pouch can be easily removed without also having to remove the skin barrier/wafer.

¹ Rolstad BS, Boarini J. Principles and techniques in the use of convexity. Ostomy/Wound Management; January/February 1996; (42)1: 26.

² Colwell JC. Principles of stoma management. In: Colwell JC, Goldberg MT, Carmel JE (eds). Fecal & Urinary Diversions: Management Principles. St.Louis: MO. Mosby 2004:240-262.

GLOSSARY

WOUND & SKIN CARE TERMS

DEFINITION

Bactericidal	Bacteria will be killed and therefore growth will not restart when contact with the antimicrobial agent is stopped.
Bacteriostatic	Bacterial growth will be inhibited, but will start again when contact with the antimicrobial agent is stopped.
Biofilm	Biofilm describes surface-adherent bacteria, encased in and protected by an extracellular polymeric substance (EPS) composed of polysaccharides, DNA and proteins.
Contamination	Soiling by contact or introduction of organisms into a wound. ¹
Hydrocolloid	A substance that yields a gel with water. ²
Hydrofiber® Technology	The Hydrofiber® Technology is exclusive to ConvaTec. It is composed of 100% pure sodium carboxymethylcellulose. It is the underlying technology behind AQUACEL®, AQUACEL® Ag, AQUACEL® Surgical, AQUACEL® Ag Surgical dressings.
Infection	An overgrowth of microorganisms capable of tissue invasion and destruction. Local and systemic symptoms ¹ may be observed.
Microbicidal	Bacteria, yeasts and fungi will be killed and therefore growth will not restart when contact with the antimicrobial agent is stopped.
Partial-thickness burn	(Second degree): A burn that damages the epidermal and some dermal tissues but does not damage the lower-lying hair follicles, sweat or sebaceous glands. ³

¹ Barton P, Parslow N. Wound Care-A Comprehensive Guide for Community Nurses. Saint Elizabeth Health Care; 1996: 86-87.

² Webster's Ninth New Collegiate Dictionary. Philippines Copyright by Merriam-Webster Inc.; 1984: 589.

³ Taber's Cyclopedic Medical Dictionary. Edition 19: F.A. Davis Company; 2001: 295.

ALPHABETICAL INDEX

Description	Page
A	
Accuseal® Adapter	40
Active Life® Closed-End Pouch with Filter and Stomahesive® Flexible Skin Barrier	26
Active Life® Closed-End Pouch with Filter and Stomahesive® Solid Skin Barrier	26
Active Life® Drainable Pouch with Convex Durahesive® Skin Barrier	27
Active Life® Drainable Pouch with Flexible Durahesive Plus® Skin Barrier	25
Active Life® Drainable Pouch with Stomahesive® Solid Skin Barrier	26
Active Life® Stoma Cap with Filter and Stomahesive® Solid Skin Barrier	26
Active Life® Urostomy Pouch with Bendable Tap and Stomahesive® Skin Barrier	27
AllKare® Adhesive Remover Wipe	36
AllKare® Protective Barrier Wipe	36
Aloe Vesta® Body Wash & Shampoo	57
Aloe Vesta® Cleansing Foam	57
Aloe Vesta® Daily Moisturizer	58
Aloe Vesta® Perineal / Skin Cleanser	57
Aloe Vesta® Protective Ointment	58
AQUACEL® Ag Burn Dressing	53
AQUACEL® Ag+ Extra™ Dressing	50
AQUACEL® Ag+ Ribbon Dressing	50
AQUACEL® Ag Extra Dressing with Strengthening Fiber	52
AQUACEL® Extra Dressing with Strengthening Fiber	52
AQUACEL® Foam Hydrofiber® Dressing	47
AQUACEL® Ag Foam Hydrofiber® Dressing	48
AQUACEL® Ag Hydrofiber® Wound Dressing	51
AQUACEL® Ag Hydrofiber® Ribbon Dressing	51
AQUACEL® Hydrofiber® Ribbon Dressing with Strengthening Fiber	51
AQUACEL® Hydrofiber® Wound Dressing	51
AQUACEL® Ag Surgical Cover Dressing	49
AQUACEL® Surgical Cover Dressing	49
Avelle™ Negative Pressure Wound Therapy System NEW	59
C	
CarboFlex® Odor Control Dressing	49
ConvaTec Loop Ostomy Rod	43
ConvaTec Moldable Technology™	12
ConvaTec Tail Closure	40
D	
Diamonds® Gelling and Odor Control Sachets	37
Dispensing Pump	57-58
Drain-Fix™ Securement Device	45
DuoDERM® CGF® Border Hydrocolloid Dressing	55
DuoDERM® CGF® Hydrocolloid Dressing	55
DuoDERM® Extra Thin CGF® Hydrocolloid Dressing	54
DuoDERM® Hydroactive® Gel	56
DuoDERM® Hydroactive® Hydrocolloid Dressing	55
DuoDERM Signal® A Tapered-Edge CGF® Hydrocolloid Dressing	54
DuoLock® Curved Tail Closure	40
E	
Eakin Cohesive® Fistula & Wound Pouch	44
Eakin Cohesive® Paste NEW	38
Eakin Cohesive® Seal	38
Eakin Cohesive® Skin Barrier	38
Esteem® ⊕ Accuseal® Tap, Urostomy Pouches with Durahesive® Skin Barrier, Cut-to-Fit NEW	31
Esteem® ⊕ Convex Durahesive® drainable pouches with InvisiClose®, Lock-It Pocket™	32
Esteem® ⊕ Durahesive®, drainable pouches with InvisiClose®, Lock-It Pocket™	31
Esteem® ⊕ Flex Convex	32
Esteem® ⊕ Moldable One-piece System	29

ALPHABETICAL INDEX

Description	Page
Esteem®  Stomahesive® Closed-end pouches with filter.....	30
Esteem®  Stomahesive®, drainable pouches with InvisiClose®, Lock-It Pocket™.....	31
Esteem®  Stomahesive® Post-Operative Kits with drainable pouches, InvisiClose®, Lock-It Pocket™.....	43
Esteem synergy® Durahesive® Moldable Convex Flexible Skin Barrier.....	22
Esteem synergy® Professional Kit.....	43
Esteem synergy® Stomahesive® Flexible Skin Barrier.....	21
Esteem synergy® Stomahesive Skin Barrier with ConvaTec Moldable Technology™.....	22
Esteem synergy® Stomahesive® Solid Skin Barrier.....	21
Esteem synergy® Straight Drainable Pouch.....	22
Esteem synergy® Urostomy Pouch with Accuseal® Tap.....	24
Esteem synergy®  Drainable pouches with InvisiClose® and Lock-It Pocket™.....	23
Esteem synergy®  Closed-end pouches.....	24
F	
Flexi-Seal® FMS (Fecal Management System).....	45
Flexi-Seal® FMS Collection Bags.....	45
Flexi-Seal® PROTECT FMS.....	45
Flexi-Seal® SIGNAL™ FMS.....	45
FLEXI-TRAK® Anchoring Device.....	45
FoamLite™ ConvaTec Dressing.....	48
G	
GentleCath™ Intermittent Urinary Catheter™.....	46
K	
KALTOSTAT® Calcium Alginate Wound Dressing.....	53
L	
Little Ones® Closed-end Pouch with Filter Adhesive Coupling Technology™ (ACT).....	35
Little Ones® Drainable Pouch with Filter and InvisiClose® (Clipless) Tail Closure Adhesive Coupling Technology™ (ACT) ..	35
Little Ones® One-Piece Drainable Pouch with InvisiClose® (Clipless) Tail Closure.....	33
Little Ones® One-Piece Drainable Pouch with Stomahesive® Skin Barrier.....	33
Little Ones® One-Piece Urostomy Pouch with Stomahesive® Skin Barrier.....	33
Little Ones® Stomahesive® Flexible Skin Barrier with Adhesive Coupling Technology™ (ACT).....	35
Little Ones® Stomahesive® Solid Skin Barrier with Adhesive Coupling Technology™ (ACT).....	35
Little Ones® Two-Piece Drainable Pouch.....	34
Little Ones® Two-Piece Mini-Pouch.....	34
Little Ones® Two-Piece Urostomy Pouch.....	34
Little Ones® Two-Piece Stomahesive® Flexible Skin Barrier.....	34
N	
Natura® Closed-End Pouch.....	18
Natura® Closed-End Pouch with Filter.....	18
Natura® Disposable Convex Inserts.....	7
Natura® Drainable Pouch.....	15
Natura® Drainable Pouch with Filter.....	16
Natura® Drainable Pouch with Filter and InvisiClose® (Clipless) Tail Closure.....	17
Natura® Drainable Pouch with InvisiClose® (Clipless) Tail Closure.....	17
Natura® Durahesive® Flexible Skin Barrier.....	10
Natura® Durahesive® Convex Cut-to-Fit Skin Barrier with ACCORDION flange (AC) or (HC).....	14
Natura® Durahesive® Moldable Skin Barrier.....	12
Natura® Durahesive® Moldable Skin Barrier with Accordion Flange.....	13
Natura® Durahesive® Cut-to-Fit Skin Barrier with Accordion Flange.....	14
Natura® Durahesive® Moldable Convex Flexible Skin Barrier.....	12
Natura® Durahesive® Flexible Skin Barrier with Convex-It®.....	11
Natura® Durahesive® Solid Skin Barrier.....	10
Natura® Flange Cap with Filter.....	18
Natura® High Output Drainable Pouch with Filter.....	16
Natura® Low-Pressure Adaptor.....	7
Natura® Moldable Skin Barrier.....	12

ALPHABETICAL INDEX

Description	Page
Natura® Mini-Pouch.....	18
Natura® Postoperative Kit	41
Natura® Durahesive® Convex Cut-to-Fit Skin Barrier with ACCORDION flange	14
Natura® Stomahesive® Flexible Skin Barrier	8
Natura® Stomahesive® Moldable Skin Barrier with Accordion flange.....	13
Natura® Stomahesive® Cut-to-Fit Skin Barrier with Accordion flange.....	14
Natura® Stomahesive® Solid Skin Barrier.....	8
Natura® Urostomy Pouch with Accuseal® Tap	19
Natura® Urostomy Pouch with Bendable Tap.....	19
Natura® Visi-Flow® Irrigation Adapter Faceplate	39
Natura® Visi-Flow® Irrigation Sleeve	39
Natura® Visi-Flow® Irrigation Sleeve Tail Closure	40
Natura® ☞ Drainable pouches with InvisiClose® and Lock-It Pocket™.....	17
Natura® ☞ Closed-End Pouch.....	18
Natura® ☞ Postoperative Kit	41
Natura® ☞ Postoperative Kit including Moldable Skin Barrier with ACCORDION Flange	42
Natura® ☞ Loop Ostomy Postoperative Kit	43
Natura® ☞ Urostomy Pouch with Soft Tap.....	19
Night Drainage Container Set	40
Night Drainage Container Tube	40
Niltac™ Sting Free Adhesive Remover Wipes.....	37
Niltac™ Sting Free Adhesive Releaser Spray.....	37
O	
Orabase® Oral Protective Paste	37
Orahesive® Denture Adhesive (Powder)	37
Ostomy Appliance Belt	40
S	
Sensi-Care® Perineal / Skin Cleanser	57
Sensi-Care® Protective Barrier	58
Sillesse™ Sting Free Protective Skin Barrier Wipes.....	37
Sillesse™ Sting Free Protective Skin Barrier Spray	37
Stomahesive® Paste	36
Stomahesive® Protective Powder	36
Stomahesive® Skin Barrier	36
Stomahesive® Strip	37
SurePress® Absorbent Padding	56
SurePress® High Sustained Compression Bandage.....	56
V	
Visi-Flow® Irrigator	39
Visi-Flow® Stoma Cone	39

NUMERICAL INDEX

Product Code	UPC Codes	Page	Product Code	UPC Codes	Page
680M	5060038510310	.45	125342	7 68455 10203 7	.25
685M	5060038511997	.45	125343	7 68455 10204 4	.25
20917	7 68455 10869 5	.33	125352	7 68455 10206 8	.27
20922	7 68455 10139 9	.33	125353	7 68455 10207 5	.27
21712	7 68455 10694 3	.36	125354	7 68455 10208 2	.27
21715	7 68455 10696 7	.36	125355	7 68455 10209 9	.27
21915	3 00030 21915 2	.40	125356	7 68455 10210 5	.27
22355	7 68455 11550 1	.43	125357	7 68455 10211 2	.27
22356	7 68455 11551 8	.43	125358	7 68455 10212 9	.27
22736	7 68455 10708 7	.39	125359	7 68455 10213 6	.27
22751	7 68455 10141 2	.25	125361	7 68455 10214 3	.27
22764	7 68455 10154 2	.25	125362	7 68455 10215 0	.27
22765	7 68455 10155 9	.25	125363	7 68455 10216 7	.27
22766	7 68455 10156 6	.25	125364	7 68455 10217 4	.27
22767	7 68455 10157 3	.25	125365	7 68455 10218 1	.27
22768	7 68455 10158 0	.25	125366	7 68455 10219 8	.27
22769	7 68455 10159 7	.25	125367	7 68455 10220 4	.27
22771	7 68455 10161 0	.25	125368	7 68455 10221 1	.27
25510	7 68455 10826 8	.36	125369	7 68455 10222 8	.27
25542	7 68455 10702 5	.37	168117	7 68455 08116 5	.53
25545	7 68455 10986 9	.37	168210	7 68455 08114 1	.53
27060	7 68455 10913 5	.40	168212	7 68455 08115 8	.53
27062	7 68455 10709 4	.40	168214	7 68455 08111 0	.53
037436	7 68455 10767 4	.36	175507	7 68455 10718 6	.40
037439	7 68455 10766 7	.36	175611	7 68455 10236 5	.26
37443	7 68455 10760 5	.36	175635	7 68455 10719 3	.40
37444	7 68455 10761 2	.36	175650	7 68455 10707 0	.40
37449	7 68455 10762 9	.45	175652	7 68455 10867 1	.40
62436	7 68455 10876 3	.37	175769	7 68455 10245 7	.26
64927	7 68455 10168 9	.27	175770	7 68455 10247 1	.26
100543	0 40000 10054 3	.57-58	175771	7 68455 10248 8	.26
125257	7 68455 10170 2	.8	175772	7 68455 10249 5	.26
125258	7 68455 10171 9	.8	175773	7 68455 10250 1	.26
125259	7 68455 10172 6	.8	175777	7 68455 10251 8	.27
125260	7 68455 10173 3	.8	175778	7 68455 10252 5	.27
125261	7 68455 10174 0	.8	175779	7 68455 10253 2	.27
125262	7 68455 10175 7	.8	175780	7 68455 10254 9	.27
125263	7 68455 10176 4	.8	177901	7 68455 10935 7	.51
125264	7 68455 10177 1	.8	177902	7 68455 10936 4	.51
125265	7 68455 10178 8	.8	177903	7 68455 10937 1	.51
125266	7 68455 10179 5	.8	183910	7 68455 10690 5	.36
125267	7 68455 10180 1	.9	187611	7 68455 10693 6	.55
125268	7 68455 10181 8	.9	187643	7 68455 10269 3	.55
125269	7 68455 10182 5	.9	187644	7 68455 10670 7	.55
125270	7 68455 10183 2	.9	187658	7 68455 10270 9	.55
125271	7 68455 10184 9	.9	187659	7 68455 10271 6	.55
125272	7 68455 10185 6	.9	187660	7 68455 10697 4	.55
125273	7 68455 10186 3	.9	187661	7 68455 10272 3	.55
125274	7 68455 10187 0	.9	187662	7 68455 10698 1	.55
125275	7 68455 10188 7	.9	187900	7 68455 10277 8	.54
125276	7 68455 10189 4	.9	187901	7 68455 10278 5	.54
125277	7 68455 10190 0	.9	187902	7 68455 10279 2	.54
125330	7 68455 10191 7	.25	187903	7 68455 10280 8	.54
125331	7 68455 10192 4	.25	187955	7 68455 10691 2	.54
125333	7 68455 10194 8	.25	187957	7 68455 10692 9	.54
125334	7 68455 10195 5	.25	187961	7 68455 10282 2	.54
125335	7 68455 10196 2	.25	187970	7 68455 10283 9	.55
125336	7 68455 10197 9	.25	187971	7 68455 10284 6	.55
125338	7 68455 10199 3	.25	187972	7 68455 10285 3	.55
125339	7 68455 10200 6	.25	187973	7 68455 10286 0	.55
125340	7 68455 10201 3	.25	187974	7 68455 10287 7	.55
125341	7 68455 10202 0	.25	187987	7 68455 10706 3	.56

NUMERICAL INDEX

Product Code	UPC Codes	Page	Product Code	UPC Codes	Page
187990	7 68455 10864 0	.56	401575	7 68455 10377 5	.8
324612	7 68455 10857 2	.57	401576	7 68455 10378 2	.8
325204	7 68455 10808 4	.57	401577	7 68455 10379 9	.8
325614	7 68455 08848 5	.58	401865	7 68455 09980 1	.57
400598	7 68455 10289 1	.25	401866	7 68455 06300 0	.57
400599	7 68455 10290 7	.25	401867	7 68455 06301 7	.57
401500	7 68455 10319 5	.15	401868	7 68455 09981 8	.57
401501	7 68455 10321 8	.15	401869	7 68455 06303 1	.57
401502	7 68455 10323 2	.15	401870	7 68455 06304 8	.57
401503	7 68455 10325 6	.15	401871	7 68455 06305 5	.57
401504	7 68455 10326 3	.15	401872	7 68455 06306 2	.58
401505	7 68455 10327 0	.15	401873	7 68455 06307 9	.58
401506	7 68455 10328 7	.15	401882	7 68455 09979 5	.57
401507	7 68455 10329 4	.15	401883	7 68455 06317 8	.57
401508	7 68455 10330 0	.15	401898	7 68455 06310 9	.58
401509	7 68455 10331 7	.15	401899	7 68455 06309 3	.58
401510	7 68455 10332 4	.15	401905	7 68455 10393 5	.10
401511	7 68455 10333 1	.15	401906	7 68455 10394 2	.8
401512	7 68455 10334 8	.15	401909	7 68455 10397 3	.18
401513	7 68455 10335 5	.15	401911	7 68455 10399 7	.39
401514	7 68455 10336 2	.15	401912	7 68455 10400 0	.39
401515	7 68455 10337 9	.15	401913	7 68455 10401 7	.39
401521	7 68455 10338 6	.18	401914	7 68455 10402 4	.39
401522	7 68455 10339 3	.18	401918	7 68455 10405 5	.39
401523	7 68455 10340 9	.18	401919	7 68455 10406 2	.39
401524	7 68455 10341 6	.18	401920	7 68455 10407 9	.41
401525	7 68455 10342 3	.18	401921	7 68455 10408 6	.41
401526	7 68455 10343 0	.18	401922	7 68455 10409 3	.41
401527	7 68455 10344 7	.18	401923	7 68455 10410 9	.41
401528	7 68455 10345 4	.18	401925	7 68455 10412 3	.34
401529	7 68455 10346 1	.18	401927	7 68455 10414 7	.34
401530	7 68455 10347 8	.18	401928	7 68455 10415 4	.34
401531	7 68455 10348 5	.18	401929	7 68455 10416 1	.34
401532	7 68455 10349 2	.18	401930	7 68455 10417 8	.34
401533	7 68455 10350 8	.19	401989	7 68455 10425 3	.39
401534	7 68455 10351 5	.19	401993	7 68455 10861 9	.7
401535	7 68455 10352 2	.19	401994	7 68455 10860 2	.7
401536	7 68455 10353 9	.19	401995	7 68455 10859 6	.7
401537	7 68455 10354 6	.19	401996	7 68455 10858 9	.7
401538	7 68455 10671 4	.19	403202	7 68455 10428 4	.49
401539	7 68455 10355 3	.19	403203	7 68455 10429 1	.49
401540	7 68455 10356 0	.19	403204	7 68455 10430 7	.49
401541	7 68455 10357 7	.19	403326	7 68455 10725 4	.54
401542	7 68455 10358 4	.19	403327	7 68455 10723 0	.54
401543	7 68455 10359 1	.19	403328	7 68455 10726 1	.54
401544	7 68455 10360 7	.19	403332	7 68455 10724 7	.54
401545	7 68455 10361 4	.19	403333	7 68455 10727 8	.54
401546	7 68455 10362 1	.19	403706	7 68455 10946 3	.51
401547	7 68455 10363 8	.19	403708	7 68455 10947 0	.51
401548	7 68455 10364 5	.19	403710	7 68455 10948 7	.51
401549	7 68455 10365 2	.19	403711	7 68455 10951 7	.51
401550	7 68455 10366 9	.19	403770	7 68455 10677 6	.51
401551	7 68455 10367 6	.19	403771	7 68455 10674 5	.51
401552	7 68455 10368 3	.19	403786	7 68455 11051 3	.53
401553	7 68455 10369 0	.19	403787	7 68455 11054 4	.53
401554	7 68455 10370 6	.19	403788	7 68455 11057 5	.53
401555	7 68455 10371 3	.19	403789	7 68455 11060 5	.53
401557	7 68455 10372 0	.16	403791	7 68455 11066 7	.53
401558	7 68455 10373 7	.16	403792	7 68455 11069 8	.53
401559	7 68455 10374 4	.16	403793	7 68455 11072 8	.53
401573	7 68455 10375 1	.8	403794	7 68455 11075 9	.53
401574	7 68455 10376 8	.8	403795	7 68455 11078 0	.53

NUMERICAL INDEX

Product Code	UPC Codes	Page	Product Code	UPC Codes	Page
403948	7 68455 10434 5	.21	411310	7 68455 10571 7	.17
404006	7 68455 10710 0	.7	411311	7 68455 10572 4	.17
404007	7 68455 10711 7	.7	411312	7 68455 10573 1	.17
404008	7 68455 10712 4	.7	411360	7 68455 10574 8	.17
404009	7 68455 10713 1	.7	411361	7 68455 10575 5	.17
404010	7 68455 10714 8	.7	411362	7 68455 10570 0	.17
404011	7 68455 10715 5	.7	411490	7 68455 10581 6	.16
404012	7 68455 10716 2	.7	411491	7 68455 10582 3	.16
404013	7 68455 10717 9	.7	411492	7 68455 10583 0	.16
404014	7 68455 10437 6	.15	411493	7 68455 10584 7	.16
404015	7 68455 10438 3	.15	411633	7 68455 10587 8	.33
404016	7 68455 10439 0	.15	411637	7 68455 10588 5	.35
404017	7 68455 10440 6	.15	411639	7 68455 10590 8	.35
404018	7 68455 10441 3	.15	411640	7 68455 10591 5	.35
404025	7 68455 10442 0	.15	411641	7 68455 10592 2	.35
404026	7 68455 10443 7	.15	411642	7 68455 10593 9	.35
404027	7 68455 10444 4	.15	411657	7 68455 11733 8	.22
404028	7 68455 10445 1	.15	411658	7 68455 11734 5	.22
404029	7 68455 10446 8	.15	411659	7 68455 11735 2	.22
404030	7 68455 10447 5	.15	411800	7 68455 10594 6	.12
404031	7 68455 10448 2	.15	411801	7 68455 10595 3	.12
404032	7 68455 10449 9	.15	411802	7 68455 10596 0	.12
404033	7 68455 10450 5	.15	411803	7 68455 10597 7	.12
404034	7 68455 10451 2	.15	411804	7 68455 10598 4	.12
404592	7 68455 10452 9	.12	411805	7 68455 10599 1	.12
404593	7 68455 10453 6	.12	411806	7 68455 10600 4	.12
404594	7 68455 10454 3	.12	411807	7 68455 10601 1	.12
405406	7 68455 10477 2	.22	412009	7 68455 11116 9	.49
405451	7 68455 10491 8	.24	412010	7 68455 11119 0	.49
405456	7 68455 10496 3	.21	412011	7 68455 11122 0	.49
405457	7 68455 10497 0	.21	412012	7 68455 11125 1	.49
405458	7 68455 10498 7	.21	412017	7 68455 11128 2	.49
405467	7 68455 10499 4	.21	412018	7 68455 11131 2	.49
405468	7 68455 10500 7	.21	412019	7 68455 11134 3	.49
405474	7 68455 10505 2	.21	412020	7 68455 11137 4	.49
405475	7 68455 10507 6	.21	413142	7 68455 10603 5	.26
405476	7 68455 10508 3	.21	413143	7 68455 10604 2	.26
405477	7 68455 10509 0	.21	413144	7 68455 10605 9	.26
405478	7 68455 10510 6	.21	413145	7 68455 10606 6	.26
405479	7 68455 10511 3	.21	413146	7 68455 10607 3	.26
405480	7 68455 10512 0	.21	413153	7 68455 10614 1	.10
405482	7 68455 10514 4	.21	413154	7 68455 10615 8	.10
405483	7 68455 10515 1	.21	413155	7 68455 10616 5	.10
409268	7 68455 10523 6	.22	413156	7 68455 10617 2	.10
409269	7 68455 10524 3	.22	413157	7 68455 10618 9	.10
409270	7 68455 10525 0	.22	413159	7 68455 10619 6	.10
409275	7 68455 10526 7	.43	413160	7 68455 10620 2	.10
410500	7 68455 10747 6	.54	413161	7 68455 10621 9	.10
410501	7 68455 10748 3	.54	413162	7 68455 10622 6	.10
410510	7 68455 10750 6	.54	413163	7 68455 10623 3	.10
411100	7 68455 10856 5	.45	413164	7 68455 10624 0	.10
411101	7 68455 10790 2	.45	413165	7 68455 10625 7	.10
411102	7 68455 10862 6	.45	413166	7 68455 10626 4	.10
411108	7 68455 12572 2	.45	413167	7 68455 10627 1	.10
411264	7 68455 10561 8	.16	413168	7 68455 10628 8	.10
411265	7 68455 10562 5	.16	413177	7 68455 10637 0	.11
411266	7 68455 10563 2	.16	413178	7 68455 10638 7	.11
411267	7 68455 10564 9	.16	413179	7 68455 10639 4	.11
411289	3 00034 11289 3	.16	413180	7 68455 10640 0	.11
411290	7 68455 10567 0	.16	413181	7 68455 10641 7	.11
411291	7 68455 10568 7	.16	413182	7 68455 10642 4	.11
411309	7 68455 10576 2	.17	413183	7 68455 10643 1	.11

NUMERICAL INDEX

Product Code	UPC Codes	Page	Product Code	UPC Codes	Page
413184	7 68455 10644 8	11	416420	7 68455 11338 5	17
413185	7 68455 10645 5	11	416421	7 68455 11339 2	17
413186	7 68455 10646 2	11	416422	7 68455 11340 8	17
413187	7 68455 10647 9	11	416423	7 68455 11341 5	17
413188	7 68455 10648 6	11	416424	7 68455 11342 2	17
413312	7 68455 10650 9	17	416472	7 68455 11572 3	17
413313	7 68455 10651 6	17	416700	7 68455 11170 1	30
413314	7 68455 10652 3	17	416701	7 68455 11171 8	30
413332	7 68455 10653 0	17	416702	7 68455 11172 5	30
413334	7 68455 10655 4	17	416703	7 68455 11173 2	30
413335	7 68455 10656 1	17	416704	7 68455 11174 9	30
413336	7 68455 10657 8	17	416705	7 68455 11175 6	30
413417	7 68455 12391 9	12	416706	7 68455 11176 3	30
413418	7 68455 12392 6	12	416707	7 68455 11177 0	30
413419	7 68455 12393 3	12	416709	7 68455 11179 4	30
413420	7 68455 12394 0	12	416710	7 68455 11180 0	30
413421	7 68455 12395 7	12	416711	7 68455 11181 7	30
413422	7 68455 12396 4	12	416712	7 68455 11182 4	30
413423	7 68455 12397 1	12	416713	7 68455 11183 1	30
413424	7 68455 12398 8	12	416714	7 68455 11184 8	30
413435	7 68455 12881 5	19	416715	7 68455 11185 5	30
413436	7 68455 12883 9	19	416718	7 68455 11188 6	31
413437	7 68455 12885 3	19	416719	7 68455 11189 3	31
413438	7 68455 12887 7	19	416721	7 68455 11191 6	31
413439	7 68455 12889 1	19	416724	7 68455 11194 7	31
413488	7 68455 12879 2	54	416725	7 68455 11195 4	31
413509	7 68455 13005 4	29	416728	7 68455 11198 5	31
413510	7 68455 13006 1	29	416729	7 68455 11199 2	31
413511	7 68455 13007 8	29	416732	7 68455 11202 9	31
413515	7 68455 13017 7	29	416733	7 68455 11203 6	31
413516	7 68455 13018 4	29	416736	7 68455 11206 7	31
413517	7 68455 13019 1	29	416737	7 68455 11207 4	31
413518	7 68455 13020 7	29	416738	7 68455 11208 1	32
413519	7 68455 13021 4	29	416739	7 68455 11209 8	32
413521	7 68455 13023 8	29	416740	7 68455 11210 4	32
413522	7 68455 13024 5	29	416741	7 68455 11211 1	32
413523	7 68455 13025 2	29	416742	7 68455 11212 8	32
413566	7 68455 13206 5	50	416743	7 68455 11213 5	32
413567	7 68455 13207 2	50	416744	7 68455 11214 2	32
413568	7 68455 13208 9	50	416745	7 68455 11215 9	32
413569	7 68455 13209 6	50	416746	7 68455 11216 6	32
413570	7 68455 13210 2	50	416747	7 68455 11217 3	32
413571	7 68455 13261 4	50	416748	7 68455 11218 0	32
416400	7 68455 11318 7	18	416749	7 68455 11219 7	32
416401	7 68455 11319 4	18	416750	7 68455 11220 3	32
416402	7 68455 11320 0	18	416751	7 68455 11221 0	32
416403	7 68455 11321 7	18	416774	7 68455 11460 3	24
416404	7 68455 11322 4	18	416775	7 68455 11461 0	24
416405	7 68455 11323 1	18	416776	7 68455 11462 7	24
416406	7 68455 11324 8	18	416777	7 68455 11463 4	24
416407	7 68455 11325 5	18	416778	7 68455 11464 1	24
416408	7 68455 11326 2	18	416779	7 68455 11465 8	24
416409	7 68455 11327 9	18	416780	7 68455 11466 5	24
416410	7 68455 11328 6	18	416781	7 68455 11467 2	24
416411	7 68455 11329 3	18	416782	7 68455 11468 9	24
416412	7 68455 11330 9	18	416783	7 68455 11469 6	24
416413	7 68455 11331 6	18	416795	7 68455 11481 8	23
416415	7 68455 11333 0	17	416796	7 68455 11482 5	23
416416	7 68455 11334 7	17	416797	7 68455 11483 2	23
416417	7 68455 11335 4	17	416798	7 68455 11484 9	23
416418	7 68455 11336 1	17	416799	7 68455 11485 6	23
416419	7 68455 11337 8	17	416900	7 68455 11486 3	23

NUMERICAL INDEX

Product Code	UPC Codes	Page	Product Code	UPC Codes	Page
416901	7 68455 11487 0	.23	420787	7 68455 12983 6	.37
416902	7 68455 11488 7	.23	420788	7 68455 12984 3	.37
416903	7 68455 11489 4	.23	420789	7 68455 12985 0	.37
416908	7 68455 11812 0	.31	420790	7 68455 12986 7	.37
416915	7 68455 11649 2	.43	420791	5 06014 8230085	.37
416916	7 68455 11587 7	.43	420804	7 68455 12909 6	.47
416917	7 68455 11588 4	.41	420805	7 68455 12913 3	.48
416918	7 68455 11589 1	.41	420806	7 68455 12917 1	.48
416919	7 68455 11590 7	.41	420807	7 68455 12921 8	.48
416920	7 68455 11591 4	.42	421033	7 68455 13307 9	.13
416921	7 68455 11592 1	.41	421034	7 68455 13308 6	.13
416922	7 68455 11593 8	.41	421035	7 68455 13309 3	.13
416923	7 68455 11594 5	.41	421039	7 68455 13319 2	.13
416924	7 68455 11595 2	.42	421040	7 68455 13320 8	.13
416925	7 68455 11596 9	.41	421041	7 68455 13321 5	.13
416926	7 68455 11597 6	.41	421042	7 68455 13322 2	.42
416927	7 68455 11598 3	.43	421043	7 68455 13323 9	.42
416928	7 68455 11599 0	.43	421044	7 68455 13324 6	.42
416929	7 68455 11600 3	.43	421045	7 68455 13325 3	.42
416930	7 68455 11601 0	.43	421046	7 68455 13326 0	.42
416938	7 38455 11609 6	.42	421047	7 68455 13327 7	.42
416946	7 68455 11617 1	.42	421150	7 68455 13641 4	.47
416954	7 68455 11625 6	.42	421152	7 68455 13643 8	.47
416975	7 68455 11533 4	.31	421154	7 68455 13645 2	.47
416976	7 68455 11535 8	.31	421156	7 68455 13647 6	.47
418000	7 68455 11717 8	.45	421454	7 68455 13826 5	.14
420127	7 68455 11821 2	.51	421456	7 68455 13830 2	.14
420128	7 68455 11817 5	.51	421458	7 68455 13834 0	.14
420619	7 68455 12287 5	.47	421460	7 68455 13838 8	.14
420621	7 68455 12291 2	.47	421462	7 68455 13842 5	.14
420623	7 68455 12294 3	.47	421464	7 68455 13846 3	.14
420624	7 68455 12298 1	.47	421466	7 68455 13850 0	.14
420625	7 68455 12302 5	.47	421468	7 68455 13854 8	.14
420626	7 68455 12306 3	.47	421551	7 68455 13919 4	.59
420627	7 68455 12311 7	.48	421552	7 68455 13920 8	.59
420628	7 68455 12315 5	.48	421553	7 68455 13921 5	.59
420629	7 68455 12318 6	.48	421554	7 68455 13922 2	.59
420631	7 68455 12331 5	.47	421555	7 68455 13923 9	.59
420633	7 68455 12336 0	.47	421557	7 68455 14227 9	.48
420635	7 68455 12341 4	.47	421559	7 68455 14230 9	.48
420636	7 68455 12345 2	.47	421561	7 68455 14233 0	.48
420637	7 68455 12349 0	.47	421563	7 68455 14236 1	.48
420639	7 68455 12354 4	.48	421612	7 68455 14315 3	.32
420642	7 68455 12359 9	.48	421613	7 68455 14317 7	.32
420645	7 68455 12364 3	.48	421614	7 68455 14319 1	.32
420646	7 68455 12368 1	.48	421615	7 68455 14321 4	.32
420647	7 68455 12322 3	.48	421619	7 68455 14329 0	.32
420648	7 68455 12326 1	.48	421622	7 68455 14335 1	.32
420669	7 68455 12507 4	.49	421623	7 68455 14337 5	.32
420670	7 68455 12511 1	.49	421626	7 68455 14343 6	.32
420671	7 68455 12559 3	.52	421628	7 68455 14347 4	.32
420672	7 68455 12665 1	.52	421630	7 68455 14301 6	.45
420673	7 68455 12564 7	.52	421631	7 68455 14609 3	.31
420675	7 68455 12738 2	.52	421638	7 68455 14374 0	.14
420676	7 68455 12741 2	.52	421639	7 68455 14376 4	.14
420678	7 68455 12746 7	.52	421640	7 68455 14378 8	.14
420679	7 68455 12749 8	.52	421641	7 68455 14380 1	.14
420680	7 68455 12712 2	.47	421642	7 68455 14382 5	.14
420681	7 68455 12716 0	.48	421643	7 68455 14384 9	.14
420695	7 68455 12761 0	.16	421927	7 68455 15112 7	.48
420696	7 68455 12762 7	.16	421928	7 68455 15115 8	.48
420697	7 68455 12763 4	.16	501002	7 68455 11140 4	.46

NUMERICAL INDEX

Product Code	UPC Codes	Page	Product Code	UPC Codes	Page
501003	7 68455 11141 1	.46			
501004	7 68455 11142 8	.46			
501005	7 68455 11143 5	.46			
501006	7 68455 11144 2	.46			
501011	7 68455 11149 7	.46			
501012	7 68455 11150 3	.46			
501013	7 68455 11151 0	.46			
501014	7 68455 11152 7	.46			
501015	7 68455 11153 4	.46			
501016	7 68455 11154 1	.46			
501019	7 68455 11157 2	.46			
501020	7 68455 11158 9	.46			
501021	7 68455 11159 6	.46			
501022	7 68455 11160 2	.46			
501023	7 68455 11161 9	.46			
650408	7 68455 10659 2	.26			
650422	7 68455 10660 8	.26			
650803	7 68455 10875 6	.26			
650828	7 68455 10661 5	.27			
650829	7 68455 10662 2	.27			
650830	7 68455 10663 9	.27			
650831	7 68455 10664 6	.27			
650832	7 68455 10665 3	.27			
650833	7 68455 10666 0	.27			
650947	7 68455 10782 7	.56			
650948	7 68455 10896 1	.56			
839001	7 68455 10872 5	.38			
839002	7 68455 10873 2	.38			
839003	7 68455 10933 3	.38			
839004	7 68455 10934 0	.38			
839005	7 68455 11825 0	.38			
839006	7 68455 13917 0	.38			
839010	(01)05060267416704	.38			
839250	7 68455 12687 3	.44			
839251	7 68455 12688 0	.44			
839252	7 68455 12689 7	.44			
839253	7 68455 12679 8	.44			
839260	7 68455 12690 3	.44			
839261	7 68455 12691 0	.44			
839262	7 68455 12692 7	.44			
839263	7 68455 12680 4	.44			
839264	7 68455 12682 8	.44			
839265	7 68455 12681 1	.44			
839266	7 68455 12686 6	.44			
839267	7 68455 12683 5	.44			
839268	7 68455 12684 2	.44			
839280	7 68455 12685 9	.44			

USEFUL RESOURCES

Ostomy Canada Society

5800 Ambler Dr., Suite 210,
Mississauga, ON L4W 4J4
Tel.: 1 905 212-7111
Toll Free: 1 888 969-9698
Fax: 1 905 212-9002
E-mail: info1@ostomycanada.ca
www.ostomycanada.ca

Canadian Cancer Society

55 St. Clair Avenue West, Suite 300
Toronto, ON M4V 2Y7
Tel.: (416) 961-7223
Toll Free: 1 888 939-3333
Fax: (416) 961-4189
Email: info@cis.cancer.ca
www.cancer.ca

Crohn's & Colitis Foundation of Canada

600-60 St. Clair Avenue East
Toronto, ON M4T 1N5
Tel.: (416) 920-5035
Toll Free: 1 800 387-1479
Fax: (416) 929-0364
E-mail: ccfc@ccfc.ca
www.ccfc.ca

Canadian Association for Enterostomal Therapy (CAET)

66 Leopolds Drive
Ottawa, Ontario
K1V 7E3
Tel: 888-739-5072
Fax: 613-834-6351
E-mail: office@caet.ca
www.caet.ca

Colorectal Cancer Association of Canada

Montreal

1350 Sherbrooke St West, Suite 300
Montreal, QC H3G 1J1
E-mail: information@colorectal-cancer.ca
Tel: (514) 875-7745
Fax: (514) 875-7746

Toronto

2 Bloor St. West, Suite 700
Toronto, ON M4W 3R1
E-mail: information@colorectal-cancer.ca
Tel: 1 877 502-6566
www.colorectal-cancer.ca

Friends of Ostomates Worldwide (FOW)

Suite 609, 344 Bloor St. W.,
Toronto, ON M5S 3A7
Phone: 416-618-5564
Email: lornea@fowc.ca
www.fowc.ca/

Intestinal Disease Education and Awareness Society (IDEAS)

3885 203A Street
Langley, BC V3A 1W6
(604) 315-3371
E-mail: info@weneedideas.com
www.weneedideas.com

Wounds Canada

(The Canadian Association of Wound Care)
P.O. Box 35569
York Mills Plaza,
North York, ON M2L 2Y4
Tel.: (416) 485-2292
Fax: 1 888 710-2289
Email: info@woundscanada.ca
www.woundscanada.ca

Canadian Home Care Association

2000 Argentinia Road
Plaza 3, Suite 302
Mississauga, ON L5N 1W1
Tel.: (905) 567-7373
Email: chca@cdnhomecare.ca
www.cdnhomecare.ca

Infection Prevention and Control Canada (IPAC)

P.O. Box 46125 RPO Westdale
Winnipeg, MB R3R 3S3
Tel.: 1 866 999-7111 or (204) 897-5990
Fax: (204) 895-9595
Email: info@ipac.canada.org
www.ipac.canada.org

Canadian Dermatology Association

1385 Bank Street, Suite 425
Ottawa, ON K1H 8N4
Tel.: 1 800 267-3376 or (613) 728-1748
Fax: 1 866 267-2178 or (613) 738-4695
www.dermatology.ca

Diabetes Canada

National Office
1400-522 University Ave
Toronto, ON M5G 2R5
Tel.: 1 800 BANTING or (416) 363-3373
Fax: (416) 363-3393
Email: info@diabetes.ca
www.diabetes.ca

Health Canada

Address Locator 0900C2
Ottawa, ON K1A 0K9
Tel.: 1 866 225-0709 or (613) 957-2991
Fax: (613) 941-5366
www.hc-sc.gc.ca
www.canada.ca

**ConvaTec****Autorisation de retour de marchandise**
Authorization to return merchandiseNuméro d'autorisation
OBLIGATOIRE

N°

Authorization number
MANDATORY

Les retours seront acceptés seulement si ce numéro d'autorisation aura été obtenu préalablement du bureau des commandes de ConvaTec Canada au 1-866-987-1471. Returns will only be accepted if this Authorization Number has been obtained from ConvaTec Canada Order Desk at 1-866-987-1471.

NOM DU CLIENT
CUSTOMER NAME _____ N° DE COMPTE
ACCOUNT No. _____ADRESSE
ADDRESS _____CODE POSTAL
POSTAL CODE _____

No DE CODE CODE N°.	QUANT.	DESCRIPTION	DATE DE PÉREMPTION EXP. DATE	N° DE LOT LOT N°.

REMARQUE: TOUTES LES MARCHANDISES RETOURNÉES DOIVENT ÊTRE PORT PAYÉES, SINON LES FRAIS DE TRANSPORT SERONT DÉDUITS DE LA NOTE DE CRÉDIT.
NOTE: ALL RETURNED SHIPMENTS MUST BE PREPAID OTHERWISE TRANSPORTATION CHARGES WILL BE DEDUCTED FROM THE CREDIT NOTE.

COMMENTAIRES:
COMMENTS: _____APPROUVÉ PAR REP: _____ N° DE TERRITOIRE:
APPROVED BY REP: _____ TERRITORY No. _____ DATE: _____
(LETTRES MOULÉES) / (PLEASE PRINT)

SIGNATURE DU CLIENT / CUSTOMER SIGNATURE _____

ADRESSE DE RETOUR: ConvaTec Canada
RETURN ADDRESS: Bureau des commandes / Order Desk
1375, Newton
Boucherville, QC
J4B 5H2

À REMPLIR PAR LE SERVICE DE LA RÉCEPTION / TO BE COMPLETED BY RECEIVING DEPT.		
DATE	NOMBRE DE BOÎTES: NO. OF CTNS:	
PAR: VIA:		
☐ PORT PAYÉ PREPAID	☐ PORT DÛ COLLECT	FRAIS DE TRANSPORT: TRANSPORT CHARGES:
ÉTAT DE LA MARCHANDISE: CONDITION OF GOODS:		
SIGNATURE: SIGNED:		N° DE LA NOTE: MEMO No.:

Explanation of CONVATEC Package Code Numbering System

LOT NUMBERS

Every ConvaTec product carries a lot number that serves to indicate when the item was produced, and from which a complete history as well as expiration date of the product can be derived.

Most products use the following format:

Example - Lot number **9J00978**

9 = Year (2019)

J = Month (September)

00978 = Lot number assignment for positive identification of each lot produced

Expiration Designation:

- If an expiration date is noted on the box, it will appear below the lot number.
- If an expiration date is *not* noted on the box, it can be derived by adding five years to date of manufacture of product. See above for details on manufacture dates.

Example: If date of manufacture is determined as September 2009 (ex. 9J), the product's shelf life will be until September 2014; however, limit to return the product is *six months prior* to this date which would be: March 1, 2014 .

A - January	B - February	C - March	D - April
E - May	F - June	G - July	H - August
J - September	K - October	L - November	M - December

Please contact your ConvaTec Territory Manager or Customer Relations Center at 1 800 465-6302 for clarification of control numbers or expiration dates.

Ostomy Care

ConvaTec ostomy products and services are designed to help people live with a renewed sense of confidence, comfort and peace of mind. It's all about enhancing quality of life.



Natura®/ Natura® ⊕

eakin®

Esteem synergy®

Esteem synergy® ⊕

Active Life®

Little Ones®

AllKare®

Esteem® ⊕

Wound and Skin Care

ConvaTec supports healthcare professionals with a comprehensive approach to healing and patient management by offering a range of life-enhancing solutions.



AQUACEL®
Foam

AQUACEL®
Dressings

AQUACEL®
Surgical

AQUACEL® Ag.
Foam

AQUACEL® Ag.
Dressings

AQUACEL® Ag.
Surgical

AQUACEL®
Extra™

AQUACEL® Ag.
Burn

AQUACEL® Ag+
Extra™

AQUACEL® Ag.
Extra™

FoamLite®
ConvaTec

DuoDERM®
Signal®

DuoDERM®
Extra Thin

KALTOSTAT®

CarboFlex®
Odor Control Dressing

ALOE VESTA®

SENSI-CARE®

SurePress®

Continence and Critical Care

ConvaTec provides innovative nurse-oriented solutions for the management of fecal incontinence for enhanced patient outcomes.

FlexiSeal™ PROTECT

FlexiSeal™
FMS

FlexiSeal™ SIGNAL™
The seal of innovation™

Negative Pressure Wound Therapy

The first system to combine negative pressure wound therapy with Hydrofiber® Technology.

Avelle™
Negative pressure, positive power



For more information, contact our Customer Relations Center (Registered Nurses on staff) at **1-800-465-6302**, Monday through Friday, 8:00 AM to 6:00 PM (EST), or visit our Web Site at **www.convatec.ca**